



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110033 (INDIA)

Phone: 91-8130973844 | E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC

TESTING EXCELLENCE
IN EVERY RESULT

ACCELERATED / COOLING CURING TANK

PRECISION TEMPERATURE CONTROL. RELIABLE RESULTS.

Designed for laboratory curing of concrete and cement specimens under controlled temperature conditions.



**WIDE TEMPERATURE
RANGE**
20°C to Boiling
Temperature



HIGH ACCURACY
±1°C



**DURABLE
CONSTRUCTION**
Built to Last



EASY TO OPERATE
User Friendly
Control Panel



TECHNICAL SPECIFICATIONS

Inside Working Chamber Size	36 inch (L) x 30 inch (W) x 30 inch (H)
Temperature Range	20°C to Boiling Temperature
Temperature Accuracy	±1°C
Temperature Controller	Digital PID Controller
Interior Material	Stainless Steel (SS 304)
Exterior Material	Mild Steel (MS) duly Powder Coated
Insulation	High Density Glass Wool
Heating System	Immersion Heater (SS)
Cooling System	Forced Air / Water Cooling
Water Circulation	High Efficiency Pump for Uniform Temp.
Power Supply	230V AC, 50Hz (Single Phase)
Safety Features	Over Temperature Cutoff, Low Water Level Protection
Applications	Concrete, Cement & Mortar Specimens Curing as per IS / ASTM Standards

INSIDE CHAMBER



APPLICATIONS

- Concrete Cube & Cylinder Curing
- Cement & Mortar Specimen Curing
- Accelerated Curing Tests
- Heat Curing Studies
- Laboratory Research & Quality Control

STANDARD COMPLIANCE

- IS 516, IS 4031 (Part 2)
- ASTM C31, ASTM C192
- ASTM C511
- And other relevant International Standards



SCAN TO VISIT
OUR WEBSITE

ARCHIT ASSOCIATES

A-242, Jahangir Puri, Delhi - 110033

Phone: 8130973844

Email: architasst@gmail.com

Web: www.architassociates.com



PRECISE
TEMPERATURE
CONTROL



UNIFORM
CURING
ENVIRONMENT



RUST RESISTANT
& DURABLE



RELIABLE
PERFORMANCE

BUILT FOR LABS.
ENGINEERED FOR
EXCELLENCE.



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844

E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



SCAN TO VIEW
PRODUCT RANGE

COMPLETE RANGE OF CONCRETE TESTING EQUIPMENT

ACCURATE • RELIABLE • DURABLE

Designed for quality control laboratories, construction sites & research institutions.



ACCURATE
RESULTS



RUGGED
CONSTRUCTION



EASY TO
OPERATE



LOW
MAINTENANCE



IDEAL FOR
LABORATORIES

DIGITAL COMPRESSION TESTING MACHINE



CAPACITY: 2000 kN (200 TON)

HAND OPERATED COMPRESSION TESTING MACHINE



CAPACITY: 2000 kN (200 TON)

FLEXURE TESTING MACHINE



CAPACITY: 200 kN

DIGITAL COMPRESSION TESTING MACHINE

- Suitable for testing cubes, cylinders & beams
- Digital load indicator with high accuracy
- Automatic peak load hold
- High rigidity two column frame
- Safety door with limit switch
- Supplied complete with control unit

TECHNICAL SPECIFICATIONS

Capacity	2000 kN (200 TON)
Least Count	1 kN
Pace Rate	140 kg/cm ² /min (Approx.)
Daylight	330 mm (Approx.)
Test Space	300 mm (Approx.)
Power Supply	230V, 50Hz, Single Phase

HAND OPERATED COMPRESSION TESTING MACHINE

- Suitable for cubes up to 150 mm
- Compact, simple & rugged design
- Accurate & easy to operate
- Fitted with pressure gauge
- Ideal for site & laboratory use

TECHNICAL SPECIFICATIONS

Capacity	2000 kN (200 TON)
Least Count	10 kN
Platen Size	165 mm (Dia)
Daylight	220 mm (Approx.)
Piston Stroke	50 mm (Approx.)
Weight (Approx.)	120 kg

FLEXURE TESTING MACHINE

- Suitable for testing concrete beams (100 x 100 x 500 mm)
- Accurate load measurement with dial gauge
- Two point / Three point loading
- Sturdy & stable construction
- Easy to operate

TECHNICAL SPECIFICATIONS

Capacity	200 kN
Least Count	0.1 kN
Load Indicator	Dial Gauge
Roller Span	400 mm
Roller Dia.	38 mm
Dimensions (LxWxH)	900 x 450 x 950 mm (Approx.)

CONCRETE MOULDS

CUBE MOULD



Made from Cast Iron / Steel
As per IS: 10086

AVAILABLE SIZES

- 70.6 mm
- 100 mm
- 150 mm
- 200 mm

CYLINDRICAL MOULD



Made from Cast Iron / Steel
As per IS: 10086

AVAILABLE SIZES

- 100 x 200 mm
- 150 x 300 mm
- 100 x 100 mm

BEAM MOULD



Made from Cast Iron / Steel
As per IS: 10086

AVAILABLE SIZES

- 100 x 100 x 500 mm
- 150 x 150 x 700 mm
- 150 x 150 x 750 mm

Heavy Duty • Precise Dimensions • Long Life

APPLICATIONS

- ✓ Cement Plants
- ✓ Construction Laboratories
- ✓ Quality Control Labs
- ✓ Research & Development
- ✓ Infrastructure Projects
- ✓ Educational Institutions

ARCHIT ASSOCIATES

A-242, Jahangir Puri,
Near Kali Mandir, Delhi-110 033 (INDIA)

91-8130973844

✉ architasst@gmail.com

🌐 www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



SCAN TO VIEW
PRODUCT RANGE



ARCHIT ASSOCIATES

📍 A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

☎ Phone: 91-8130973844

✉ E-mail: architasst@gmail.com

🌐 Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F22C



ACCURATE
RESULTS



RELIABLE
PERFORMANCE



QUALITY
ASSURED



ASPHALT CONTENT TESTER

300 x 340 x 400 mm

Two Chamber Combustion & Ignition System

Advanced fully automatic Asphalt Content Tester for determination of asphalt content in bituminous mixtures (Hot Extraction Method).



TWO CHAMBER SYSTEM

Combustion chamber and ignition chamber for efficient and complete burning of sample.



HIGH TEMPERATURE CAPABILITY

Ignition chamber temperature 540°C and 750°C for accurate results.



TOUCH DISPLAY

User friendly touch display for easy operation and program selection.



BUILT-IN PRINTER

Automatic print slip after completion of test.



AUTOMATIC WEIGHING BALANCE

Integrated automatic weight balance for precise measurement.



DATA RECORDING

USB port for recording and transferring data.



RUGGED & RELIABLE DESIGN

Robust construction for long life and consistent performance.



TECHNICAL SPECIFICATIONS

Chamber Type	Two Chamber (Combustion & Ignition)
Overall Dimensions (L x W x H)	300 x 340 x 400 mm (approx.)
Ignition Chamber Temperature	540°C & 750°C
Temperature Control	PID controlled
Display	7" Touch Screen Display
Print	Automatic print slip after test
Weighing System	Automatic Weight Balance (built-in)
Power Supply	230V AC, 50 Hz, Single Phase
Power Consumption	1500 W (approx.)
Data Storage	USB for data record
Alarm	Visual & Audible Alarm
Sample Trays	2 Sample Trays + 1 Tray Handle
Body Material	Mild Steel with Powder Coating
Application	Determination of Asphalt Content in bituminous mixes by Ignition Method
Standards	IS: 2386 (Part 2), AASHTO T 308, ASTM D6307 (Ignition Method)

APPLICATIONS

- Bitumen & Asphalt Testing Laboratories
- Quality Control in Road Construction
- Asphalt Mix Design & Evaluation
- Research & Development Institutes
- AMC & Field Testing Laboratories



HIGH ACCURACY
& PRECISION



540°C & 750°C
HIGH TEMPERATURE



EASY OPERATION
TOUCH DISPLAY



AUTO PRINT
RESULTS



AUTOMATIC
WEIGHING



USB PORT
DATA RECORD



SAFE & DURABLE
CONSTRUCTION



ONE YEAR WARRANTY

Against manufacturing defects
(excluding consumables)



AMC & REPAIRING FACILITY AVAILABLE

Comprehensive Annual Maintenance Contract
and repairing services available

ARCHIT ASSOCIATES

📍 A-242, Jahangir Puri,
Near Kali Mandir, Delhi-110 033 (INDIA)

☎ 91-8130973844

✉ architasst@gmail.com

GSTIN NO. 07A0UPG0300F22C

ARCHIT ASSOCIATES – YOUR PARTNER IN TESTING EXCELLENCE



ARCHIT ASSOCIATES

Trusted Supplier of
MATERIAL TESTING & LABORATORY INSTRUMENTS



A-242, Jahangir Puri, Near Kali Mandir,
Delhi-110 033 (INDIA)

91-8130973844, 9871406036,
9810052732, 9810052758

architassst@gmail.com

www.architassociates.com

COMPLETE LAB SOLUTIONS FOR AAC BLOCK MANUFACTURERS IN INDIA



QUALITY
ASSURED



ACCURATE &
RELIABLE



LABORATORY
SOLUTIONS



CUSTOMIZED
EQUIPMENT

- ✓ Advanced Instruments for
- ✓ Quality, Consistency &
- ✓ Compliance in AAC Production



COMPLETE RANGE OF LABORATORY INSTRUMENTS FOR AAC BLOCK MANUFACTURING PLANTS

LAB OVEN



Uniform Heating
Accurate Temp. Control

MOISTURE BALANCE



Fast & Accurate
Moisture Determination

WEIGHING BALANCES



High Precision
Reliable Weighing

MUFFLE FURNACE



High Temperature
Reliable Performance

COMPRESSION TESTING MACHINE



Accurate Strength Testing
of AAC Blocks

FLEXURE TESTING MACHINE



Flexural Strength
Determination

VICAT APPARATUS



Setting Time
Determination

LE-CHATLIER MOULD



Soundness Test
(Expansion Check)

WATER BATH



For Curing &
Temperature Control

SLAKE REACTIVITY APPARATUS



Lime Reactivity Testing
As per Standards

JAW CRUSHER



Crushing of Raw Materials
(Lime, Fly Ash, Gypsum, Sand)

SIEVE SHAKER



Efficient Sieving for
Particle Analysis

BRASS SIEVES 20 cm DIA



High Quality Sieves for
Accurate Analysis

RIFLE SAMPLE DIVIDER



Uniform & Representative
Sample Division

PH METER



Accurate pH Measurement
for Slurry & Mix Water

LAB GLASSWARE



High Quality Glassware
for Laboratory Use

LAB CHEMICALS



Lab Chemicals

MEASURING CYLINDER & BEAKERS



Accurate Volume
Measurement

THERMOMETER (DIGITAL)



- Raw Material
Testing
- Precise Control
& QA Assurance
- Product Quality
Evaluating



Moisture Free
Sample Storage

APPLICATIONS IN AAC BLOCK MANUFACTURING

- ✓ Raw Material Testing (Lime, Fly Ash, Gypsum, Sand)
- ✓ Slurry & Mix Quality Control
- ✓ Green Cake Density & Strength Testing
- ✓ Autoclaved Block Quality Testing
- ✓ Process Control & Quality Assurance
- ✓ Research & Development Laboratories



KEY BENEFITS

- ✓ High Accuracy & Precision
- ✓ Durable & Long Lasting
- ✓ Easy to Operate & Maintain
- ✓ Compliance with International Standards
- ✓ Ideal for AAC Quality Control Laboratories



Raw Material
Testing



Slurry & Mix
Quality Control



Green & Hardened
Block Testing



Process Control
& QA Assurance



Product Quality
Evaluation



R & D
Laboratories

A-242, Jahangir Puri,
Near Kali Mandir,
Delhi-110 033 (INDIA)

91-8130973844, 9871406036,
9810052732, 9810052758

architassst@gmail.com

www.architassociates.com

ALL PRODUCTS COMPLY WITH RELEVANT INDUSTRY STANDARDS



ARCHIT ASSOCIATES

YOUR TRUSTED PARTNER IN LABORATORY & MATERIAL TESTING

📍 A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

☎ +91-8130973844, 9871406036 ✉ architasst@gmail.com

🌐 www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



ACCURATE
RESULTS



RELIABLE
PERFORMANCE



QUALITY
ASSURED

TESTING EXCELLENCE IN EVERY RESULT

ACCELERATED POLISHING MACHINE

REF. STD. BS 812-114

Accelerated Polishing Machine is used for wear resistance testing of concrete by polishing the specimen surfaces with specified abrasive under controlled load, speed and water supply. It simulates the effect of traffic wear on concrete surfacing. The machine is robustly built for long life and consistent performance in laboratories.

FEATURES

- Complies with BS 812-114
- Dual Abrasive Feeders For Corn Emery & Flour Emery
- Controlled Abrasive Feed Rate
- Controlled Water Supply
- Spring Loaded Rubber Tyre
- Standard Road Wheel Speed
- Applied Load as per Standard
- Large Removable Tray For collection of used water & abrasives
- Thermal Cutout Protection
- Emergency Stop Button
- Robust & Heavy Duty Construction
- Operates on 220V, Single Phase, AC Supply

STANDARD ACCESSORIES

- Road Wheel – 1 No.
- Rubber Tyre – 1 No.
- Abrasive Feeders (Corn & Flour) – 2 Nos.
- Water Container with Valve – 1 No.
- Measuring Cylinder – 1 No.
- Spanners & Tools – 1 Set
- Instruction Manual – 1 No.



TECHNICAL SPECIFICATIONS

Standard:	BS 812-114
Road Wheel Speed	315 to 325 rpm
Applied Load on Specimen	715 to 735 N
Abrasive	Corn Emery & Flour Emery
Feed Rate (Corn Emery)	20 – 34 g/min
Feed Rate (Flour Emery)	2 – 4 g/min
Water Supply (Corn Emery)	20 – 34 ml/min
Water Supply (Flour Emery)	5 – 8 ml/min
Specimen Size	As per BS 812-114
Tyre	Solid Rubber, Spring Loaded
Electrical Supply	220 V, 50 Hz, Single Phase, AC
Motor	0.5 HP, 1440 rpm
Safety	Thermal Cutout & Emergency Stop
Overall Dimensions (L x W x H)	750 x 600 x 1100 mm (Approx.)
Net Weight (Approx.)	120 – 140 kg (Approx.)

PRINCIPLE OF OPERATION

Specimens are manufactured in accurately machined and matched moulds; they are then removed from the moulds and located on the 'Road Wheel'. The wheel is rotated and comes in contact with the solid rubber tyre, spring loaded. Abrasive charges are continuously fed by mechanical feeders at fixed speed. Feed rates for corn emery are 20–34 g/min and for flour emery 2–4 g/min. Water is supplied at a controlled rate of 20–34 ml/min (corn emery) and 5–8 ml/min (flour emery) through a water container. Used water and abrasives are collected in a large removable tray. Electrical supply is protected by thermal cutouts and an emergency stop button.

APPLICATIONS



CONCRETE LABS



R & D CENTRES



QUALITY CONTROL LABORATORIES



CIVIL ENGINEERING LABORATORIES

OPTIONAL ACCESSORIES

- Digital Tachometer
- Additional Road Wheel
- Additional Rubber Tyre
- Abrasive (Corn Emery & Flour Emery)
- Stopwatch
- Extra Measuring Cylinder
- Spare Water Valve
- Calibration Certificate



PRECISE RESULTS
Consistent & accurate
testing performance



HEAVY DUTY
CONSTRUCTION
For long life



EASY TO OPERATE
User friendly design
& controls



TEST CERTIFICATES
Supplied with
calibration certificates



AFTER SALES SUPPORT
Prompt service &
technical support



SCAN TO VISIT
OUR WEBSITE

ARCHIT ASSOCIATES

📍 A-242, Jahangir Puri, Delhi - 110033 (INDIA)

☎ +91-8130973844, +91 9871406036

✉ architasst@gmail.com

28
YEARS OF
EXPERIENCE



BUILT FOR LABS.
ENGINEERED FOR EXCELLENCE.
Precision, Performance, Trust.



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)
Phone: 91-8130973844 E-mail: architasst@gmail.com
Web: www.architassociates.com
GSTIN NO. 07A0UPG0300F22C



AGGREGATE IMPACT TESTER

IS: 2386 (Part 4)

The Aggregate Impact Tester is used to determine the impact value (toughness) of aggregates. The instrument consists of a metal cup mounted in a base, a hammer of specified weight which is released from a fixed height and an arrangement to raise the hammer. The test results are used in the design of road surfacing and other construction works.

FEATURES



STURDY CONSTRUCTION

Made from high quality materials for durability and long service life.



PRECISION HAMMER

Hammer of 13.5 – 14.0 kg released from a height of 380 ± 5 mm.



STANDARD COMPLIANCE

Designed and manufactured as per IS: 2386 (Part 4).



EASY OPERATION

Simple mechanism for quick and accurate testing.



ACCURATE RESULTS

Provides reliable and repeatable impact value of aggregates.



TECHNICAL SPECIFICATIONS

PARAMETER	SPECIFICATION
Standard	IS: 2386 (Part 4)
Hammer Weight	13.5 kg to 14.0 kg
Hammer Diameter	102 ± 1 mm
Height of Fall	380 ± 5 mm
Diameter of Cup	102 ± 1 mm
Depth of Cup	50 ± 1 mm
Base Diameter	150 mm (Min.)
Overall Dimensions (Approx.)	450 mm (L.) x 300 mm (W) x 750 mm (H)
Net Weight (Approx.)	25 – 30 kg

STANDARD ACCESSORIES

- Metal Cup
- Hammer
- Base Plate
- Tamping Rod
- Measure Scale
- Lifting Guide
- Instruction Manual

OPTIONAL ACCESSORIES

- Extra Metal Cup
- Digital Height Gauge
- Carrying Case



APPLICATIONS



Road Construction



Quality Control Laboratories



Research & Development



Educational Institutions

KEY BENEFITS



Reliable & Accurate Testing



Rugged & Long Lasting



Easy to Operate & Maintain



Compliance with Indian Standards



Ideal for Laboratory & Field Use

ARCHIT ASSOCIATES

A-242, Jahangir Puri,
Near Kali Mandir, Delhi-110 033 (INDIA)
91-8130973844, 9871406036

✉ architasst@gmail.com

🌐 <https://sites.google.com/view/architassociates>

GSTIN NO. 07A0UPG0300F22C

ARCHIT ASSOCIATES – YOUR PARTNER IN TESTING EXCELLENCE



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844 | E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC

TESTING
EXCELLENCE IN
EVERY RESULT


COMPLETE RANGE OF CIVIL MATERIAL TESTING EQUIPMENT

ONE STOP SOLUTION FOR ALL YOUR TESTING NEEDS

For Soil, Cement, Concrete, Bitumen, Paint, Water, Wood & General Lab Testing



PREMIUM
QUALITY

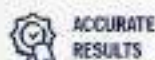

ACCURATE
RESULTS



RELIABLE &
DURABLE



EASY TO USE
& MAINTAIN



COMPLIANT TO
IS / ASTM STANDARDS

SOIL TESTING EQUIPMENT



- Standard Proctor Compaction Test Apparatus
- CBR Testing Machine
- Direct Shear Test Apparatus
- Consolidation Test Apparatus
- Sieve Shaker, Sieve Set
- Field Density, Sand Replacement Apparatus
- Adair's Limit, Permeability Test Apparatus
- Moisture Content, Core Cutter, etc.

CEMENT TESTING EQUIPMENT



- Vicat Apparatus
- Le Chatelier Apparatus
- Cement Compression Testing Machine
- Fineness of Cement Apparatus
- Soundness Test Apparatus
- Setting Time Test Apparatus
- Specific Gravity, Autoclave, etc.

CONCRETE TESTING EQUIPMENT



- Compression Testing Machine (CTM)
- Flexure Testing Machine
- Slump Test Apparatus
- Cube, Cylinder & Beam Molds
- Vibrating Table
- Compacting Factor Apparatus
- Rebound Hammer, Core Cutter, etc.

BITUMEN TESTING EQUIPMENT



- Ring and Ball Apparatus (Semi Automatic & Digital)
- Bitumen Penetrometer
- Softening Point Apparatus
- Ductility Testing Machine
- Flash & Fire Point Apparatus
- Marshall Stability Test Apparatus
- Bitumen Extraction Apparatus, etc.

PAINT & COATING TESTING EQUIPMENT



- Viscosity Cup & Viscometer
- Drying Time Recorder
- Gloss Meter
- Paint Shear
- Wet Film Thickness Gauge
- Scratch Hardness Tester
- Elcometer, Adhesion Tester, etc.

WATER TESTING EQUIPMENT



- pH Meter
- TDS Meter, Conductivity Meter
- Turbidity Meter
- Jar Test Apparatus
- BOD Incubator
- Water Distillation Unit
- Water Bath, COD Apparatus, etc.

WOOD TESTING EQUIPMENT



- Universal Testing Machine for Wood
- Moisture Meter
- Density Determination Kit
- Hardness Tester (Brinell)
- Thickness Planner
- Compression & Bending Test Fixtures
- Adhesion & Bond Strength Tester, etc.

GENERAL LAB TESTING EQUIPMENT



- Hot Air Oven & Muffle Furnace
- Electronic Balance
- Magnetic Stirrer
- Mechanical Shaker
- Hot Plate
- Incubator
- Desiccator, pH Paper, Glassware, etc.

FEATURED PRODUCTS

COMPRESSION TESTING MACHINE (CTM)



For Concrete Cube, Cylinder, Beam & Block Testing

RING & BALL APPARATUS - DIGITAL



For Softening Point Determination of Bitumen Materials

CBR TESTING MACHINE



For Subgrade & Base Course Evaluation

STANDARD PROCTOR COMPACTION APPARATUS



For Maximum Dry Density & Optimum Moisture Content

FLEXURE TESTING MACHINE



For Concrete Beams Bending Strength

VICAT APPARATUS



For Setting Time & Consistency of Cement

DIGITAL pH & TDS METER



For Water Quality Analysis

INDUSTRIES WE SERVE



CONFORMS TO IS / ASTM STANDARDS



SOIL
IS 2720 (Part I) to XI
ASTM D422, D598
D4318, D1533



CEMENT
IS 4031 (Part 4,5,6)
ASTM C187, C191
C204, C349



CONCRETE
IS 516, IS 1199
ASTM C39, C78
C293, D699



BITUMEN
IS 1205, IS 73
ASTM D5, D36
D133, D4402



WATER
APHA, IS 3025
ASTM D1067
D1888



SCAN TO VISIT
OUR WEBSITE

ARCHIT ASSOCIATES

A-242, Jahangir Puri, Delhi - 110033

8130973844, 9871406036

architasst@gmail.com



WIDE RANGE
OF PRODUCTS



ACCURATE
TEST RESULTS



AFTER SALES
SUPPORT



TRUSTED BY
PROFESSIONALS

ACCURATE TESTING. RELIABLE RESULTS. BUILT FOR LABS. ENGINEERED FOR EXCELLENCE.



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844 | E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC

TESTING
EXCELLENCE IN
EVERY RESULT

LABORATORY BALL MILL

5KG, 10KG & 20KG CAPACITY

Used for grinding and homogenizing of bituminous materials, minerals, ores, chemicals and other lab samples to the required fineness.

FEATURES & BENEFITS



Robust Construction
Heavy duty MS fabricated body for vibration free performance.



Powerful Motor
High torque motor ensures smooth and efficient operation.



Efficient Grinding
Uniform grinding action for consistent and repeatable results.



Safety Features
Secure lid with safety interlock & overload protection.



User Friendly Operation
Easy to operate with minimum maintenance.



Wide Application
Ideal for bitumen, minerals, ores, chemicals, soils and other lab materials.



AVAILABLE CAPACITIES

5KG CAPACITY



Drum Volume (Approx.)
12 Liters
Ball Charge (Approx.)
5 kg
Suitable for small batch laboratory testing.

10KG CAPACITY



Drum Volume (Approx.)
20 Liters
Ball Charge (Approx.)
10 kg
Ideal for medium batch laboratory operations.

20KG CAPACITY



Drum Volume (Approx.)
40 Liters
Ball Charge (Approx.)
20 kg
Suitable for large batch laboratory applications.

APPLICATION

- Bitumen & Asphalt Testing
- Minerals & Ores
- Cement & Concrete Materials
- Soils & Aggregates
- Chemical & Pharmaceutical Labs
- Research & Development



CONFORMS TO RELEVANT IS / ASTM STANDARDS

- IS 1205 Part II
- ASTM D36
- AASHTO T53



**DURABLE
CONSTRUCTION**



**HIGH
PERFORMANCE**



**PRECISE
RESULTS**



**RELIABLE
& EFFICIENT**



SCAN TO VISIT
OUR WEBSITE

ARCHIT ASSOCIATES

A-242, Jahangir Puri, Delhi - 110033

Phone: 8130973844, 9871406036

Email: architasst@gmail.com

QUALITY INSTRUMENTS

ACCURATE TESTING | RELIABLE RESULTS
BUILT FOR LABS, ENGINEERED FOR EXCELLENCE



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844

E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



MIXER BITUMEN

CAPACITY : 5 LITERS, 10 LITERS, 20 LITERS

For Bituminous Mix Design & Laboratory Use

Mixer Bitumen is designed for uniform mixing of bitumen with aggregates in bituminous mix design. The machine ensures thorough and homogeneous mixing with variable speed control to meet different testing standards.

It is widely used in road construction laboratories, quality control labs and research institutions.

FEATURES



RUGGED CONSTRUCTION

Heavy duty structure for long life and reliable performance.



VARIABLE SPEED CONTROL

Variable speed for proper and uniform mixing of bitumen and aggregates.



EFFICIENT MIXING

Ensures homogeneous mixing for accurate and reproducible test results.



SAFETY FEATURES

Fitted with safety guard and limit switch for safe operation.



STAINLESS STEEL BOWL

High grade stainless steel bowl for corrosion resistance and easy cleaning.



EASY OPERATION

User friendly operation with reliable performance.



TECHNICAL SPECIFICATIONS

Parameter	Specification
Capacity Options	5 Liters, 10 Liters, 20 Liters
Bowl Material	Stainless Steel
Mixing Type	Planetary
Speed Control	Variable Speed
Speed Range	60 to 150 rpm (Approx.)
Motor	0.5 HP / 1 HP (Approx.) (Depending on Capacity)
Power Supply	220V, 50 Hz, Single Phase, AC Supply
Application	Bituminous Mix Design, Laboratory Testing

STANDARD ACCESSORIES



Stainless Steel Bowl



Mixing Paddle



Scraper

* Accessories may vary with capacity model.

KEY BENEFITS

- Uniform and consistent mixing
- Improves test accuracy and reliability
- Saves time and manual effort
- Low maintenance and easy cleaning
- Suitable for laboratory and research use

APPLICATIONS

- Bituminous Mix Design
- Road Construction Laboratories
- Quality Control Laboratories
- Research & Development
- Educational Institutions



RUGGED & DURABLE



UNIFORM MIXING



VARIABLE SPEED



EASY TO OPERATE



CORROSION RESISTANT



SAFE & RELIABLE



LABORATORY USE

ARCHIT ASSOCIATES

A-242, Jahangir Puri,
Near Kali Mandir, Delhi-110 033 (INDIA)
91-8130973844, 9871406036

architasst@gmail.com

<https://sites.google.com/view/architassociates>

GSTIN NO. 07A0UPG0300F2ZC

ARCHIT ASSOCIATES – YOUR PARTNER IN TESTING EXCELLENCE



ARCHIT ASSOCIATES

YOUR TRUSTED PARTNER IN LABORATORY & MATERIAL TESTING

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

91-8130973844

architasst@gmail.com

www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



ACCURATE
RESULTS



RELIABLE
PERFORMANCE



QUALITY
ASSURED

DIGITAL BOMB CALORIMETER

ACCURATE DETERMINATION OF CALORIFIC VALUE OF SOLID & LIQUID FUELS

Digital Bomb Calorimeter is used for the precise determination of calorific value (Gross Calorific Value) of solid and liquid fuels. It is designed in compliance with the specifications of Institute of Petroleum (IP) and British Standard Institution (BSI).



MICROPROCESSOR
BASED DIGITAL
TEMPERATURE
INDICATOR



BUILT-IN PRINTER
FOR DIRECT
PRINTOUT



STIRRER FOR
UNIFORM TEMPERATURE
DISTRIBUTION



EASY TO OPERATE
& MAINTAIN



APPLICATIONS



Coal & Coke



Biomass &
Wood



Liquid Fuels



Petroleum &
Chemical Industry



Research &
Laboratories



HIGH ACCURACY
& RESOLUTION



LONG LASTING
PERFORMANCE



EASY OPERATION
& MAINTENANCE



TEST CERTIFICATES
PROVIDED



COMPLIANCE WITH
IP & BSI STANDARDS

SALIENT FEATURES

- Microprocessor based digital temperature indicator
- Built-in printer for direct printout
- High accuracy and resolution
- Stirrer for uniform temperature distribution
- Easy to operate and maintain
- Robust construction for long life
- Supplied with test certificates

TECHNICAL SPECIFICATIONS

Temperature Range	: Ambient to 40°C
Temperature Accuracy	: $\pm 0.1^\circ\text{C}$
Resolution	: 0.01°C
Bomb Volume	: 300 ml (approx.)
Calorimeter Capacity	: 4.5 Liters (approx.)
Power Supply	: 230V AC, 50Hz, Single Phase
Standard	: IP & BSI

SUPPLIED COMPLETE WITH

- Calorimeter with Insulation
- Bomb Assembly (300 ml)
- Bomb Jacket
- Stirrer with Motor
- Temperature Probe
- Digital Temperature Indicator
- Control Unit with Printer
- Pressure Gauge
- Firing Wire, Ignition Coil
- Accessories & Tools



SCAN TO VISIT
OUR WEBSITE

ARCHIT ASSOCIATES

A-242, Jahangir Puri, Delhi - 110033 (INDIA)

91-8130973844, 9871406036

architasst@gmail.com

<https://sites.google.com/view/architassociates>

28
YEARS OF
EXCELLENCE

DELIVERING QUALITY EQUIPMENT
FOR ACCURATE RESULTS AND
BETTER PERFORMANCE.



ARCHIT ASSOCIATES

📍 A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

☎ Phone: 91-8130973844

✉ E-mail: architasst@gmail.com

🌐 Web: www.architassociates.com

GSTIN NO. 07AOUPG0300F2ZC



ACCURATE
RESULTS



RELIABLE
PERFORMANCE



QUALITY
ASSURED



CARBONATION CHAMBER

TEMPERATURE, HUMIDITY & CO₂ CONTROLLED

STANDARD CHAMBER SIZE: 600 x 600 x 600 mm
(CUSTOM SIZES AVAILABLE)

Carbonation chamber is designed to accelerate the carbonation process in concrete specimens under controlled temperature, humidity and CO₂ concentration. It is widely used in research and quality control laboratories to evaluate the durability and performance of concrete.

KEY FEATURES



PRECISE ENVIRONMENT CONTROL

Accurate control of temperature, humidity and CO₂ concentration for consistent results.



WIDE RANGE CONTROL

Temperature 10°C to 50°C, Humidity 50% to 95% and CO₂ 1% to 10%.



RUGGED & DURABLE DESIGN

Heavy-duty insulated chamber with corrosion resistant construction for long life.



DIGITAL CONTROLLER

Microprocessor based PID controller with digital display for easy monitoring and operation.



UNIFORM AIR CIRCULATION

Forced air circulation ensures uniform distribution of temperature, humidity and CO₂.



CUSTOMIZABLE

Chamber size and specifications can be modified as per customer requirement.



TECHNICAL SPECIFICATIONS

Temperature Range	10°C to 50°C
Temperature Accuracy	±1°C
Humidity Range	50% to 95% RH
Humidity Accuracy	±3% RH
CO ₂ Concentration Range	1% to 10%
CO ₂ Accuracy	±0.1%
Chamber Size (Standard)	600 x 600 x 600 mm (L x W x H) (Custom sizes available)
Controller	Microprocessor based PID controller
Display	Digital LED Display (Temperature, Humidity, CO ₂)
Sensor	High precision Temperature, Humidity & CO ₂ sensor
Material (Interior)	Stainless Steel
Material (Exterior)	Mild Steel with Powder Coating
Insulation	High density PUF insulation
Air Circulation	Forced air circulation with blower
Door	Single door with viewing window & gasket seal
Power Supply	230V AC, 50 Hz, Single Phase
Power Consumption	1500 W (approx.)
Gas Supply	CO ₂ from cylinder with pressure regulator
Safety Features	Over temperature protection, Alarm, Gas leak protection (optional)

APPLICATIONS

- Durability testing of concrete structures
- Research & development laboratories
- Evaluation of concrete mixes
- Quality control in construction projects
- Academic & testing institutes



TEMPERATURE
CONTROL
10°C to 50°C



HUMIDITY
CONTROL
50% to 95%



CO₂ CONTROL
1% to 10%



UNIFORM AIR
CIRCULATION



DIGITAL DISPLAY
& PID CONTROL



SAFE & RELIABLE
OPERATION



CUSTOM OPTIONS
AVAILABLE



ONE YEAR WARRANTY

Against manufacturing defects
(excluding consumables)



AMC & REPAIRING FACILITY AVAILABLE

Comprehensive Annual Maintenance Contract
and repairing services available

ARCHIT ASSOCIATES

📍 A-242, Jahangir Puri,
Near Kali Mandir, Delhi-110 033 (INDIA)
☎ 91-8130973844

✉ architasst@gmail.com

GSTIN NO. 07AOUPG0300F2ZC

ARCHIT ASSOCIATES – YOUR PARTNER IN TESTING EXCELLENCE



ARCHIT ASSOCIATES

YOUR TRUSTED PARTNER IN LABORATORY & MATERIAL TESTING

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

+91-8130973844, 9871406036 | architasst@gmail.com

www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



ACCURATE
RESULTS



RELIABLE
PERFORMANCE



QUALITY
ASSURED

TESTING EXCELLENCE IN EVERY RESULT

CALIFORNIA BEARING RATIO (CBR) TEST APPARATUS

50 kN CAPACITY

California Bearing Ratio (CBR) Test Apparatus is designed for performing laboratory evaluation of the California Bearing Ratio value of highway sub bases and subgrade and for the determination of strength of cohesive materials which have maximum particle sizes less than 19 mm (3/4").

The machine is designed to be mounted on a suitable bench and comprises of a robust and compact two column frame with adjustable upper cross beam. The frame is 50 kN capacity.

Three test Rate of Strain speeds are provided 1.25 mm/min, 1.50 mm/min & 2.5 mm/min tests.

KEY FEATURES

- Robust and compact two column frame
- Adjustable upper cross beam
- 50 kN capacity
- Three test rate of strain speeds 1.25 mm/min, 1.50 mm/min & 2.5 mm/min
- Suitable for subgrade, sub-base & cohesive soils
- Easy to operate and maintain
- Precision and reliable performance

STANDARDS

IS:9669, IS: 2720 (Part XVI),
BS 1377; 1924; ASTM D 1883

TECHNICAL SPECIFICATIONS

Capacity	50 kN
Frame Type	Two Column
Upper Cross Beam	Adjustable
Test Rate of Strain	1.25 mm/min, 1.50 mm/min & 2.5 mm/min
Proving Ring Capacity	50 kN
Dial Gauge	0 - 25 mm x 0.01 mm
Power Supply	230V AC, 50 Hz
Motor Power	0.5 HP
Overall Dimensions (L x W x H)	750 x 500 x 1250 mm (Approx.)
Net Weight (Approx.)	120 kg (Approx.)

STANDARD ACCESSORIES

- Mould (152 mm ID x 177.8 mm Height) - 1 No.
- Perforated Base Plate - 1 No.
- Spacer Disc - 1 No.
- Compaction Rammer (4.89 kg) - 1 No.
- Dial Gauge (0 - 25 mm x 0.01 mm) - 1 No.
- Proving Ring (50 kN) - 1 No.
- Tripod for Dial Gauge - 1 No.
- Surcharge Weight - 2.5 kg x 1 No.
- Instruction Manual - 1 No.

OPTIONAL ACCESSORIES

- Penetration Plunger (50 mm, 100 mm)
- Extension Rods
- Digital Dial Gauge
- Data Acquisition System



PRECISE RESULTS
Consistent & accurate
testing performance



SUPERIOR QUALITY
High grade components
for reliability



EASY TO OPERATE
User friendly design
& simple controls



TEST CERTIFICATES
Supplied with
calibration certificates



AFTER SALES SUPPORT
Prompt service &
technical support



SCAN TO VISIT
OUR WEBSITE

ARCHIT ASSOCIATES

A-242, Jahangir Puri, Delhi - 110033 (INDIA)

+91-8130973844, +9871406036

architasst@gmail.com | <https://bit.ly/google.com/view/architassociates>

28
YEARS OF
EXPERIENCE



BUILT FOR LABS.
ENGINEERED FOR EXCELLENCE.
Precision. Performance. Perfection.



ARCHIT ASSOCIATES

INSTRUMENTATION | MEASUREMENT | SOLUTIONS



A-242, Jahangir Puri,
Near Kali Mandir, Delhi-110 033 (INDIA)



+91 8130973844



architasst@gmail.com



www.architassociates.com

GSTIN NO. 07AOUPG0300F2ZC

PRECISION • PERFORMANCE • TRUST

POSITECTOR® 6000 COATING THICKNESS GAUGE

Accurate. Durable. Dependable.

PosiTector 6000 is a robust, easy-to-use coating thickness gauge designed for professionals who demand accuracy and reliability in the most challenging environments.

PRODUCT HIGHLIGHTS



High Accuracy

$\pm(2 \mu\text{m} + 3\% \text{ of reading})$
 $\pm(0.1 \text{ mils} + 3\% \text{ of reading})$



Wide Measuring Range

50 - 3800 microns (Cabled Probe)
2 to 150 mils



Non-Metallic Coating Thickness Gauge

Ideal for measuring coating on non-metallic substrates



Supplied with Manufacturers
Calibration Certificate
[DeFelsko USA]



Durable Design

Impact resistant, dustproof and
water-resistant – built for tough jobs



Built-in Memory

Stores up to 1,000 readings – view,
download or transfer easily



User Friendly

Large 2.8" color touchscreen with
redesigned keypad for quick navigation



WHAT'S IN THE BOX

- PosiTector 6000 Standard Body (Model: BDYSTD-E)
- PosiTector 200 Coating Thickness Gauge – Cabled Probe (Model: PRB200C)
- Calibration Certificate (DeFelsko USA)
- Protective Rubber Holster with Belt Clip
- 3 AAA Batteries
- Wrist Strap
- Instruction Manual
- Carrying Case



Automotive



Paint & Coatings



Manufacturing



Inspection



Quality Control

Whether in the field or in the lab, trust
PosiTector 6000 for accurate coating
measurements every time.

MEASURE WITH CONFIDENCE.

Authorized Distributor & Solution Partner for DeFelsko USA Products

DeFelsko®
The Measure of Quality



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844

E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



COD DIGESTER

SUITABLE FOR 15 GLASS TUBES

For Water & Wastewater Analysis

The COD Digester is designed for the safe and efficient digestion of water and wastewater samples in accordance with Standard Methods for the Examination of Water and Wastewater. The unit provides uniform heating and precise temperature control for accurate and repeatable results.

FEATURES



PRECISE TEMPERATURE CONTROL

Digital PID temperature controller ensures accurate temperature control from ambient +5°C to 150°C.



DIGITAL TIMER

Built-in digital timer for setting digestion time with audible alarm.



SUITABLE FOR 15 TUBES

Comes complete with 15 glass tubes, condensers and stainless steel rack.



SAFE & RELIABLE

Heater with overload protection and stainless steel top for durability and safety.



NABL CALIBRATED

Supplied duly NABL calibrated for temperature accuracy.



TECHNICAL SPECIFICATIONS

Parameter	Specification
Model	COD Digester - 15 Tubes
Capacity	15 Glass Tubes with Condensers
Temperature Range	Ambient +5°C to 150°C
Temperature Controller	Digital PID Controller-cum-indicator
Temperature Accuracy	± 0.5°C
Timer	Digital Timer with Alarm
Heating System	High quality immersion heaters
Rack Material	Stainless Steel
Top Plate Material	Stainless Steel
Included Accessories	15 Glass Tubes, 15 Condensers, Stainless Steel Rack
Power Supply	220V, 50 Hz, Single Phase, AC Supply

STANDARD ACCESSORIES



Stainless Steel Rack



15 Glass Tubes



15 Condensers

KEY BENEFITS

- Uniform heating for consistent results
- Accurate temperature & time control
- Easy operation and maintenance
- Robust construction for long life
- Ideal for laboratory use

APPLICATIONS

- Chemical Oxygen Demand (COD) analysis
- Water & Wastewater Treatment Plants
- Environmental Laboratories
- Industrial Effluent Testing
- Educational & Research Institutions



ACCURATE RESULTS



MICROPROCESSOR CONTROLLED



DIGITAL DISPLAY & TIMER



NABL CALIBRATED



ROBUST & DURABLE



EASY TO USE



LABORATORY USE

ARCHIT ASSOCIATES

A-242, Jahangir Puri,
Near Kali Mandir, Delhi-110 033 (INDIA)
91-8130973844

architasst@gmail.com

<https://sites.google.com/view/architassociates>

GSTIN NO. 07A0UPG0300F2ZC

ARCHIT ASSOCIATES – YOUR PARTNER IN TESTING EXCELLENCE



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844 | E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F22C

TESTING TODAY
BUILDING TOMORROW



ACCURATE



RELIABLE



DURABLE



TRUSTED

COMPLETE RANGE OF CIVIL MATERIAL TESTING EQUIPMENT

ONE STOP SOLUTION FOR ALL YOUR TESTING NEEDS

High quality, precision engineered equipment for testing of

SOIL • CEMENT • CONCRETE • BITUMEN • PAINT • WATER • WOOD & MORE



SOIL TESTING EQUIPMENT



TRIAxIAL SHEAR
TEST APPARATUS



COMPACTION
APPARATUS



SIEVE ANALYSIS
SET



CBR TEST
APPARATUS



LIQUID LIMIT
APPARATUS

CEMENT TESTING EQUIPMENT



VICAT
APPARATUS



LE CHATELIER
MOULD



AUTOMATIC
CONSISTENCY
APPARATUS



CEMENT COMPRESSIVE
STRENGTH MACHINE

CONCRETE TESTING EQUIPMENT



CONCRETE COMPRESSION
TESTING MACHINE



FLEXURAL STRENGTH
TESTING MACHINE



SLUMP CONE
APPARATUS



CONCRETE VIBRATING
TABLE



CURE MOULD &
CYLINDER MOULD

BITUMEN TESTING EQUIPMENT



RING & BALL
APPARATUS
(DIGITAL)



RING & BALL
APPARATUS
(SEMI AUTOMATIC)



LABORATORY
BALL MILL



BITUMEN EXTRACTION
APPARATUS

PAINT TESTING EQUIPMENT



VISCOSITY CUP



FINENESS OF
GRIND GAUGE



ELOCOMETER



PULL OFF
ADHESION TESTER



DRYING TIME
RECORDER

WATER TESTING EQUIPMENT



pH METER



TURBIDITY
METER



TDS METER



WATER BATH



ODI
INCUBATOR

WOOD TESTING EQUIPMENT



WOOD MOISTURE
METER



UNIVERSAL TESTING
MACHINE (WOOD)



WOOD COMPRESSION
TEST APPARATUS



WOOD CORE
CUTTER

OTHER ESSENTIAL TESTING EQUIPMENT



AGGREGATE CRUSHING
VALUE APPARATUS



LOS ANGELES
ABRASION MACHINE



SPECIFIC GRAVITY
APPARATUS



DIGITAL WEIGHING
BALANCE



MUFFLE FURNACE



HOT AIR OVEN



PLATE LOAD
TEST APPARATUS



MARSHALL STABILITY
TEST APPARATUS



FIELD DENSITY
TEST KIT

CUSTOMIZED SOLUTIONS

We also offer customized testing solutions as per client requirements and international standards.

- Custom Design
- Logo Branding
- Special Size & Capacity
- Turnkey Lab Setup

OUR PROMISE




PREMIUM QUALITY
PRODUCTS




ACCURATE &
RELIABLE RESULTS



LONG LASTING &
DURABLE

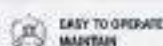

EASY TO OPERATE &
MAINTAIN

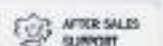

AFTER SALES
SUPPORT

APPLICATION AREAS

- Research & Development Labs
- Consulting & Testing Laboratories
- Roads, Bridges, Dams & Buildings
- Quality Control Laboratories
- Government & Private Institutions
- Industrial & Manufacturing Units



SCAN TO VISIT
OUR WEBSITE

ARCHIT ASSOCIATES

A-242, Jahangir Puri, Delhi - 110033

Phone: 8130973844

Email: architasst@gmail.com

QUALITY EQUIPMENT. TRUSTED RESULTS.
STRONGER INFRASTRUCTURE.

For Better Testing Today, For Better Tomorrow.



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844

E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



LABORATORY CONCRETE MIXER

DURABLE • EFFICIENT • RELIABLE

Electrically operated concrete mixer designed for laboratory use. Suitable for mixing concrete, mortar and other construction materials in laboratories and quality control departments.



FEATURES



Electrically Operated
Smooth and efficient performance.



Swivel Arrangement
Drum can be tilted for easy charging and discharge.



Standard Make Motor
Fitted with 1 HP motor for reliable operation.



Compact & Sturdy Design
Robust construction for long service life.



Easy to Operate
User friendly and low maintenance.

TECHNICAL SPECIFICATIONS

Drum Volume (Approx.)	1 cu.ft.
Operation	Electrically Operated
Drum Type	Swivel Arrangement
Motor	1 HP, Standard Make
Power Supply	230V AC, 50Hz, Single Phase
Drum Material	Heavy Duty Steel
Frame	Robust Mild Steel Construction
Mobility	2 Fixed Wheels for Easy Movement
Note	Supplied without Calibration Certificate

IDEAL FOR



Civil Engineering
Laboratories



Construction
Laboratories



Quality Control
Departments



Research &
Development



Educational
Institutions



**ROBUST
CONSTRUCTION**



**EFFICIENT
PERFORMANCE**



**EASY
OPERATION**



**LONG LASTING
DURABILITY**



**RELIABLE
& SAFE**

ARCHIT ASSOCIATES – YOUR PARTNER IN TESTING EXCELLENCE



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kail Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844

E-mail: architassst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F22C



ACCURATE
RESULTS



RELIABLE
PERFORMANCE



QUALITY
ASSURED



SCAN TO VIEW
PRODUCT RANGE

CONCRETE PENETROMETER SPRING TYPE

AS PER IS:1199 (PART-7)

Used to determine the approximate strength of concrete by measuring the resistance to penetration of a spring loaded plunger. It is used for determining the initial and final setting times of cement and concrete.

CONSTRUCTION & FEATURES

- ✓ Cylindrical spring housing with plunger on the top
- ✓ Penetration needle attached to the other end
- ✓ Plunger graduated in 1 kg divisions (to a maximum capacity of 60 kg)
- ✓ Set of six needle points with areas of 645, 323, 161, 65, 32 & 16 mm²
- ✓ Supplied complete in a wooden carrying case
- ✓ All parts made from corrosion resistant material

ACCURATE
& EASY TO
USE

APPLICATION

A representative sample of fresh concrete is sieved and the mortar is placed in a container. At regular time intervals, the resistance of the mortar to penetration is measured to estimate the initial and final setting times.



TECHNICAL SPECIFICATIONS

Standard	IS:1199 (Part-7)
Type	Spring Type
Plunger Graduated	1 kg divisions up to 60 kg capacity
Spring Housing	Cylindrical type
Penetration Needle	Attached to other end of spring housing
Needle Points (Areas)	645, 323, 161, 65, 32 & 16 mm ² (Set of 6 Nos.)
Maximum Capacity	60 kg
Graduation	1 kg
Material	Stainless Steel / Hardened Steel
Supplied With	Penetrometer, Penetration Needle, 6 Nos. Needle Points, Plunger, Spring Housing & Wooden Carrying Case
Weight (Approx.)	2.5 - 3.0 kg

SUPPLIED WITH



Spring Housing
(With Plunger)



Penetration
Needle



Plunger
(Graduated)



Needle Points
(Set of 6 Nos.)



Wooden
Carrying Case



Precise
Measurement



High Accuracy
& Reliability



Portable
& Durable



Easy to Use
& Maintain

+91-8130973844

architassst@gmail.com



www.architassociates.com



ONE YEAR WARRANTY

Against manufacturing defects
(excluding consumables)



AMC & REPAIRING FACILITY AVAILABLE

Comprehensive Annual Maintenance Contract
and repairing services available

Conforms to
IS:1199 (Part-7)

ARCHIT ASSOCIATES - YOUR PARTNER IN TESTING EXCELLENCE



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir,
Delhi-110 033 (INDIA)

91-8130973844, 9871406036

architasst@gmail.com

www.architassociates.com

GSTIN NO. 07A0UPG0300F22C

TRUSTED SUPPLIER OF
**MATERIAL TESTING &
LABORATORY INSTRUMENTS**



QUALITY
ASSURED
PRODUCTS



ACCURATE &
RELIABLE



LABORATORY &
SITE SOLUTIONS



CUSTOMIZED
TESTING
EQUIPMENT

COOLING CURING TANK

TEMPERATURE FROM 10°C TO BOILING TEMPERATURE

SUITABLE FOR CUBE MOULDS

Available Capacity: 6 Nos, 12 Nos, 18 Nos, 24 Nos. etc.

PRODUCT DESCRIPTION

Cooling Curing Tank is specially designed for curing of concrete cubes, mortar cubes, cylinders and other concrete specimens as per IS standards.

It is fitted with a high precision digital temperature controller with an inbuilt cooling system and heating arrangement to maintain water temperature from 10°C to boiling temperature.

The tank is made of thick M.S. sheet with corrosion resistant powder coated finish for long lasting performance. The unit is highly reliable, accurate and easy to operate.



TECHNICAL SPECIFICATIONS

Parameter	Specification
Temperature Range	10°C to Boiling Temperature
Temperature Controller	Digital PID Controller
Temperature Accuracy	± 1°C
Cooling System	Compressor Based
Heating System	Immersion Heater
Tank Material	M.S. Sheet (Powder Coated)
Insulation	High Density PUF / Glass Wool
Drain Valve	Provided
Power Supply	230 V AC, 50 Hz
Capacity (Cube Moulds)	6 Nos, 12 Nos, 18 Nos, 24 Nos. etc.

SALIENT FEATURES

- Temperature range from 10°C to boiling temperature.
- Digital PID temperature controller for accurate control.
- Efficient cooling with powerful compressor.
- Uniform temperature distribution.
- Rust proof tank with powder coated finish.
- Easy operation and low maintenance.
- Suitable for IS standard cube moulds.

AVAILABLE MODELS (AS PER CAPACITY)



WIDELY USED IN



CEMENT
LABORATORIES



CIVIL ENGINEERING
LABORATORIES



QUALITY CONTROL
DEPARTMENTS



RESEARCH &
INSTITUTES



CONSTRUCTION
INDUSTRY



ARCHIT ASSOCIATES
A-242, Jahangir Puri,
Near Kali Mandir,
Delhi-110 033 (INDIA)

91-8130973844

architasst@gmail.com

www.architassociates.com

ALL PRODUCTS COMPLY WITH RELEVANT INDUSTRY STANDARDS



ARCHIT ASSOCIATES

YOUR TRUSTED PARTNER IN LABORATORY & MATERIAL TESTING

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

+91-8130973844 9871406036 architasst@gmail.com

www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



SCAN TO VIEW
PRODUCT RANGE

DEEP ABRASION TESTING MACHINE

REF. STD. IS 13630 (PART 12), ISO 10545 (PART 6)

Used to determine the resistance of unglazed tiles to deep abrasion. The machine presses the tile against a rotating steel disc while abrasive is dropped at a controlled rate.

SALIENT FEATURES

- Conforms to IS : 13630 (Part 12), ISO 10545 (Part 6)
- Adjustable Load : 5 to 20 kg (in steps of 1 kg)
- Rotating Steel Disc : 75 ± 5 rpm
- Preset Electronic Counter with Auto Motor Stop
- Abrasive Feed with Flow Control Valve from Storage Hopper
- Robust Construction with Corrosion Resistant Finish
- Collection Tray for Used Abrasive
- Supplied with Calibration Glass Plates (150 x 150 x 10 mm) - 2 Nos.

OPERATING PRINCIPLE

The ability of unglazed tiles to resist deep abrasion is determined by holding the tile under test in a vertical plane and pressing it against a steel disc rotating at a specified speed. The abrasive material of specified grit size is dropped from a funnel at fixed rate so as to fall between the face of the tile under test and the steel disc. The length of impression produced on the surface of the tile under test is measured. The volume of material abraded away from the tile is calculated to evaluate resistance to deep abrasion of the tile.

TECHNICAL SPECIFICATIONS

Diameter of Steel Disc	200 $\pm$ 0.2 mm
Thickness of Steel Disc	10 $\pm$ 0.1 mm
Material of Steel Disc	Mild Steel (Commercially Grade)
Speed of Rotation of Disc	75 $\pm$ 5 rpm
Force Applied on Test Specimen	Adjustable between 5 and 20 kg in steps of 1 kg
Motor	1/4 HP Single Phase 230 Volts AC
Counter	3 Digit Pre Set Type Electronic Counter with Memory Back Up and Push Button Reset
Capacity of Storage Hopper	5 Liter
Power Supply	230 Volts, 50 Hz, Single Phase, AC Supply
Weight	90 - 100 Kgs Approx.

CONSTRUCTION DETAILS

- Metal disc held with its axis in a horizontal plane.
- Motorized arrangement to rotate the disc
- Pre set type electronic counter to record the number of rotations and auto stop the motor
- Frame to hold the test specimen in a vertical plane tangential to the disc
- Trolley on which the test specimen holder is mounted
- Dead weights to pull the trolley and press the specimen against the disc

APPLICATIONS

- Tile Manufacturing Industry
- Quality Control Laboratories
- Research & Development Centres
- Government Testing Laboratories
- Building Material Testing
- Educational Institutions

DETAILS

- The metal disc is made of commercially available mild steel. It is rotated about a horizontal axis at the specified speed with the help of an electric motor and a two stage V-belt and pulley arrangement. A pre set type electronic counter with reset and memory backup is provided to record the total number of rotations and also to switch off the motor after the desired number of rotations.
- The test specimen is held in a clamp such that its upper face lies in a vertical plane and can be pressed against the outer face of the metal disc. The clamp is mounted on a wheeled trolley moving over guide rods. The trolley can be pulled with dead weights to press the test specimen against the metal disc.
- Abrasive material in the form of fused aluminum oxide powder of 80 grit is dropped between the test specimen and the outer face of metal disc from a funnel. A storage hopper over the funnel supplies abrasive continuously. The rate of flow can be controlled manually with the ball valve fitted below the hopper.
- All components are mounted on a rigid metallic frame with a collection tray for used abrasive. Transparent silica glass plates (150 x 150 x 10 mm) - 2 Nos. are supplied for calibration purpose. The machine is finished in grey hammerstone stoving painting and chrome/zinc plated parts for corrosion resistance.

STANDARD ACCESSORIES



OPTIONAL ACCESSORIES



ACCURATE
RESULTS



RUDD & DURABLE
CONSTRUCTION



ADJUSTABLE
LOAD



PRESET COUNTER
WITH AUTO STOP



CORROSION RESISTANT
FINISH



EASY TO OPERATE &
MAINTAIN



RELIABLE &
LONG LASTING

ARCHIT ASSOCIATES - YOUR PARTNER IN TESTING EXCELLENCE



ARCHIT ASSOCIATES

YOUR TRUSTED PARTNER IN LABORATORY & MATERIAL TESTING

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

91-8130973844 ✉ architassst@gmail.com

www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



ACCURATE
RESULTS



RELIABLE
PERFORMANCE



QUALITY
ASSURED

TESTING EXCELLENCE IN EVERY RESULT

DIGITAL WRAP REEL



3 Lea



6 Lea



AUTOMATIC WRAP REEL

- ◆ Digital Wrap Reel Manufactured by us Is Of Finest Quality And Covered By Warranty For Manufacturing Defects.
- ◆ Digital Wrap Reel Is A Precision Made, Robust And Reliable Reading Machine Which Is Used For Rapid Preparation Of Accurate Lengths Of Yarn. In Hank Form For Subsequent Weighing And Determination Of Count And Strength Of The Yarn .It Gives Accurate And Perfect Of Final Count / Denier Of The Yarn Sample / Specimen.



FEATURES

- ◆ The Equipment Consists Of Solidly Constructed Swift / Wheel.
- ◆ It Has Traverse Mechanism
- ◆ Wrap Reel Is Fitted With A Digital Counter.
- ◆ It Is Powered By A Special Geared Motor.



SPECIFICATION

- ◆ Swift Perimeter : 1 Meter Number Of Lead Made At Time: 3/6/10/12/16
- ◆ Spread Of Each Lea : 20 mm
- ◆ Net Weight Of The Main Unit : 20 Kg
- ◆ Controls : Main , Zero / Reset Controller.



HOW DIGITAL WRAP REEL WORKS?

- ◆ Bobbins / yarn Cones Are Placed On The Yarn Cone Stand.
- ◆ The end of the Yarn Has To Be Passed Through The Tensioning Device Into The Ceramic Eye Lids, Through To The Ceramic Pig Tail and Finally Onto The Plate With Screw Heads. The Free End Of The Yarn Has To Be Tied On The Screw Head By Making Few Rounds.
- ◆ The Mains Is Put On And The Counter Is Set At Desired Setting (90 For Denier) With The Help Of Thumb Wheel Setting.
- ◆ After The Required Revolutions , Wrap Reel Will Stop Automatically.
- ◆ The Free end of Each Of The Lead has to be United With The Other End Of The Same Lead. The Bobbin Side The Knot Has To Be Cut.



HIGH ACCURACY

Precision made for accurate length & count measurement.



ROBUST & DURABLE

Sturdy construction for long life and reliable performance.



DIGITAL COUNTER

Built-in digital counter for precise and easy reading.



POWERFUL MOTOR

Equipped with special geared motor for smooth operation.



EASY OPERATION

User friendly design for quick and efficient testing.



QUALITY ASSURED

Manufactured under strict quality control with warranty coverage.



PRECISE RESULTS

Consistent & accurate testing performance



SUPERIOR QUALITY

High grade components for long life



EASY TO OPERATE

Intuitive controls with advanced technology



TEST CERTIFICATES

Supplied with calibration certificates



AFTER SALES SUPPORT

Prompt service & technical support



SCAN TO VISIT
OUR WEBSITE

ARCHIT ASSOCIATES

A-242, Jahangir Puri, Delhi - 110033 (INDIA)

91-8130973844

✉ architassst@gmail.com

28

YEARS OF
EXPERIENCE



BUILT FOR LABS.
ENGINEERED FOR EXCELLENCE.
Precision. Performance. Perfection.



ARCHIT ASSOCIATES

YOUR TRUSTED PARTNER IN LABORATORY & MATERIAL TESTING

📍 A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

☎ +91-8130973844, 9871406036 ✉ architasst@gmail.com

🌐 www.architassociates.com

GSTIN NO. 07A0UPG0300FZZC



ACCURATE
RESULTS



RELIABLE
PERFORMANCE



QUALITY
ASSURED

TESTING EXCELLENCE IN EVERY RESULT

14" DISC POLISHING MACHINE

HEAVY DUTY • DUST COLLECTION • THREE PHASE

14" Disc Polishing Machine is designed for grinding and polishing of metallographic specimens to achieve mirror finish. The machine is fitted with heavy duty motors, dust collection system and user friendly control panel. Ideal for laboratories, research centres and industrial quality control departments.

FEATURES

- 14 Inch MS Polishing Disc**
Heavy duty MS disc for efficient polishing
- Dual Motor System**
2800 RPM Motor for polishing
1440 RPM Motor for disc rotation
- Dust Collection System**
Built-in outlet for effective dust extraction
- Three LED Indicators**
R-Y-B LED for power & operation status
- Start / Stop Push Button**
Easy & safe operation
- Three Phase MC**
Reliable magnetic contactor for safe control
- Heavy Duty MS Body**
Rugged construction for long life
- Low Maintenance**
Simple design for easy maintenance
- Stable & Vibration Free Operation**
Ensures uniform polishing & accuracy



COMPLETE WITH

- 14 Inch MS Polishing Disc – 1 No.
- Dust Collection Outlet – 1 No.
- Start / Stop Push Buttons – 1 Set
- Three LED Indicators (R-Y-B) – 1 Set
- Three Phase MC with Overload – 1 Set
- Instruction Manual – 1 No.



TECHNICAL SPECIFICATIONS

Disc Diameter	14 Inch (355 mm)
Disc Material	Mild Steel (MS)
Disc Type	Fixed Disc
Polishing Motor	2 HP, 2800 RPM, Three Phase
Disc Motor	1 HP, 1440 RPM, Three Phase
Power Supply	415V, Three Phase, 50Hz
Control	Start / Stop Push Button
LED Indicators	3 Nos. (Red, Yellow, Blue)
Dust Collection	Yes, with Outlet
Electrical Control	Three Phase MC with Overload
Body Material	Mild Steel (MS)
Body Finish	Powder Coated
Noise Level	Low Noise
Overall Dimensions (L x W x H)	750 x 600 x 950 mm (Approx.)
Net Weight (Approx.)	110 – 120 kg (Approx.)

IDEAL FOR

- METALLOGRAPHY LABS
- RESEARCH & DEVELOPMENT
- QUALITY CONTROL LABORATORIES
- STEEL PLANTS & FOUNDRIES
- EDUCATIONAL INSTITUTES

OPTIONAL ACCESSORIES

- Silicon Carbide Papers
- Alumina Powder
- Diamond Polishing Cloth
- Velvet Cloth
- Diamond Suspension



PRECISE RESULTS
Consistent & accurate
testing performance



SUPERIOR QUALITY
High grade components
for reliability



EASY TO OPERATE
User friendly design
& simple controls



TEST CERTIFICATES
Supplied with
calibration certificates



AFTER SALES SUPPORT
Prompt service &
technical support



SCAN TO VISIT
OUR WEBSITE

ARCHIT ASSOCIATES

📍 A-242, Jahangir Puri, Delhi - 110033 (INDIA)

☎ +91-8130973844, +91 9871406036

✉ architasst@gmail.com 🌐 <http://bit.ly/architassociates>

28
YEARS OF
EXPERIENCE



**BUILT FOR LABS,
ENGINEERED FOR EXCELLENCE**
Precision. Performance. Trust.



ARCHIT ASSOCIATES

YOUR TRUSTED PARTNER IN LABORATORY & MATERIAL TESTING

📍 A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

☎ +91-8130973844 9871406036

🌐 www.architassociates.com

✉ architasst@gmail.com

GSTIN NO. 07A0UPG0300F2ZC



SCAN TO VIEW
PRODUCT RANGE

DORRY ABRASION TESTING MACHINE

REF. STD. BS 812

Used for determining the resistance to abrasion of natural stone and other building materials.

SALIENT FEATURES

- ✓ Conforms to BS 812
- ⚙ Robust & Heavy Duty Construction
- 🔄 Uniform & Consistent Abrasion Test
- 🌀 Rotating Disc with Precise Speed Control
- ⌚ Dual Sand Hopper with Flow Control
- ⚖ Specimen Holder with Dead Weight
- 🔧 Easy Operation & Maintenance
- 🛡 Corrosion Resistant Finish
- 🕒 Long Life & Reliable Performance



STANDARD ACCESSORIES



ABRASION DISC
(1 No.)



SAND HOPPER
(2 Nos.)



SPECIMEN WEIGHT
CONTAINER (2 Nos.)

TECHNICAL SPECIFICATIONS

Standards	BS 812
Rotating Disc Diameter	710 ± 5 mm
Disc Rotating Speed	28 - 30 RPM
Motor	0.5 HP, Single Phase
Power Supply	220 V, 50 Hz, Single Phase, AC Supply
Sand Hopper	Two Conical Hoppers with Flow Control
Specimen Holder	Two with Dead Weight
Abrasive	Standard Silica Sand
Material of Disc	Cast Iron / Mild Steel
Overall Dimensions (L x W x H)	Approx. 1100 x 700 x 1050 mm
Net Weight (Approx.)	75 - 85 Kgs

WORKING PRINCIPLE

A specimen is pressed against the rotating abrasion disc by means of a loaded container. Abrasive sand is fed uniformly on the disc through two conical hoppers. The disc rotates at 28-30 RPM and the abrasive action wears the surface of the specimen. After the specified number of cycles, the specimen is cleaned and the loss in thickness/weight is measured to evaluate the abrasion resistance.

APPLICATIONS



Building
Materials



Natural Stone
Industry



Quality Control
Laboratories



Research &
Development



Testing
Laboratories

OPTIONAL ACCESSORIES



Thickness Gauge
(Dial / Digital)



Weighing Balance
(Top Pan)



Silica Sand
(As per Std.)



Spare Abrasion Disc



Calibration Certificate



ACCURATE
TEST RESULTS



RUGGED & DURABLE
CONSTRUCTION



EASY TO OPERATE
& MAINTAIN



PRECISE DISC
ROTATION



QUALITY ASSURED
EQUIPMENT



RELIABLE & LONG
SERVICE LIFE



SCAN TO VISIT
OUR WEBSITE

ARCHIT ASSOCIATES

📍 A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

☎ +91-8130973844 | architasst@gmail.com | www.architassociates.com

28+
YEARS OF
EXPERIENCE



**BUILT FOR LABS.
ENGINEERED FOR EXCELLENCE.**
Precision. Performance. Trust.



ARCHIT ASSOCIATES

YOUR TRUSTED PARTNER IN LABORATORY & MATERIAL TESTING

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

+91-8130973844, 9871406036 | architasst@gmail.com

www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



ACCURATE
RESULTS



RELIABLE
PERFORMANCE



QUALITY
ASSURED

TESTING EXCELLENCE IN EVERY RESULT

DOUBLE DISC POLISH MACHINE

TWO DISC • VARIABLE SPEED • WATER INLET / OUTLET

Double Disc Polish Machine is designed for polishing metallographic specimens to obtain a scratch-free mirror finish. The machine is widely used in laboratories, research centres and quality control departments.

FEATURES

- Aluminium Polishing Disc**
Durable & corrosion resistant 8 inch disc
- Dual Drive System**
Independent drive for both discs
- RPM Display (2 Nos.)**
Digital RPM display for both discs
- Water Inlet / Outlet**
Easy water connection for cooling & cleaning during polishing
- Start / Stop Push Button**
User friendly control operation
- 0.5 HP Motor**
High performance & energy efficient
- Single Phase MC**
Reliable & safe electrical control panel
- Robust Construction**
Heavy duty body for long life & maintenance free operation



COMPLETE WITH

- Aluminium Disc 8 inch (2 Nos.)
- Water Inlet / Outlet Pipes
- Power Cable & Plug
- Instruction Manual



TECHNICAL SPECIFICATIONS

Type	Double Disc Polish Machine
Disc Material	Aluminium
Disc Diameter	8 inch (200 mm)
Number of Discs	2
Drive	Dual Drive (Individual)
Motor	0.5 HP
Speed (RPM)	0 - 1500 RPM (Variable)
RPM Display	2 Nos. (Digital)
Water Connection	Inlet / Outlet
Control	Start / Stop Push Button
Electrical Supply	Single Phase, 230V, 50Hz
Control Panel	Single Phase MC
Construction	Heavy Duty Mild Steel Body
Overall Dimensions (L x W x H)	700 x 450 x 450 mm (Approx.)
Net Weight (Approx.)	70 kg (Approx.)

IDEAL FOR

METALLOGRAPHY
LABS

RESEARCH &
DEVELOPMENT

QUALITY CONTROL
LABORATORIES

EDUCATIONAL
INSTITUTES

OPTIONAL ACCESSORIES

- Polishing Cloths
- Diamond Suspensions
- Alumina Suspensions
- Specimen Holder



PRECISE RESULTS
Consistent & accurate
testing performance



SUPERIOR QUALITY
High grade components
for reliability



EASY TO OPERATE
User friendly design
& simple controls



TEST CERTIFICATES
Supplied with
calibration certificates



AFTER SALES SUPPORT
Prompt service &
technical support



SCAN TO VISIT
OUR WEBSITE

ARCHIT ASSOCIATES

A-242, Jahangir Puri, Delhi - 110033 (INDIA)

+91-8130973844, +91 9871406036

architasst@gmail.com | <https://bit.ly/google.com/view/architassociates>

28
YEARS OF
EXPERIENCE



BUILT FOR LABS.
ENGINEERED FOR EXCELLENCE.
Precision. Performance. Trust.



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir,
Delhi-110 033 (INDIA)

91-8130973844 architasst@gmail.com

www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC

TRUSTED SUPPLIER OF
MATERIAL TESTING &
LABORATORY INSTRUMENTS



QUALITY
ASSURED
PRODUCTS



ACCURATE &
RELIABLE



LABORATORY &
SITE SOLUTIONS



CUSTOMIZED
TESTING
EQUIPMENT

FLOW TABLE MOTORISED

Flow Table is used for determining the flow and workability of cement mortar and other related materials. It is widely used in cement laboratories, civil engineering laboratories and quality control departments.

FEATURES



Motor driven for smooth and consistent operation



High accuracy and reliability



Sturdy and compact design



Low maintenance and easy to operate



Better aesthetics and compact look

FLOW TABLE – MOTORISED



TECHNICAL SPECIFICATIONS

Parameter	Details
Table Top Material	Brass
Table Top Diameter	250 ± 2.5 mm
Drop Height	12 mm
Cam Mechanism	Ground & Hardened
Motor Speed (for Motorised)	100 RPM
Flow Mould Base Diameter	-
Flow Mould Top Diameter	-
Flow Mould Height	-
Construction	Heavy-duty
Finish	Powder Coated
Models Available	Motorised

WIDELY USED IN



CEMENT
LABORATORIES



CIVIL ENGINEERING
LABORATORIES



QUALITY CONTROL
DEPARTMENTS



RESEARCH &
INSTITUTES



CONSTRUCTION
INDUSTRY

MODELS AVAILABLE



FLOW TABLE
HAND OPERATED



FLOW TABLE
MOTORISED



QUALITY ASSURED
PRODUCTS



ACCURATE &
RELIABLE



CUSTOMIZED
TESTING EQUIPMENT



LABORATORY &
SITE SOLUTIONS



ARCHIT ASSOCIATES
A-242, Jahangir Puri,
Near Kali Mandir,
Delhi-110 033 (INDIA)



91-8130973844



architasst@gmail.com

ALL PRODUCTS COMPLY WITH RELEVANT INDUSTRY STANDARDS



ARCHIT ASSOCIATES

📍 A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

☎ Phone: 91-8130973844

✉ E-mail: architasst@gmail.com

🌐 Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



FLOW TABLE FOR SELF COMPACTING CONCRETE

Ref. Standards - IS: 9103:1999

Flow table is used to determine the flow of fresh mixed super plasticized concrete to high working.

The slump cone is placed centrally on the table to be held in position by standing on the two foot pieces. A wooden tamping bar is provided for lightly tamping for each layer.



STANDARD

IS: 9103 : 1999

Specification for flow table for self compacting concrete.

APPLICATION

Used in laboratories and at site to determine the flow/slump flow of self compacting concrete (SCC) to assess its workability and filling ability.

KEY FEATURES

- ✓ Robust and rigid construction for long life
- ✓ Uniform table surface for accurate results
- ✓ Easy to operate and maintain
- ✓ Provided with slump cone and tamping bar
- ✓ Corrosion resistant and powder coated finish
- ✓ Conforms to IS: 9103:1999


TECHNICAL SPECIFICATIONS

Table Top	900 mm x 900 mm
Table Top Thickness	12 mm
Table Height	200 mm
Drop Height	12 ± 2.0 mm
Slump Cone Top Diameter	100 ± 2 mm
Slump Cone Bottom Diameter	200 ± 2 mm
Slump Cone Height	300 ± 2 mm
Tamping Bar	40 mm dia x 300 mm long
Finish	Powder Coated

PACKAGE / SYSTEM INCLUDES

- | | |
|------------------------|----------------------|
| ✓ Flow Table | ✓ Two Foot Pieces |
| ✓ Slump Cone (Brass) | ✓ Handle |
| ✓ Tamping Bar (Wooden) | ✓ Instruction Manual |

DIMENSIONS (Approx.)

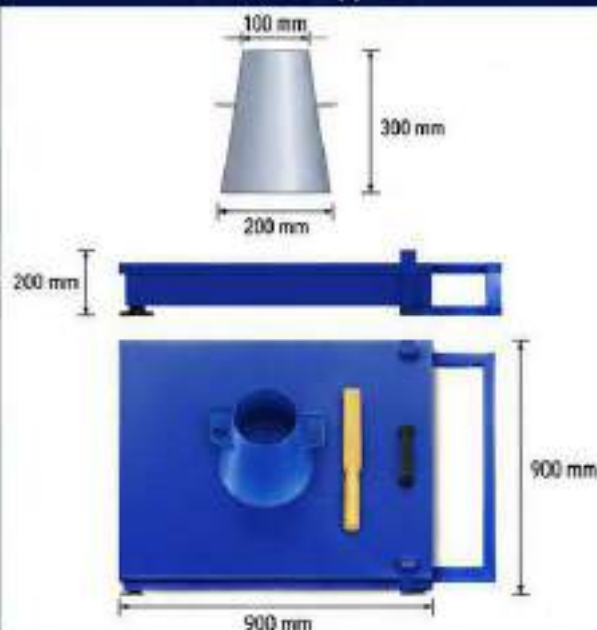

ARCHIT ASSOCIATES

📍 A-242, Jahangir Puri,
Near Kali Mandir, Delhi-110 033 (INDIA)

☎ 91-8130973844

✉ architasst@gmail.com

ARCHIT ASSOCIATES – YOUR PARTNER IN TESTING EXCELLENCE



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kaili Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844

E-mail: architassst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07AOUPG0300F2ZC



ACCURATE
RESULTS



RELIABLE
PERFORMANCE



QUALITY
ASSURED



FREEZE THAW CHAMBER

TEMPERATURE RANGE: -20°C TO 50°C

Storage conditions must be distributed in a manner that ensures that the product quality will not be adversely affected. With the exception of short transit times within the same climatic zone, it is virtually impossible to validate a shipping method against all environmental conditions.

Therefore, distribution stability studies are designed to generate additional data set to complement preformulation development and routine (ICH Q1A/R2) studies.



SPECIFICATIONS

- Temperature Range : -20°C to 50°C

SELECTION GUIDE

Capacity		Internal Size (W×H×D) mm	External Size (W×H×D) mm	Trays Qty	Trays Size (W×D) mm	Approx. Weight (Kg)	
Cu.ft.	Liter					Emp.	Ship.
7	200	510 x 580 x 700	1500 x 720 x 880	2	540 x 605	250	300
11	300	910 x 580 x 580	1855 x 720 x 770	3	540 x 495	270	425
14	400	910 x 580 x 785	1855 x 720 x 965	3	540 x 800	280	445
25	700	1200 x 800 x 730	2120 x 800 x 1020	4	780 x 635	310	500
34	1050	1200 x 1200 x 730	2110 x 1300 x 1020	6	515 x 835	380	600
54	1500	1500 x 1200 x 830	2310 x 1500 x 1120	8	515 x 735	490	700

Note: Custom sizes and tray configurations available on request.



PRECISE
TEMPERATURE
CONTROL



HUMIDITY
CONTROL
50% to 95%



CO₂ CONTROL
1% to 10%



UNIFORM AIR
CIRCULATION



VERTICAL ROTATION
15 ±0.2 RPM



AIR FLOW
4000 mL/min



SAFE & RELIABLE
OPERATION



RUGGED & DURABLE
CONSTRUCTION



ONE YEAR WARRANTY
Against manufacturing defects
(excluding consumables)



AMC & REPAIRING FACILITY AVAILABLE
Comprehensive Annual Maintenance Contract
and repairing services available

ARCHIT ASSOCIATES

A-242, Jahangir Puri,
Near Kaili Mandir, Delhi-110 033 (INDIA)
91-8130973844

architassst@gmail.com

GSTIN NO. 07AOUPG0300F2ZC

ARCHIT ASSOCIATES – YOUR PARTNER IN TESTING EXCELLENCE

STANDARD EQUIPMENT

- Double walled with PUF insulation
- SS304 Interior & Exterior of Mild Steel
- Stainless Steel trays
- Metal door with gasket & lock on the outside
- Forced air circulation for uniform temperature
- PU wheels for easy movement
- U shaped SS tubular heaters
- CFC free cooling system consisting of hermetically sealed compressor coupled with evaporation coil & condenser
- Complete Operations driven through PLC with data logging capacity of upto 1000 readings
- High speed Ethernet based communication
- All operations through a single 7" Color Touch Screen
- High-quality PT-100 Temperature sensor
- Safety devices: High Temperature cut-off settable thermostats
- Audio/Visual alarms for deviations
- Complete with DO, IQ, OQ, PQ & SOP & Operational Manual

Power Supply: 230 volts 50 Hz main supply. Also available with 60 Hz Cycle.

MODE OF OPERATION

- User can program a specific profile needed for their application to enable rapid simulation conditions & the need to replicate freezing & thawing during product testing
- Upto 4 temperatures set point can be set in cyclic mode.
(eg. -20°C set for 48 hrs, 5°C for 24 hrs, 25°C for 72 hrs and 40°C for 120 hrs) No. of such temperature cycles can be repeated

OPTIONAL ACCESSORIES

- ICDAS - A Single Software Installation for all Newtronic make equipments
- Door Access system (Numeric or Biometric)
- Remote alarm system suitable for upto 4 nos. Newtronic make equipments
- Mobile Alert System



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844 | E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



LABORATORY SPARES & ESSENTIAL TOOLS

QUALITY YOU CAN TRUST, PERFORMANCE YOU CAN RELY ON

A wide range of laboratory spares, tools and accessories for civil engineering, material testing and general laboratory applications. Built for accuracy, durability and long service life.



ACCURATE
RESULTS



DURABLE &
RELIABLE



WIDE RANGE OF
APPLICATIONS



QUALITY
ASSURED



LAB GLASSWARE



Beaker



Conical Flask



Measuring Cylinder



Volumetric Flask



Pipette



Burette



Funnel



Watch Glass



Test Tube

LABORATORY SPARES & ACCESSORIES



Enamel Tray



Standard Sand Bag
25 kg



Trowel



Plumb Bob



Masonry Hammer



Scoop



Brush

MEASURING & TESTING TOOLS



Vernier Caliper



Digital Vernier Caliper



Steel Scale



Weighing Scale



Thermometer



Dial Gauge



Stop Watch

POWER TOOLS & EQUIPMENT



Drill Machine



Concrete Mixer



Vibrating Table



Hand Rammer



Sieves
(Various Sizes)



Wheel Barrow

OTHER ESSENTIALS

- Spatula
- Glass Rod
- Mortar & Pestle
- Rubber Mallet
- Sample Containers
- Filter Paper



IDEAL FOR

- Civil Engineering Laboratories
- Construction Sites
- Quality Control Labs
- Research & Development
- Educational Institutions

QUALITY YOU CAN RELY ON

All products are carefully selected for their quality, precision and durability to ensure consistent and accurate performance in every test.



ARCHIT ASSOCIATES

A-242, Jahangir Puri,
Near Kali Mandir, Delhi-110 033 (INDIA)
91-8130973844

architasst@gmail.com

GSTIN NO. 07A0UPG0300F2ZC

YOUR PARTNER IN
LABORATORY EXCELLENCE



ARCHIT ASSOCIATES

Trusted Supplier of
MATERIAL TESTING & LABORATORY INSTRUMENTS

A-242, Jahangir Puri, Near Kali Mandir,
Delhi-110 033 (INDIA)

architasst@gmail.com

91-8130973844, 9871406036,
9810052732, 9810052758

www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC

HOT PLATE

Reliable Heating. Precise Results.

FITTED WITH ENERGY REGULATOR &
DIGITAL TEMPERATURE CONTROLLER



QUALITY
ASSURED



ACCURATE &
RELIABLE



LABORATORY
SOLUTIONS



CUSTOMIZED
EQUIPMENT

PRODUCT RANGE



RECTANGULAR MODEL



ROUND MODEL

KEY FEATURES



Uniform
Heating



Energy
Regulator



Digital Temperature
Controller (Optional)



Safe &
Reliable



Easy Temperature
Control



Low
Maintenance



Energy
Efficient



Special Sizes
Available for
Customization

STANDARD MODELS – FITTED WITH ENERGY REGULATOR & ON/OFF SWITCH

S. No.	Size (Top Plate)	Top Plate Type	Dimensions (L x W x H)	Power Supply
1	200 mm dia (Round)	MS	Dia 200 x 180 mm	230V AC, 50 Hz
2	250 x 300 mm	MS	250 x 300 x 180 mm	230V AC, 50 Hz
3	250 x 400 mm	MS	250 x 400 x 180 mm	230V AC, 50 Hz
4	300 x 450 mm	MS	300 x 450 x 180 mm	230V AC, 50 Hz
5	450 x 600 mm	MS	450 x 600 x 180 mm	230V AC, 50 Hz

STANDARD MODELS – FITTED WITH DIGITAL TEMPERATURE CONTROLLER-CUM-INDICATOR

S. No.	Size (Top Plate)	Top Plate Type	Dimensions (L x W x H)	Power Supply
1	200 mm dia (Round)	MS	Dia 200 x 180 mm	230V AC, 50 Hz
2	250 x 300 mm	MS	250 x 300 x 180 mm	230V AC, 50 Hz
3	250 x 400 mm	MS	250 x 400 x 180 mm	230V AC, 50 Hz
4	300 x 450 mm	MS	300 x 450 x 180 mm	230V AC, 50 Hz
5	450 x 600 mm	MS	450 x 600 x 180 mm	230V AC, 50 Hz



**SPECIAL SIZES
FOR CUSTOMIZATION**

Special sizes of hot plates can be
manufactured as per customer
requirements.
Please contact us for customization.



Customized
Solutions



High Quality
Construction



Consistent
Performance



Ideal for
Laboratories

A-242, Jahangir Puri,
Near Kali Mandir,
Delhi-110 033 (INDIA)



91-8130973844
9871406036
9810052732
9810052758



architasst@gmail.com

ALL PRODUCTS COMPLY WITH RELEVANT INDUSTRY STANDARDS



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844 E-mail: architasst@gmail.com

Web: www.architassociates.com



HUMIDITY & TEMPERATURE CABINETS

PRECISE CONTROL. RELIABLE PERFORMANCE.

HUMIDITY & TEMPERATURE CABINET



Ideal for testing the effect of humidity and temperature on materials, components and products in research, industrial and quality control laboratories.

- Temperature Range
Ambient to 60°C $\pm$ 1°C
- Humidity Range
Ambient to 95% $\pm$ 3%
- Inner Chamber
SS 304 Grade
- Digital Controller
Cum Indicator
- Steam Based
Humidity Generation
- Inside Acrylic Door &
SS Water Tank Provided
- Safe, Reliable &
Easy to Operate

HUMIDITY CABINET COOLING / ENVIRONMENTAL TEST CHAMBER (ECO FRIENDLY)



Suitable for environmental simulation, stability studies, material testing and accelerated life testing applications.

- Cooling Temp. Range
10°C to 60°C $\pm$ 1°C
- Humidity Range
50% to 95% $\pm$ 4%
- Eco Friendly
Cooling System
- Inner Chamber
SS 304 Grade
- Digital Controller
Cum Indicator
- High Performance
& Durable Design

TECHNICAL SPECIFICATIONS (HUMIDITY & TEMPERATURE CABINET)

Temperature Range	Ambient to 60°C $\pm$ 1°C
Humidity Range	Ambient to 95% $\pm$ 3%
Temperature Controller	Digital Controller Cum Indicator
Humidity Controller	Digital Controller Cum Indicator
Display	Digital Display for Temperature & Humidity (% RH)
Chamber Material	Stainless Steel 304 Grade
Outer Body	Double Walled Mild Steel Duty Painted
Insulation	High Density PUF Insulation
Humidification	By Steam Generator
Door	Inside Acrylic Door
Water Tank	Stainless Steel
Power Supply	230 Volts AC, 50 Hz

TECHNICAL SPECIFICATIONS (COOLING / ENVIRONMENTAL TEST CHAMBER)

Cooling Temperature Range	10°C to 60°C $\pm$ 1°C
Humidity Range	50% to 95% $\pm$ 4%
Temperature Controller	Digital Controller Cum Indicator
Humidity Controller	Digital Controller Cum Indicator
Display	Digital Display for Temperature & Humidity (% RH)
Chamber Material	Stainless Steel 304 Grade
Outer Body	Double Walled Mild Steel Duty Painted
Insulation	High Density PUF Insulation
Refrigeration System	Eco Friendly (CFC / HCFC Free)
Compressor	Hermetically Sealed
Power Supply	230 Volts AC, 50 Hz

AVAILABLE SIZES (STANDARD)

Model (Approx.)	Inner Chamber Size (W x D x H)
HTC-150	450 x 450 x 750 mm
HTC-250	600 x 500 x 850 mm
HTC-400	600 x 600 x 1000 mm
HTC-600	750 x 750 x 1100 mm
HTC-1000	1000 x 800 x 1250 mm

* Other customized sizes available on request.

APPLICATIONS

- Pharmaceutical Testing
- Electronics & Components
- Automotive Testing
- Material Testing
- Research & Development
- Quality Control Laboratories

KEY FEATURES

- Precise Temperature & Humidity Control
- Uniform Air Circulation for Consistent Results
- High Quality SS 304 Inner Chamber
- Digital Control with High Accuracy
- Energy Efficient & Eco Friendly Operation
- Safe, Reliable & Easy Maintenance

OPTIONAL ACCESSORIES



DATA LOGGER



PRINTER



RS-485 COMMUNICATION



ADDITIONAL SHELVES



CASTOR WHEELS



WATER PURIFICATION UNIT

ARCHIT ASSOCIATES

A-242, Jahangir Puri,
Near Kali Mandir, Delhi-110 033 (INDIA)

91-8130973844

architasst@gmail.com

GSTIN NO. 07A0UPG0300F2ZC

ARCHIT ASSOCIATES - YOUR PARTNER IN TESTING EXCELLENCE



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir,
Delhi-110 033 (INDIA)
91-8130973844 architasst@gmail.com
www.architassociates.com
GSTIN NO. 07A0UPG0300F2ZC

TRUSTED SUPPLIER OF
**MATERIAL TESTING &
LABORATORY INSTRUMENTS**



QUALITY
ASSURED
PRODUCTS



ACCURATE &
RELIABLE



LABORATORY &
SITE SOLUTIONS



CUSTOMIZED
TESTING
EQUIPMENT

IN-SITU VANE SHEAR TEST APPARATUS

FOR DETERMINATION OF IN-PLACE SHEARING
RESISTANCE OF SOFT COHESIVE SOILS

CAT. NO. AA-138
(HAND OPERATED / STANDARD MODEL)

CAT. NO. AA-138-1
(MOTORISED VERSION)

Designed in accordance with IS:4434 for conducting in-situ Vane Shear test from the bottom of bore hole in saturated cohesive deposits to determine the undrained shear strength and sensitivity of soil.



ACCURATE
RESULTS



ROBUST &
DURABLE



EASY TO
OPERATE



IDEAL FOR
FIELD USE



CONFORMS TO
IS:4434

DESCRIPTION

The In-Situ Vane Shear Test Apparatus is designed in accordance with IS:4434 for conducting vane shear tests in saturated cohesive soils directly from bore holes. The equipment is used to determine the undrained shear strength and sensitivity of soft clay deposits.

- Torque applicator assembly mounted on a rigid base with graduated gear wheel and crank mechanism.
- Torque ring deforms as torque is applied and deformation is measured through a precision dial gauge.
- Pointer provided for registering rotation of vane.
- Detachable stand for anchoring and attachment to securely hold the rod string.
- Calibration curve supplied to convert dial gauge readings to torque in kg-cm.

TECHNICAL SPECIFICATIONS

Parameter	Details
Standard	IS:4434
Torque Applicator Assembly (Capacity)	2000 kg-cm
Torque Measurement	Torque Ring with Dial Gauge
Dial Gauge	0.002 mm ÷ 5 mm
Torque Rod	Square Section, 600 mm Long
Rods (Quick Coupling Type)	1 m Long – 25 Nos.
Additional Rods	0.5 m Long – 10 Nos.
Dummy Rod	1 No.
Vanes Supplied	37.5 mm dia × 75 mm high 50 mm dia × 100 mm high
Calibration	Calibration Curve Supplied (kg-cm of torque)
Stand	Detachable Type
Rod Holding Attachment	Provided
Weight (Approx.)	15 – 18 kg (without rods & vanes)

SALIENT FEATURES

- Designed specifically for in-situ vane shear testing in bore holes.
- Torque ring with section cut for accurate measurement.
- Dial gauge to measure torque ring deformation.
- Pointer for registering rotation of vane.
- Gear wheel marked in degrees.
- Smooth operation through crank mechanism.
- Detachable stand for easy anchoring and handling.
- Attachment provided to hold rod string securely.
- Calibration curve supplied to convert dial gauge readings to kg-cm of torque.
- Sturdy construction for field use in tough conditions.

STANDARD SUPPLY INCLUDES

- Torque Applicator Assembly, Capacity 2000 kg-cm
- Vane (with Vane Rod), 37.5 mm dia × 75 mm high
- Vane (with Vane Rod), 50 mm dia × 100 mm high
- Torque Rod (Square Cross Section), 600 mm Long
- Rods (Quick Coupling Type), 1 m Long
- Rods (Quick Coupling Type), 0.5 m Long
- Dummy Rod
- Dial Gauge 0.002 mm ÷ 5 mm
- Calibration Curve (kg-cm of torque)

* Other ranges available on request.

- 1 No.
- 1 No.
- 1 No.
- 1 No.
- 25 Nos.
- 10 Nos.
- 1 No.
- 1 No.
- 1 Set

OPTIONAL ACCESSORIES

- Vane (with Vane Rod), 65 mm dia × 130 mm high
- Vane (with Vane Rod), 75 mm dia × 150 mm high
- Vane (with Vane Rod), 100 mm dia × 200 mm high
- Guide with Ball Bearing Arrangement for use with 100 mm Casing Pipe
- Guide with Ball Bearing Arrangement for use with 150 mm Casing Pipe



MODELS AVAILABLE

Model No.	Description
AA-138	In-Situ Vane Shear Test Apparatus – Standard Hand Operated
AA-138-1	In-Situ Vane Shear Test Apparatus – Motorised Version

Standard: IS:4434 Compliant

ORDERING INFORMATION

Cat. No.	Description
AA-138	In-Situ Vane Shear Test Apparatus Complete with Standard Accessories
AA-138-1	In-Situ Vane Shear Test Apparatus, Motorised Version Complete with Standard Accessories

APPLICATIONS

- Foundation Engineering
- Geotechnical Investigations
- Soil Exploration in Bore Holes
- Design of Piles, Caissons & Foundations
- Research & Development
- Quality Control & Field Testing

STANDARDS

Conforms to
IS:4434
relevant IS & ISO standards
for in-situ vane shear testing



IDEAL FOR

Geotechnical
Consultants

Soil Investigation
Agencies

Construction
Projects

**BUILT FOR ACCURACY.
DESIGNED FOR FIELD PERFORMANCE.**

- Reliable Results
- Durable & Rugged Construction
- Easy to Operate
- Field Proven Design



ARCHIT ASSOCIATES

A-242, Jahangir Puri,
Near Kali Mandir, Delhi-110 033 (INDIA)
91-8130973844

architasst@gmail.com

GSTIN NO. 07A0UPG0300F2ZC



SCAN TO VIEW
PRODUCT RANGE

CAT. NO. AA-138
HAND OPERATED

CAT. NO. AA-138-1
MOTORISED VERSION



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844 E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



LABORATORY CRUSHER & PULVERISER

FOR CEMENT & MINERAL TESTING

Engineered for sample preparation in cement industries and laboratories. Ensures efficient size reduction with uniform output for accurate test results.



RUGGED
CONSTRUCTION



HIGH
PERFORMANCE



ACCURATE &
REPRODUCIBLE
RESULTS



EASY TO
OPERATE



IDEAL FOR
LABORATORIES

LABORATORY JAW CRUSHER



LABORATORY PULVERISER



LABORATORY JAW CRUSHER

- Designed for primary crushing of hard, medium hard, brittle & tough materials.
- Suitable for cement, clinker, limestone, slag, coal, ores & other minerals.
- Provides uniform output size for reliable laboratory analysis.



TECHNICAL SPECIFICATIONS

Feed Size	Up to 50 mm
Output Size	2 - 10 mm (Adjustable)
Crushing Capacity	100 - 150 kg/hr (Approx.)
Jaw Size	100 mm x 60 mm / 150 mm x 100 mm
Power	1 HP / 1.5 HP, 1440 RPM, 230V, 50Hz
Construction	Mild Steel Fabrication
Features	Replaceable Manganese Steel Jaws, Adjustable Discharge, Easy to Clean
Applications	Cement Plants, R&D Labs, Quality Control Labs

LABORATORY PULVERISER

- Designed for fine grinding of hard, medium hard & brittle materials.
- Ideal for cement, clinker, gypsum, limestone, slag & other minerals.
- Produces fine powder for chemical & physical testing.



TECHNICAL SPECIFICATIONS

Feed Size	Up to 10 mm
Output Size	90 - 200 Micron (Adjustable)
Grinding Capacity	20 - 30 kg/hr (Approx.)
Grinding Chamber	Hardened Steel / Carbide Lined
Power	1.5 HP / 2 HP, 1440 RPM, 230V, 50Hz
Construction	Mild Steel Fabrication
Features	High Speed Rotor, Replaceable Liners, Easy to Clean Collection Drawer
Applications	Cement Plants, R&D Labs, Quality Control Labs

IDEAL FOR SAMPLE PREPARATION IN



CEMENT PLANTS



QUALITY CONTROL
LABORATORIES



R & D
LABORATORIES



MINERAL TESTING
LABORATORIES

SALIENT FEATURES

- ✓ Compact & Space Saving Design
- ✓ High Crushing & Grinding Efficiency
- ✓ Low Maintenance & Easy Operation
- ✓ Long Life & Reliable Performance
- ✓ Suitable for Continuous Laboratory Use

ARCHIT ASSOCIATES

A-242, Jahangir Puri,
Near Kali Mandir, Delhi-110 033 (INDIA)

91-8130973844

✉ architasst@gmail.com

GSTIN NO. 07A0UPG0300F2ZC



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844 E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



JOLTING APPARATUS

REF. STANDARD - IS:10078, EN-196-1, 13454-2

For Cement & Concrete Testing

The Jolting Apparatus is used to compact cement mortar and concrete in moulds by jolting action. It is designed as per IS:10078, EN-196-1 and 13454-2 standards.

The apparatus ensures uniform compaction and consistency in test specimens.



FEATURES



Robust Construction
Heavy duty construction for long life and reliable performance.



Standard Compliance
Designed as per IS:10078, EN-196-1 and 13454-2 standards.



Digital Counter
Provided with sensor-based digital counter with memory retention facility.



Locating Pins
Locating pins provided for easy and accurate mounting of moulds.



Mould & Hopper
Mould with hopper can be clamped rigidly to the table.

TECHNICAL SPECIFICATIONS

Table	Rectangular table, rigidly connected by two support arms to spindle.
Horizontal Distance	800 mm from the centre of the table.
Jolting Motion	Table can be raised and allowed to fall freely 15 mm by a cam.
Cam Speed	Cam rotates at a rate of 60 rpm.
Counter	Digital counter with sensor to stop the machine after 60 set jolts.
Memory Retention	Digital counter has memory retention facility.
Locating Pins	Provided for mounting the mould compartments on the table.
Mould Clamping	Mould, surrounded by the hopper, can be clamped rigidly to the table.
Power Supply	220V, 50 Hz, Single Phase, AC Supply.

STANDARD ACCESSORIES

MOULD WITH HOPPER



LOCATING PINS



DIGITAL COUNTER



KEY BENEFITS



Ensures uniform compaction of test specimens



Improves accuracy & repeatability of test results



Easy to operate & maintain



Rugged and durable design for laboratory use

APPLICATIONS



CEMENT TESTING



CONCRETE TESTING



AGGREGATE TESTING



SOIL TESTING



RESEARCH & QUALITY CONTROL LABS



PRECISE & RELIABLE



ACCURATE RESULTS



EASY TO USE



STANDARD COMPLIANT



DURABLE CONSTRUCTION



AFTER SALES SUPPORT

ARCHIT ASSOCIATES

A-242, Jahangir Puri,
Near Kali Mandir, Delhi-110 033 (INDIA)
91-8130973844

architasst@gmail.com

GSTIN NO. 07A0UPG0300F2ZC

ARCHIT ASSOCIATES - YOUR PARTNER IN TESTING EXCELLENCE



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)
Phone: 91-8130973844 E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07AOUPG0300F2ZC



LE - CHATELIER WATER BATH DIGITAL

FOR CEMENT TESTING

Designed for carrying out Le-Chatelier soundness test on cement in accordance with IS: 4031 (Part 3). The water bath maintains uniform temperature from 27°C to boiling in 27 ± 3 minutes. Inside and outside chamber made of Stainless Steel 304 grade.

FEATURES



WIDE TEMPERATURE RANGE

Temperature range from 27°C to boiling in 27 ± 3 minutes.



PREMIUM STAINLESS STEEL

Inside and outside chamber made of Stainless Steel 304 grade.



HIGH CAPACITY

Supplied with stainless steel rack to accommodate 12 Nos of Le-Chatelier moulds.



HEAVY DUTY STIRRER

Ensures uniform water circulation and temperature consistency.



POWERFUL HEATING

2 Nos of Stainless Steel Heating Elements 2KW each (Total 4KW).



DIGITAL CONTROL

Digital Programmable PID Temperature Controller-cum-indicator.



EASY OPERATION & SAFE

Main On/Off Switch, Heater On/Off Switch for safe and convenient operation.



TECHNICAL SPECIFICATIONS

PARAMETER	SPECIFICATION
Temperature Range	27°C to boiling in 27 ± 3 minutes
Temperature Controller	Digital Programmable PID Temperature Controller-cum-indicator
Temperature Accuracy	$\pm 1^\circ\text{C}$
Chamber Material	Stainless Steel 304 Grade (Inside & Outside)
Capacity	12 Nos of Le-Chatelier Moulds (with SS Rack)
Stirrer	Heavy Duty Stirrer for uniform water circulation
Heaters	2 Nos of Stainless Steel Heating Elements 2KW each (Total 4KW)
Controls	Main On/Off Switch, Heater On/Off Switch
Power Supply	220 Volt, Single Phase, 50 Hz
Applications	Le-Chatelier Soundness Test on Cement (IS: 4031 Part 3)

STANDARD ACCESSORIES

- Stainless Steel Rack (For 12 Nos of Le-Chatelier Moulds)
- Digital Programmable PID Controller
- Heavy Duty Stirrer
- 2 Nos S.S. Heating Elements (2KW each)
- Instruction Manual



OPTIONAL ACCESSORIES

- Le-Chatelier Mould
- Thermometer
- Insulated Lid
- Water Level Indicator



DIMENSIONS (APPROX.)

Parameter	Size (mm)
Length	600
Width	400
Height	550
Weight (Approx.) : 38 - 45 Kg	



ACCURATE
& RELIABLE



DURABLE
CONSTRUCTION



UNIFORM
TEMPERATURE



EASY TO OPERATE
& MAINTAIN



IDEAL FOR CEMENT
TESTING LABS

ARCHIT ASSOCIATES

A-242, Jahangir Puri,
Near Kali Mandir, Delhi-110 033 (INDIA)
91-8130973844

✉ architasst@gmail.com

GSTIN NO. 07AOUPG0300F2ZC

ARCHIT ASSOCIATES - YOUR PARTNER IN TESTING EXCELLENCE



ARCHIT ASSOCIATES

Trusted Supplier of
MATERIAL TESTING & LABORATORY INSTRUMENTS

📍 A-242, Jahangir Puri, Delhi - 110033 INDIA

☎ 9810052732, 8130973844, 9810052758
9871406036

✉ architassst@gmail.com

🌐 www.architassociates.com



COMPLETE LAB SOLUTIONS FOR **LIME POWDER** MANUFACTURERS IN INDIA



QUALITY
ASSURED



ACCURATE &
RELIABLE



LABORATORY
SOLUTIONS



CUSTOMIZED
EQUIPMENT

PRECISION INSTRUMENTS FOR QUALITY, CONSISTENCY & COMPLIANCE

LAB OVEN



Uniform Heating
Accurate Temperature Control

MOISTURE BALANCE



Fast & Accurate
Moisture Determination

WEIGHING BALANCES



High Precision
Reliable Weighing

MUFFLE FURNACE



High Temperature
Reliable Performance

SLAKE REACTIVITY APPARATUS



For Lime Reactivity Testing
As per Standards

JAW CRUSHER



Crushing of Lumps
for Fine Grinding

SIEVES SHAKER



Efficient Sieving
for Particle Analysis

BRASS SIEVES 20 cm DIA



High Quality Brass Sieves
For Accurate Analysis

RIFFLE SAMPLE DIVIDER



Uniform & Representative
Sample Division

PH METER



Accurate pH
Measurement

LAB GLASSWARE



High Quality Glassware
for Laboratory Use

WHY CHOOSE OUR INSTRUMENTS?

- ✓ High Accuracy & Precision
- ✓ Durable & Long Lasting
- ✓ Easy to Operate & Maintain
- ✓ Compliance with International Standards
- ✓ Ideal for Lime Quality Control Laboratories

APPLICATIONS IN LIME POWDER MANUFACTURE



Raw Material
Testing



Process Control
& Quality Assurance



Product Quality
Evaluation



R & D
Laboratories



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844 E-mail: architasst@gmail.com



Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC

LABORATORY OVEN THIN FILM / LOSS ON HEATING

(FOR BITUMEN TESTING)

Designed to determine the loss in the weight of bitumen & flux oils (Loss on Heating Test) and the effect of heat & air on Asphaltic Materials (Thin Film Test). The unit is heavily insulated and has a double glass built-in door for view inside.

FEATURES



Temperature Controlled

Temperature is controlled at $163^{\circ}\text{C} \pm 1^{\circ}\text{C}$.



Rotating Platform

Inside platform rotated by external motor from 5 to 6 R.P.M.



Heavily Insulated

Superior insulation with double glass built-in door for clear view inside.



Air Circulating Fan

Ensures uniform temperature distribution inside the chamber.



Durable Construction

Inner chamber made of stainless steel sheet and outside made of Mild Steel Sheet.



COMPLETE WITH

- 2 Stainless Steel Pans (to hold cups)
- 4 Aluminum Sample Containers
- 9 Aluminum Sample Containers
- Digital Temperature Controller-Cum-Indicator
- Air Circulating Fan



TECHNICAL SPECIFICATIONS

Application	Thin Film Test / Loss on Heating Test (For Bitumen Testing)
Temperature	$163^{\circ}\text{C} \pm 1^{\circ}\text{C}$
Platform Speed	5 to 6 R.P.M.
Construction	Inner Chamber: Stainless Steel Sheet Outer Body: Mild Steel Sheet
Door	Double Glass Built-in Door for View
Insulation	Heavily Insulated
Controller	Digital Temperature Controller-Cum-Indicator
Air Circulation	Air Circulating Fan
Sample Containers	4 Aluminum Containers & 9 Aluminum Containers
Pans	2 Stainless Steel Pans (to hold cups)
Power Supply	230 Volts AC, 50Hz
Working	Designed to determine the loss in the weight of bitumen & flux oils and the effect of heat & air on Asphaltic Materials
Chamber Size (W x D x H)	450 x 450 x 450 mm

IDEAL FOR



BITUMEN
TESTING



ASPHALTIC
MATERIAL TESTING



QUALITY CONTROL
LABORATORIES



RESEARCH &
DEVELOPMENT

SAMPLE CONTAINERS INCLUDED

4 Aluminum
Containers



9 Aluminum
Containers



ACCURATE
RESULTS



HEAVY DUTY
CONSTRUCTION



UNIFORM AIR
CIRCULATION



EXTERNAL MOTOR
DRIVEN PLATFORM



PRECISE TEMPERATURE
CONTROL

ARCHIT ASSOCIATES - YOUR PARTNER IN TESTING EXCELLENCE



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844 | E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F22C



SCAN TO VISIT
OUR WEBSITE

MOBILE LAB TESTING

ACCURATE TESTING, ANYWHERE, ANYTIME

Our Mobile Laboratory is a fully equipped testing solution that brings the laboratory to your site. Ideal for on-site testing of

**SOIL • CEMENT • CONCRETE
BITUMEN • PAINT • WATER & MORE**



ON-SITE
TESTING



ACCURATE &
RELIABLE



SAVES TIME
& COST



QUALITY
ASSURANCE



AS PER IS / ASTM
STANDARDS



EXPERT TEAM
SUPPORT

EQUIPPED FOR ALL MAJOR TESTS

SOIL TESTING

- Sieve Analysis
- Moisture Content
- Proctor Compaction
- CBR Test
- Field Density Test



CEMENT TESTING

- Consistency Test
- Setting Time Test
- Fineness Test
- Soundness Test
- Compressive Strength



CONCRETE TESTING

- Slump Test
- Cube Compression Test
- Rebound Hammer Test
- Core Cutting
- NDT Tests



BITUMEN TESTING

- Flash & Fire Point
- Softening Point
- Penetration Test
- Ductility Test
- Viscosity Test



PAINT TESTING

- Viscosity
- Fineness of Grind
- Drying Time
- Adhesion Test
- Thickness Test



WATER TESTING

- pH Value
- Turbidity
- TDS / Conductivity
- Hardness
- Chloride Test



WOOD TESTING

- Moisture Content
- Density
- Compression Test
- Bending Strength
- Tensile Resistance



OTHER TESTING

- Aggregate Test
- NDT Test
- Plate Load Test
- Field Density Test
- Many More...



WHY CHOOSE OUR MOBILE LAB?

- ✓ Fully Equipped & Calibrated Lab
- ✓ On-Site Instant Results
- ✓ High Precision & Accuracy
- ✓ Cost Effective Solution
- ✓ Quick Turnaround Time
- ✓ Experienced Professionals
- ✓ Customized Solutions
- ✓ Digital Reporting



APPLICATION AREAS



ROADS & HIGHWAYS



BUILDINGS &
CONSTRUCTION



BRIDGES & FLYOVERS



DAMS & CANALS



INDUSTRIAL
PROJECTS



MINING &
QUARRIES

WE ASSURE

- QUALITY TESTING
- TIMELY DELIVERY
- COMPLETE SATISFACTION
- SAFETY & COMPLIANCE

OUR MOBILE LAB SERVICE PROCESS



REQUIREMENT
ANALYSIS



SITE VISIT &
PLANNING



ON-SITE
TESTING



DATA ANALYSIS
& REPORTING



DELIVERY OF
RESULTS



CALL US TODAY!

8130973844

ARCHIT ASSOCIATES

A-242, Jahangir Puri, Delhi - 110033

Email: architasst@gmail.com

**QUALITY TESTING.
STRONGER INFRASTRUCTURE.**

For Better Testing Today, For Better Tomorrow.



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: +9-8130973844 9810052732, 9810052758, 9871406036

E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



SCAN TO VIEW
PRODUCT RANGE

MOISTURE BALANCE

ACCURATE MOISTURE ANALYSIS FOR PRECISE RESULTS

AVAILABLE MODELS: 60g | 100g | 200g

Moisture Balance is designed for fast, precise and reliable moisture determination in a wide variety of materials. Ideal for use in laboratories, quality control, research, food, pharmaceuticals and industrial applications.

FEATURES

- High Accuracy**
Precise moisture measurement with readability of 0.01%.
- Halogen Heating Technology**
Fast and uniform heating ensures accurate and repeatable results.
- Digital Display**
Large backlit LCD display for easy reading of results and parameters.
- Multiple Modes**
Standard drying, fast drying, soft drying and step drying modes.
- Auto Tare & Calibration**
Built-in auto tare and external calibration for reliable performance.
- Compact & Durable Design**
Robust construction for long life and consistent results.



APPLICATIONS



FOOD & AGRICULTURE



PHARMACEUTICAL INDUSTRY



CHEMICAL & INDUSTRIAL



BUILDING MATERIALS



RESEARCH & QUALITY CONTROL

TECHNICAL SPECIFICATIONS

Parameter	MB-60G	MB-100G	MB-200G
Capacity	60g	100g	200g
Readability		0.01% (0.001g)	
Moisture Readability		0.01%	
Heating Source		Halogen Lamp	
Temperature Range		40°C to 199°C (1°C Step)	
Display		Backlit LCD with Dual Display	
Weighing Sensor		High Precision Load Cell	
Drying Modes		Standard, Fast, Soft, Step	
Timer		1 to 99 minutes (1 min. Step)	
Pan Size		Ø90 mm	
Power Supply		230V ± 10%, 50Hz	
Interface		RS-232 (Optional)	
Calibration		External Calibration	
Operating Temperature		10°C to 40°C	

OPTIONAL ACCESSORIES



ALUMINUM SAMPLE PANS



GLASS FIBER FILTERS



PRINTER



CALIBRATION WEIGHTS



RS-232 INTERFACE

STANDARD FEATURES

- Halogen Heating System
- Large Backlit LCD Display
- % Moisture, % Dry Matter, ATRO (Auto)
- Automatic Temperature Control
- Auto Power Off Function
- Overload Protection
- Pan Handle for Safe Handling
- In-built Level Indicator
- Draft Shield Ring (Optional)
- RS-232 Interface (Optional)

AVAILABLE MODELS

60g MODEL

Ideal for small sample quantity & precision applications.

100g MODEL

Ideal for general laboratory & quality control applications.

200g MODEL

Suitable for high sample volume and industrial use.



ACCURATE & RELIABLE RESULTS



FAST & EASY OPERATION



UNIFORM HEATING



RUGGED & DURABLE CONSTRUCTION



WIDE RANGE OF APPLICATIONS



QUALITY YOU CAN TRUST

ARCHIT ASSOCIATES - YOUR PARTNER IN TESTING EXCELLENCE



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir,
Delhi-110 033 (INDIA)

91-8130973844, 9871406036

architasst@gmail.com

www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC

TRUSTED SUPPLIER OF
MATERIAL TESTING &
LABORATORY INSTRUMENTS



QUALITY
ASSURED
PRODUCTS



ACCURATE &
RELIABLE



LABORATORY &
SITE SOLUTIONS



CUSTOMIZED
TESTING
EQUIPMENT

PLANETARY / MORTAR MIXER

WITH DIGITAL RPM INDICATOR

COMPLIANCE WITH FOLLOWING INTERNATIONAL STANDARDS

EN 196 1 & EN 196 3, IS : 4031, IS 10890,
IS : 1727, ASTM C451

Planetary / Mortar Mixer is used for mixing of cement mortar and construction materials in civil engineering laboratories. It is designed for uniform and consistent mixing as per standard test requirements.

CONSTRUCTION DETAILS

1. Compact and tactful design
2. Aesthetic appearance
3. Exterior body powder coated in attractive shades
4. Paddle and bowl made of Stainless Steel Material
5. An epicyclic type stainless steel paddle imparts both planetary and revolving motion.
6. The blade is about 75 mm long, 50 mm wide and tapered to a thin edge.
7. Mixer operates on two different speeds.
Low speed (140 ± 5 rpm) and Medium speed (285 ± 10 rpm).
8. Planetary movement of (62 ± 5 rpm) in low range and (125 ± 10 rpm) in medium range.
9. The mixing paddle has a planetary motion and is driven by a motor.
10. The scraper consists of semi rigid rubber blade attached to a handle about 150 mm long.
11. Suitable for operation on 230 Volts, 50 Hz, single phase, AC Supply.

BASIC SPECIFICATIONS

Parameter	Details
Product	Planetary / Mortar Mixer
Standard	As per BIS Standards
Mixing Bowl	Stainless Steel
Mixing Action	Planetary & Revolving Motion
Speeds	Low Speed : 140 ± 5 rpm Medium Speed : 285 ± 10 rpm
Planetary Movement	Low Range : 62 ± 5 rpm Medium Range : 125 ± 10 rpm
Blade Size	75 mm (L) x 50 mm (W), Tapered Edge
Motor	Heavy Duty
Power Supply	230 Volts, 50 Hz, Single Phase, AC
Operation	Motorized
Finish	Powder Coated
Application	Mixing of Cement Mortar & Construction Materials

SALIENT FEATURES

- ✓ Compact and tactful design
- ✓ Aesthetic appearance
- ✓ Stainless steel paddle and bowl
- ✓ Epicyclic type paddle for uniform mixing
- ✓ Two speed operation (140 ± 5 rpm & 285 ± 10 rpm)
- ✓ Planetary movement for better mixing consistency
- ✓ Scraper with semi-rigid rubber blade
- ✓ Digital RPM indicator for accurate speed monitoring
- ✓ Powder coated body for corrosion resistance
- ✓ Suitable for laboratory use
- ✓ Easy to operate and maintain



WIDELY USED IN



CEMENT
LABORATORIES



CIVIL ENGINEERING
LABORATORIES



QUALITY CONTROL
DEPARTMENTS



RESEARCH &
INSTITUTES



CONSTRUCTION
INDUSTRY



ARCHIT ASSOCIATES
A-242, Jahangir Puri,
Near Kali Mandir,
Delhi-110 033 (INDIA)



91-8130973844
9871406036



architasst@gmail.com

ALL PRODUCTS COMPLY WITH RELEVANT INDUSTRY STANDARDS



ARCHIT ASSOCIATES

Trusted Supplier of
MATERIAL TESTING & LABORATORY INSTRUMENTS

A-242, Jahangir Puri, Near Kali Mandir,
Delhi-110 033 (INDIA)

91-8130973844, 9871406036,
9810052732, 9810052758

architasst@gmail.com

www.architassociates.com

GSTIN NO. 07A0UPG0300F22C

TRUSTED SUPPLIER OF
MATERIAL TESTING &
LABORATORY INSTRUMENTS



QUALITY
ASSURED
PRODUCTS



ACCURATE &
RELIABLE



LABORATORY &
SITE SOLUTIONS



CUSTOMIZED
TESTING
EQUIPMENT

CONCRETE TESTING MOULDS

ISI & NON ISI MODELS AVAILABLE

CUBE MOULDS – MODELS AVAILABLE

CAST IRON CUBE MOULDS (IS:10086 / NON-ISI)



50 mm Cube
Cast Iron
NON-ISI



100 mm Cube
Cast Iron
ISI Mark / NON-ISI




150 mm Cube
Cast Iron
ISI Mark / NON-ISI

STEEL CUBE MOULDS



70.6 mm Cube
Steel Mould
ISI Mark / NON-ISI

BEAM MOULDS

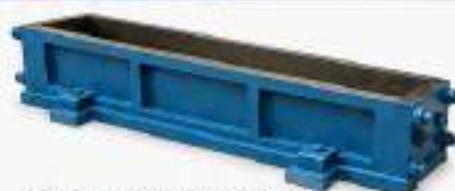

100 mm x 100 mm x 500 mm
Beam Mould (Cast Iron)



150 mm x 150 mm x 700 mm
Beam Mould (Cast Iron)

CYLINDRICAL MOULDS

CAST IRON CYLINDRICAL MOULDS (IS:10086)



150 mm Dia
x 300 mm High
(IS:10086)



100 mm Dia
x 100 mm High
(Split)



100 mm Dia
x 200 mm High
(Split)



150 mm Dia
x 150 mm High
(Split)



300 mm Dia
x 300 mm High
(Split)

PLASTIC CUBE MOULDS



100 mm Cube
Plastic Mould



150 mm Cube
Plastic Mould

KEY FEATURES

- ✓ Made from high quality cast iron / steel / plastic
- ✓ Conforms to IS: 10086 (where applicable)
- ✓ Precise dimensions for accurate test results
- ✓ Strong, durable and long service life
- ✓ Available in ISI & Non-ISI models
- ✓ Easy to use and maintain

APPLICATIONS



CEMENT
INDUSTRY



CONSTRUCTION
QUALITY CONTROL



RESEARCH
LABORATORIES



ENGINEERING
COLLEGES



TESTING
LABORATORIES

ALL PRODUCTS COMPLY WITH RELEVANT INDIAN STANDARDS



A-242, Jahangir Puri,
Near Kali Mandir,
Delhi-110 033 (INDIA)



91-8130973844
9871406036
9810052732
9810052758



architasst@gmail.com



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir,
Delhi-1100 033 (INDIA)

+91-9810052732, 9810052758, 9871406036

architast@gmail.com

Web www.architassociates.com



TESTING EXCELLENCE IN EVERY RESULT

DIGITAL COMPRESSION TESTING MACHINE CHANNEL MODEL

WITH AUTOMATIC PACE RATE CONTROL SYSTEM

AVAILABLE CAPACITIES

- 100 kN
- 250 kN
- 500 kN
- 1000 kN
- 2000 kN
- 3000 kN



SUPPLIED WITH
NABL CALIBRATION
CERTIFICATE



Digital Display with
Microprocessor Based Controller



Automatic Pace Rate Control for
Uniform & Accurate Loading



High Precision, Repeatability
& Reliability



Robust Construction for
Long Lasting Performance



TECHNICAL SPECIFICATIONS

Capacity (kN)	100, 250, 500, 1000, 2000 & 3000 kN
Model	: Channel Model
Display	: Digital
Pace Rate Control	: Automatic
Control Unit	: Microprocessor Based
Load Accuracy	: $\pm 1\%$ of Indicated Load
Load Resolution	: 0.1 kN
Platen Size	: 165 mm (Standard) (or as per requirement)
Power Supply	: 230V AC, 50Hz, Single Phase



ACCURATE
RESULTS



AUTOMATIC PACE
RATE CONTROL



EASY TO
OPERATE



SAFE & RELIABLE
PERFORMANCE



LOW MAINTENANCE
& DURABLE



CALL US TODAY
FOR EXPERT
SOLUTIONS

9810052732
9810052758
9871406036

ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir,
Delhi-1100 033 (INDIA)

+91-9810052732, 9810052758, 9871406036

architast@gmail.com



**BUILT FOR LABS.
ENGINEERED FOR EXCELLENCE.**
Precision. Performance. Perfection.



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)
Phone: 91-8130973844 E-mail: architasst@gmail.com

Web: www.architassociates.com
GSTIN NO. 07A0UPG0300F2ZC



SCAN TO VIEW
PRODUCT RANGE

MUFFLE FURNACE

High Temperature Performance
For Laboratory Applications

AVAILABLE MODELS:

1150 – 1200°C 1350 – 1400°C 1550 – 1600°C

Designed for a wide range of laboratory applications including ashing, heat treatment, calcination, sintering and general high temperature testing. Built for accuracy, safety and durability.

FEATURES

- High Temperature Range**
Available in 1150–1200°C, 1350–1400°C & 1550–1600°C models.
- Digital PID Controller**
Microprocessor based PID controller for precise temperature control and stability.
- High Quality Insulation**
Low thermal conductivity insulation for maximum energy efficiency and safety.
- Durable Construction**
Inner chamber made of high grade refractory material & outer body of mild steel.
- Uniform Temperature**
Specially designed heating elements on three sides ensure uniform temperature distribution.
- Safety Features**
Over temperature protection, door safety interlock & indicator lights.



APPLICATIONS

- ASHING & IGNITION LOSS
- HEAT TREATMENT
- SINTERING & CALCINATION
- CERAMIC & GLASS TESTING
- RESEARCH & DEVELOPMENT

TECHNICAL SPECIFICATIONS

Model	1150 – 1200°C	1350 – 1400°C	1550 – 1600°C
Max. Temperature	1150 – 1200°C	1350 – 1400°C	1550 – 1600°C
Working Temperature	1100°C	1300°C	1500°C
Temperature Controller	Digital PID Controller – Cum – Indicator		
Temperature Accuracy	± 1°C		
Heating Elements	Kanthal A-1 / Super Kanthal		
Chamber Material	High Alumina Refractory Bricks	High Alumina Refractory Bricks	Ceramic Fiber Board / High Alumina Bricks
Insulation	High Grade Ceramic Fiber Insulation		
Outer Body	Powder Coated Mild Steel Sheet		
Power Supply	230 Volts AC, 50 Hz		
Optional	S.S. Chamber, Digital Timer, Data Logger, Over Temperature Controller, Programmable Controller		

STANDARD FEATURES

- Digital PID Temperature Controller-Cum-Indicator
- High Temperature Refractory Chamber
- High Grade Ceramic Fiber Insulation
- Heating Elements on Three Sides
- Durable Powder Coated Body
- Door Safety Interlock
- Over Temperature Protection
- Power Indicator, Heating Indicator
- Easy to Operate & Maintain

AVAILABLE MODELS

- 1150 – 1200°C MODEL**
Suitable for general ashing & heat treatment applications.
- 1350 – 1400°C MODEL**
Suitable for high temperature processing & research work.
- 1550 – 1600°C MODEL**
Suitable for advanced thermal processes & special materials.

OPTIONAL ACCESSORIES



PRECISE
TEMPERATURE
CONTROL

ROUGH & DURABLE
CONSTRUCTION

HIGH THERMAL
EFFICIENCY

SAFE & RELIABLE
OPERATION

EASY TO USE &
MAINTAIN

LONG LIFE &
LOW MAINTENANCE

ARCHIT ASSOCIATES – YOUR PARTNER IN TESTING EXCELLENCE



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)
Phone: 91-8130973844 E-mail: architasst@gmail.com

Web: www.architassociates.com
GSTIN NO. 07A0UPG0300F2ZC



SCAN TO VIEW
PRODUCT RANGE

LABORATORY OVEN DIGITAL / ANALOG TEMPERATURE

AVAILABLE MODELS: 250°C | 300°C | 400°C

Designed for general laboratory applications including drying, heating, sterilization and other temperature controlled processes. Ideal for use in educational, research, industrial and quality control laboratories.

FEATURES

-  **Wide Temperature Range**
Available in 250°C, 300°C & 400°C variants.
-  **Digital / Analog Controller**
PID Digital Controller-Cum-Indicator for accurate temperature control.
-  **Heavy Duty Construction**
Inner chamber made of Stainless Steel & outside made of Mild Steel.
-  **Uniform Heating**
Air circulating fan ensures uniform temperature distribution.
-  **Double Walled Door**
Insulated door with double wall construction for better efficiency.
-  **Safe & Reliable**
Over temperature protection & indicator for enhanced safety.



APPLICATIONS



TECHNICAL SPECIFICATIONS

Model	250°C	300°C	400°C
Temperature Range	Ambient to 250°C	Ambient to 300°C	Ambient to 400°C
Temperature Controller	Digital PID Controller-Cum-Indicator (Also available in Analog)		
Temperature Accuracy	± 1°C		
Temperature Sensor	PT-100		
Construction	Inner Chamber: Stainless Steel Sheet Outer Body: Mild Steel Sheet		
Insulation	High Grade Mineral Wool / Glass Wool		
Door	Double Walled & Insulated		
Air Circulation	Air Circulating Fan		
Shelves	Removable Wire Mesh / Perforated Shelves		
Power Supply	220 Volts AC, 50Hz		
Safety	Over Temperature Cut-off & Indicator		

OPTIONAL ACCESSORIES



STANDARD FEATURES

- Digital Temperature Controller-Cum-Indicator (Analog also available)
- Uniform Temperature Distribution
- Stainless Steel Chamber
- Insulated Double Walled Door
- Air Circulating Fan
- Powder Coated Exterior
- Stable & Vibration Free Performance

AVAILABLE MODELS

250°C MODEL

Suitable for general drying & heating applications.

300°C MODEL

Suitable for high temperature processing.

400°C MODEL

Suitable for high temperature & industrial use.



ARCHIT ASSOCIATES - YOUR PARTNER IN TESTING EXCELLENCE



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844

E-mail: architassst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



SCAN TO VIEW
PRODUCT RANGE

PAINT & COATING TESTING EQUIPMENT

ACCURATE • RELIABLE • PRECISION ENGINEERED

A comprehensive range of instruments for testing the quality, performance and durability of paints, coatings and related materials.



PRECISE
RESULTS



RUGGED &
DURABLE



ACCURATE &
RELIABLE



EASY TO USE
& MAINTAIN



IDEAL FOR LAB
& INDUSTRY



OUR PRODUCT RANGE

SALT SPRAY CHAMBER



- Corrosion testing of coated metals & components
- PVC / FRP / RP construction
- Digital PID Controller
- Complies with ASTM B117, ISO 9227, IS 10121

COLOUR VIEWING CABINET



- Assessment of colour under standard light
- D65, TUBA, UV, CWF, FVA
- Digital hour meter
- Complies with ASTM D1729

FLOW CUP



- Ford Cup 1, 2, 3, 4, 5, 6
- For viscosity measurement of paints & coatings
- Stainless Steel / Brass
- Complies with ASTM D1206, IS 7231

CONICAL MANDREL BEND TESTER



- To determine flexibility of coatings
- Conical Mandrels: 1/8" to 1"
- Complies with ASTM D622, ISO 5860, IS 101

TUBULAR IMPACT TESTER



- To determine impact resistance of coatings
- Adjustable height
- Complies with ASTM D2754, IS 101

CROSS CUT TESTER



- Adhesion testing of coatings
- 6 Teeth / 11 Teeth Options
- Complies with ASTM D3359, ISO 2409, IS 101

DFT COATING THICKNESS GAUGES



- Magnetic (Fe) & Non-Magnetic (NFe)
- Range: 0 - 2000 µm
- Digital Display
- Complies with ISO 2808

GLOSS METER



- 20° / 60° / 85° Options
- Measure gloss of coatings
- Digital Display
- Complies with ASTM D623, ISO 2813

PENCIL HARDNESS TESTER



- Determine scratch hardness of coatings
- Pencil: 6B to 6H
- Complies with ASTM D3363, IS 101

ERICHSEN SCRATCH HARDNESS TESTER



- Scratch resistance of coatings
- Adjustable load
- Complies with ISO 1518, ASTM D2197

PULL OFF ADHESION TESTER



- Measures adhesion strength
- Digital display
- Complies with ASTM D4541, ISO 4624

WET FILM THICKNESS GAUGE



- Measure wet film thickness
- Range: 2K - 2000 µm
- Stainless Steel
- Complies with ISO 2008

ACCESSORIES



Paint Shaker



Viscosity Cup



Mixing Paddle



Applicator Film (Various Types)



Elcometer Coating Gauge



Humidity Chamber



Drying Oven



Standard Panels (SS / MS / Aluminum)

IDEAL FOR

- Paint & Coating Manufacturers
- Quality Control Laboratories
- Research & Development
- Automotive Industry
- Infrastructure & Engineering
- Educational & Testing Institutes

ARCHIT ASSOCIATES

A-242, Jahangir Puri,
Near Kali Mandir, Delhi-110 033 (INDIA)

91-8130973844

architassst@gmail.com

GSTIN NO. 07A0UPG0300F2ZC



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844

E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



SCAN TO VIEW
PRODUCT RANGE

CONCRETE PAN MIXER

ROBUST • RELIABLE • EFFICIENT

Designed for uniform and rapid mixing of concrete, mortar and other construction materials. Suitable for laboratory and site use.



Heavy Duty Construction



Uniform & Homogeneous Mixing



Easy Operation &
Low Maintenance



Sturdy Pan with
High Mixing Efficiency



Ideal for Laboratory
& Field Applications



Reliable Performance
& Long Service Life

AVAILABLE CAPACITIES

CAP. 40 LTR.



Drum Capacity	40 Liters
Mixing Capacity	28 - 30 Liters
Power	0.5 HP / 1440 RPM
Weight (Approx.)	65 - 70 kg

CAP. 60 LTR.



Drum Capacity	60 Liters
Mixing Capacity	42 - 45 Liters
Power	0.75 HP / 1440 RPM
Weight (Approx.)	85 - 90 kg

CAP. 100 LTR.



Drum Capacity	100 Liters
Mixing Capacity	70 - 75 Liters
Power	1 HP / 1440 RPM
Weight (Approx.)	120 - 130 kg



Consistent
Mixing Quality



Wear Resistant
Pan



Strong & Durable
Construction



Smooth & Easy
Tilting Mechanism



Tested for
Reliable Results



IDEAL FOR

- Construction Laboratories
- Civil Engineering Projects
- Quality Control Laboratories
- RMC Plants & Sites

ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844

E-mail: architasst@gmail.com

GSTIN NO. 07A0UPG0300F2ZC



ARCHIT ASSOCIATES

YOUR TRUSTED PARTNER IN LABORATORY & MATERIAL TESTING

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

+91-8130973844, 9871406036 architasst@gmail.com

www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



ACCURATE
RESULTS



RELIABLE
PERFORMANCE



QUALITY
ASSURED

TESTING EXCELLENCE IN EVERY RESULT

BRIQUETTE MOULD (THREE GANG & SINGLE GANG MOULD)

FOR TENSILE TESTING OF CEMENT, MORTAR & BITUMEN

Highly durable, precision-engineered Briquette Moulds used for the tensile testing of cement, mortar, and bitumen.

They are designed to prevent distortion, maintain strict dimensional compliance, and feature a superior internal surface finish in accordance with IS: 10080 and IS: 10086 specifications.



THREE GANG MOULD



SINGLE GANG MOULD



THREE GANG MOULD
with grinding

KEY FEATURES

- Manufactured from high-grade cast iron / steel for maximum strength and durability
- Precision machined cavities ensure accurate and consistent briquette dimensions
- Superior internal surface finish prevents adherence and ensures easy demoulding
- Designed to prevent distortion and maintain dimensional stability
- Robust clamping mechanism for leak-proof and secure moulding
- Corrosion-resistant finish for long service life
- Ergonomic handle for safe and convenient use
- Complies with IS: 10080 and IS: 10086

TECHNICAL SPECIFICATIONS

Parameter	Three Gang Mould	Single Gang Mould
Number of Cavities	3	1
Internal Dimensions (L x W x H) (Each Cavity)	165 x 40 x 40 mm	165 x 40 x 40 mm
Overall Dimensions (L x W x H) (Approx.)	330 x 150 x 120 mm	120 x 150 x 120 mm
Material	Mild Steel	
Finish	Precision machined, smooth internal finish	
Clamping System	Wing nut with screw arrangement	
Application	Tensile testing of cement, mortar & bitumen briquettes	
Standards	IS: 10080, IS: 10086	

STANDARDS

IS: 10080

Specification for Briquettes for Tensile Strength of Cement Mortar

IS: 10086

Specification for Briquette Moulds for Cement and Mortar

SCOPE OF SUPPLY

- Briquette Mould (Three Gang / Single Gang) - 1 No.
- Base Plate - 1 No.
- Clamping Assembly - 2 Nos.
- Instruction Manual - 1 No.

THREE GANG MOULD (CAST IRON / STEEL)



Designed to cast three briquettes simultaneously, improving efficiency and ensuring uniformity in testing. Ideal for laboratories with high testing volume requirements.

Ideal For:

- Cement & mortar testing labs
- Quality control & research labs
- Educational & industrial applications

SINGLE GANG MOULD (CAST IRON / STEEL)



Compact and easy-to-handle mould for casting single briquettes. Ideal for routine testing and low-volume requirements.

Ideal For:

- Small laboratories
- Field & site testing
- Educational institutions

THREE GANG MOULD (STAINLESS STEEL)



Steel construction ensures excellent corrosion resistance, longer life and superior performance in all environments.

Ideal For:

- High-end laboratories
- Harsh environment testing
- Long-term, heavy-duty use



PRECISION ENGINEERED
Dimensional accuracy
& reliable results



HIGH DURABILITY
Strong construction
for long service life



EASY TO USE
User-friendly design
for safe operation



CORROSION RESISTANT
Protective finish for
harsh environments



TEST CERTIFICATES
Supplied with
calibration certificate



SCAN TO VISIT
OUR WEBSITE

ARCHIT ASSOCIATES

A-242, Jahangir Puri, Delhi - 110033 (INDIA)

+91-8130973844, +91-9871406036

architask@gmail.com https://www.google.com/maps/place/architassociates

28
YEARS OF
EXPERIENCE



BUILT FOR LABS.
ENGINEERED FOR EXCELLENCE.
Precision, Performance, Perfection.



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844

E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



SCAN TO VIEW
PRODUCT RANGE

ROCKBOLT PULLOUT TEST EQUIPMENT

ACCURATE • RELIABLE • RUGGED • EASY TO OPERATE

Rockbolt Pullout Test Equipment is used for the testing of Rockbolts or Rock Anchors pullout strength in rock, for the evaluation of bond between reinforcement and grout cement or grout concrete or rock.

Manufactured to suit different site requirements such as manual or motorized operation, size of rockbolt and also different length of flexible hose pipes are available.

APPLICATIONS


- Testing of Rockbolts / Rock Anchors pullout strength
- Evaluation of bond between reinforcement and grout
- Quality control in tunneling and mining projects
- Slope stabilization and underground support systems
- Research and geotechnical investigations



TUNNELING PROJECTS



SLOPE STABILIZATION

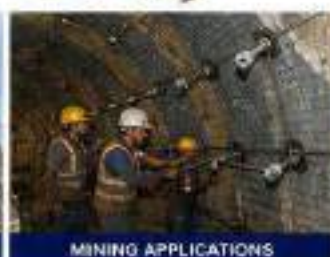

MINING APPLICATIONS



GEOTECHNICAL INVESTIGATIONS

TECHNICAL SPECIFICATIONS

PARAMETER	20T MODEL	30T MODEL	50T MODEL
Hydraulic Centre Hole Jack Range	20T	30T	50T
Centre Hole Diameter	33mm	33mm	54mm
Approx. Lift of Jack	65mm	65mm	65mm
Approx. Weight of Jack	14Kgs	16Kgs	25Kgs
Pumping Unit Type	Manual / Motorized, fitted with load gauge		

KEY FEATURES

- ✓ High accuracy and reliability
- ✓ Compact & portable design
- ✓ Easy to operate and maintain
- ✓ Heavy duty construction
- ✓ Available in manual or motorized options
- ✓ Custom hose length available
- ✓ Suitable for field and laboratory use

TYPICAL TEST PROCEDURE



ACCURATE
RESULTS



RUGGED &
DURABLE



EASY TO
OPERATE



BUILT FOR INDIAN
CONDITIONS



AFTER SALES
SUPPORT



ARCHIT ASSOCIATES

A-242, Jahangir Puri,
Near Kali Mandir, Delhi-110 033 (INDIA)
+91-8130973844

✉ architasst@gmail.com

GSTIN NO. 07A0UPG0300F2ZC

ARCHIT ASSOCIATES - YOUR PARTNER IN TESTING EXCELLENCE



ARCHIT ASSOCIATES

YOUR TRUSTED PARTNER IN LABORATORY & MATERIAL TESTING

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

91-8130973844

architasst@gmail.com

www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



ACCURATE
RESULTS



RELIABLE
PERFORMANCE



QUALITY
ASSURED

RAILWAY INSPECTION GAUGES

PRECISION GAUGES FOR SAFE, RELIABLE & EFFICIENT RAIL OPERATIONS

At ARCHIT ASSOCIATES, we provide a comprehensive range of Railway Inspection Gauges designed to meet stringent Indian Railways and International standards.

Our gauges are manufactured from high quality materials for dimensional accuracy, durability and long service life. These are essential tools for inspection, measurement and maintenance of railway tracks, rolling stock and associated components, ensuring operational safety, reliability and compliance.



OUR RANGE OF RAILWAY INSPECTION GAUGES



CHECKING GAUGES



BENCH BOLT & NUT



BENCH NUT BOLTS



END PLATE HOLE CHECKING
GAUGE WITH MASTER



GANTRY TROLLEY



HOIST HOOKS



HOIST WHEEL
ASSEMBLY



RT-381 INSERT
CHECKING GAUGE



MOTOR COUNTER WEIGHT



MOULD CHECKING
GAUGES



HOIST & TROLLEY WHEELS



PLAIN & GEARED
TROLLEY WHEEL



RT-3705 INSPECTION
GAUGE (3)



SLEEPER CHECKING
GAUGES



STENCILS



STENCILS AS PER
REQUIREMENT



TOE PIN GAP GAUGE



TOE PIN GAUGE &
MASTER GAUGE



TOE PIN GAUGE RT-381



VIBRATOR COVER 2 KW



WIRE POSITION CHECKING
GAUGE

STANDARDS COMPLIANCE



As per Indian Railways
Specifications
& International
Standards



SCAN TO VISIT
OUR WEBSITE

ARCHIT ASSOCIATES

A-242, Jahangir Puri, Delhi - 110033 (INDIA) 91-8130973844

architasst@google.com/view/architassociates

28 YEARS OF
EXPERIENCE

DELIVERING QUALITY EQUIPMENT
FOR ACCURATE RESULTS AND
BETTER PERFORMANCE.



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844

E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC

TESTING
EXCELLENCE IN
EVERY RESULT

RING AND BALL APPARATUS – SEMI AUTOMATIC & DIGITAL

For Softening Point Determination
of Bitumen Materials

FEATURES & BENEFITS



Digital Temperature Controller
Highly accurate PID based
temperature control



Precise Heating System
Uniform heating for
consistent results



Laboratory Grade Construction
Rust proof & corrosion resistant
high quality materials



Accurate Softening Point Testing
Ensures reliable and
repeatable test results



User Friendly Operation
Simple controls for
ease of operation



**Suitable for Civil Engineering
Laboratories**
Ideal for Bitumen & Asphalt
material testing



SEMI AUTOMATIC MODEL

Reliable automatic model with
PID temperature control



DIGITAL MODEL

Advanced digital model with
PID temperature controller
for high accuracy

CONFORMS TO
RELEVANT
IS / ASTM
STANDARDS

- ✓ IS 1205 Part II
- ✓ ASTM D36
- ✓ AASHTO T53



APPLICATION

Used for determining the Softening Point
of Bitumen which is an important
parameter in evaluating the performance
of bituminous materials.



ACCURATE
RESULTS



DURABLE
CONSTRUCTION



EASY TO
OPERATE



LABORATORY
GRADE



SCAN TO VISIT
OUR WEBSITE

ARCHIT ASSOCIATES

A-242, Jahangir Puri, Delhi – 110033

Phone: 8130973844

Email: architasst@gmail.com

QUALITY INSTRUMENTS

ACCURATE TESTING | RELIABLE RESULTS
BUILT FOR LABS, ENGINEERED FOR EXCELLENCE



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844

E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



LEADING SUPPLIER OF LABORATORY & MATERIAL TESTING EQUIPMENT FOR RMC PLANTS

Complete BIS Standard Testing Solutions for
Concrete, Cement, Aggregates, Soil, Asphalt &
Construction Materials



ABOUT COMPANY

Archit Associates is a trusted supplier of high-quality laboratory and material testing equipment for RMC Plants, Civil Engineering Laboratories, Construction Projects, Infrastructure Companies and Quality Control Laboratories across India. We supply reliable Indian market available equipment strictly as per BIS Standards.



BIS STANDARDS COMPLIANCE

Equipment supplied strictly as per
relevant BIS / IS Standards applicable
in Indian Construction & RMC Industry.



BIS APPROVED
QUALITY



TESTED &
CERTIFIED



RELIABLE
& ACCURATE



BUILT FOR
PERFORMANCE

QUALITY • ACCURACY • RELIABILITY • COMPLIANCE



OUR PRODUCT CATEGORIES



Concrete Testing
Equipment



Cement Testing
Equipment



Aggregate Testing
Equipment



Soil Testing
Equipment



Asphalt & Bitumen
Testing Equipment



Compression Testing
Machines



Cube Moulds &
Slump Test Apparatus



Sieves & Laboratory
Glassware



Weighing
Balances



Non-Destructive
Testing Equipment



Laboratory Ovens &
Incubators



Survey & Civil Lab
Instruments



INDUSTRIES SERVED

- » RMC Plants
- » Infrastructure Projects
- » Highway Projects
- » Railways
- » Metro Projects
- » Cement Plants
- » Construction Companies
- » QA/QC Laboratories
- » Government Testing Laboratories



WHY CHOOSE ARCHIT ASSOCIATES?

- Genuine Quality Equipment
- Fast Delivery
- Competitive Pricing
- Calibration Support
- Technical Support
- Trusted Industrial Supplier
- PAN India Supply



YOUR TRUSTED PARTNER FOR MATERIAL TESTING & QUALITY ASSURANCE SOLUTIONS



A-242, Jahangir Puri,
Near Kali Mandir,
Delhi-110 033 (INDIA)



+91-8130973844



architasst@gmail.com

GSTIN NO. 07A0UPG0300F2ZC



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844

E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



100% INDIAN MADE
LABORATORY EQUIPMENT



SCAN TO VIEW
PRODUCT RANGE

ROCKWELL & SUPERFICIAL ROCKWELL HARDNESS TESTERS

ACCURATE • RELIABLE • PRECISE • EASY TO OPERATE • BUILT FOR LONG LIFE

AUTOMATIC MOTORIZED DIGITAL ROCKWELL & SUPERFICIAL ROCKWELL HARDNESS TESTER



Automatic Motorized Digital Rockwell and Superficial Rockwell Hardness Tester, specially designed for accurate and reliable hardness measurement in laboratories, quality control departments, production facilities, research institutes, and industrial applications.

RANGE

- Rockwell: 20-95 HRA, 10-100 HRB, 10-70 HRC
- Superficial Rockwell:
 - 70-94 HR15N, 42-86 HR30N, 20-77 HR45N
 - 67-93 HR15T, 29-82 HR30T, 10-72 HR45T

KEY FEATURES

- Multiple scales: HRA, HRB, HRC, HRD, HRF, HRE, HR15N, HR30N, HR45N, HR15T, HR30T, HR45T (with standard indenter)
- Optional scales (with optional indenters): HRE, HRH, HRK, HRL, HRM, HRP, HRH, HRS, HRV, HR15W, HR30W, HR45W, HR15X, HR30X, HR45X, HR15Y, HR30Y, HR45Y
- Preliminary test force: 98N (10kgf) for Rockwell, 29.4N (3kgf) for Superficial Rockwell
- Main test forces: 147N (15kgf), 294N (30kgf), 441N (45kgf), 588N (60kgf), 980N (100kgf), 1471N (150kgf)
- Automatic load control: Initial test force, main test force loading, holding & unloading
- Automatic stage elevation: Up, down & return
- Resolution: 0.1 HR
- Load dwell time: 2 to 60 seconds (adjustable)
- Max. workpiece height: 230 mm
- Max. testing width: 170 mm
- Power supply: 100-240V, 50/60Hz
- Dimensions (LxWxH): 475 x 200 x 700 mm
- Net weight: Approx. 70 kg

STANDARD DELIVERY



IDEAL FOR



MOTORIZED ROCKWELL HARDNESS TESTER



Motorized Rockwell Hardness Tester designed for accurate hardness measurement of metallic materials.

RANGE

- Rockwell: 20-95 HRA, 10-100 HRB, 10-70 HRC
- Scales (Standard Indenter): HRA, HRB, HRC, HRD, HRF, HRE
- Additional scales (Optional Indenters): HRE, HRH, HRK, HRL, HRM, HRP, HRR, HRS, HRV

KEY FEATURES

- Preliminary test force: 98 N (10 kgf)
- Major test forces: 588 N (60 kgf), 980 N (100 kgf), 1471 N (150 kgf)
- Stage elevation: Manual
- Load application, dwell & unloading: Automatic
- Resolution: 0.1 HR
- Load dwell time: 2 to 60 seconds (adjustable)
- Max. workpiece height: 230 mm
- Max. testing width: 170 mm
- Power supply: 100-240V, 50/60Hz
- Dimensions (LxWxH): 475 x 200 x 700 mm
- Net weight: Approx. 70 kg

STANDARD DELIVERY



NOTE:

Hardness Test Blocks (HR15N80-90, HR30N65-80, HR30T70-85) are not included in standard delivery.



ARCHIT ASSOCIATES

A-242, Jahangir Puri,
Near Kali Mandir, Delhi-110 033 (INDIA)
+91-8130973844

architasst@gmail.com

GSTIN NO. 07A0UPG0300F2ZC

ARCHIT ASSOCIATES - YOUR PARTNER IN TESTING EXCELLENCE



SCAN TO VIEW
PRODUCT RANGE



ARCHIT ASSOCIATES

📍 A-242, Jangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

☎ Phone: 91-8130973844

✉ E-mail: architasst@gmail.com

🌐 Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



ACCURATE
RESULTS



RELIABLE
PERFORMANCE



QUALITY
ASSURED



SCAN TO VIEW
PRODUCT RANGE

ROLLING THIN FILM OVEN

UNIFORM HEATING • PRECISE CONTROL • RELIABLE RESULTS

Rolling Thin Film Oven is designed for laboratory determination of rolling thin films of bituminous materials. It ensures uniform heating and controlled test conditions for accurate and repeatable results.

KEY FEATURES



DOUBLE-WALLED INSULATED CHAMBER
Ensures uniform temperature distribution and energy efficiency.



DOUBLE-PANE GLASS VIEWING WINDOW
Front door with double-pane glass for safe observation during operation.



AIR VENTS (TOP & BOTTOM)
Symmetrical top and bottom air vents to maintain proper air circulation.



AIR PLENUM SYSTEM
Integrated air plenum system for uniform heated air distribution inside the chamber.



SQUIRREL CAGE FAN
Squirrel cage type circulation fan (2800 RPM) ensures consistent air movement and temperature uniformity.



DIGITAL CONTROLLER
Digital indicating thermostat / PID controller to maintain oven temperature at $163^{\circ}\text{C} \pm 0.5^{\circ}\text{C}$.



VERTICAL ROTATING CARRIAGE
Vertical circular carriage mechanism designed to rotate sample bottles mechanically at 15 ± 0.2 RPM.



AIR JET SYSTEM
Air jet system that directs heated air into each sample bottle at the lowest point of its travel.



CALIBRATED FLOW METER
Calibrated flow meter provided to regulate air flow at 4000 mL/min.



TECHNICAL SPECIFICATIONS

Temperature Range	10°C to 50°C
Temperature Accuracy	$\pm 1^{\circ}\text{C}$
Humidity Range	50% to 95% RH
Humidity Accuracy	$\pm 3\%$ RH
CO ₂ Concentration Range	1% to 10%
CO ₂ Accuracy	$\pm 0.1\%$
Chamber Size (Standard)	600 x 600 x 600 mm (L x W x H) (Custom sizes available)
Controller	Digital Thermostat / PID Controller
Display	Digital LED Display (Temperature, Humidity, CO ₂)
Circulation Fan	Squirrel Cage Fan, 2800 RPM
Carriage Speed	15 ± 0.2 RPM
Air Flow	4000 mL/min (Calibrated Flow Meter)
Power Supply	230V AC, 50 Hz, Single Phase
Power Consumption	1500 W (approx.)
Construction	Mild Steel with Powder Coating
Insulation	High Density Ceramic Wool
Door	Front Opening with Double Pane Glass Window
Safety Features	Over Temperature Protection, Fan Failure Alarm, Door Open Alarm

APPLICATIONS

- Bitumen & Asphalt Testing Laboratories
- Quality Control in Road Construction
- Research & Development Laboratories
- Evaluation of Bituminous Materials
- Academic & Testing Institutes



PRECISE TEMPERATURE
CONTROL
10°C to 50°C



HUMIDITY
CONTROL
50% to 95% RH



CO₂ CONTROL
1% to 10%



UNIFORM AIR
CIRCULATION



DIGITAL PID
CONTROLLER



VERTICAL ROTATION
 15 ± 0.2 RPM



CALIBRATED AIR FLOW
4000 mL/min



SAFE & RELIABLE
OPERATION



ONE YEAR WARRANTY
Against manufacturing defects
(excluding consumables)



AMC & REPAIRING FACILITY AVAILABLE
Comprehensive Annual Maintenance Contract
and repairing services available

ARCHIT ASSOCIATES

📍 A-242, Jangir Puri,
Near Kali Mandir, Delhi-110 033 (INDIA)
☎ 91-8130973844

✉ architasst@gmail.com

GSTIN NO. 07A0UPG0300F2ZC

ARCHIT ASSOCIATES - YOUR PARTNER IN TESTING EXCELLENCE



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir,
Delhi-110 033 (INDIA)

91-8130973844, 9871406036, 9810052732, 9810052758

architasst@gmail.com

www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC

TRUSTED SUPPLIER OF
MATERIAL TESTING &
LABORATORY INSTRUMENTS



QUALITY
ASSURED
PRODUCTS



ACCURATE &
RELIABLE



LABORATORY &
SITE SOLUTIONS



CUSTOMIZED
TESTING
EQUIPMENT

ROTATIONAL VISCOMETERS (STANDARD TYPE)

- Widely used to measure the viscosity of fluids such as grease, paints, pharmaceuticals, coatings, adhesives, detergents, etc.
- For measuring the absolute viscosity of Newtonian liquids and the apparent viscosity of non-Newtonian liquids
- Measurement results show rotor number, speed, liquid viscosity, and full-scale percentage simultaneously
- 5-inch color touch screen, convenient and quick operation
- Automatic measurement mode, the rotor stops when measurements stabilize
- The temperature sensor is included for real-time sample temperature detection
- Software is included for data transfer

ROTORS (INCLUDED)



leveling bubble

rotor protection
frame (included)

temperature sensor
(included)

606

SPECIFICATION

Code	105	206	606
Range (mPa.s)	1-1×10 ⁶	1-2×10 ⁶	1-6×10 ⁶
Rotation speed (r/min)	6, 12, 30, 60	0.3, 0.6, 1.5, 3, 6, 12, 30, 60	0.1, 0.3, 0.6, 1.5, 3, 6, 12, 30, 60
Rotor No.	1#, 2#, 3#, 4#		
Resolution (mPa.s)	0.1		
Accuracy	±1.0%FS		
Repeatability	±0.5%FS		
Working environment	5°C-35°C, ≤80%		
Power supply	AC110V~220V, 50/60Hz		
Dimension (L×W×H)	460×180×420mm		
Weight	6.3kg		

STANDARD DELIVERY

Main unit	1 pc
Rotor	4 pcs (1#, 2#, 3#, 4#)
Temperature sensor	1 pc
Protection cap	1 pc
Main post for lifting	1 pc
Host connection handle	1 pc
Main engine base	1 pc
Rotor protection frame	1 pc
Power adapter	1 pc
Wrench	3 pcs

OPTIONAL DELIVERY

Printer	DRV-PRINTER
Thermostatic water bath	DRV-BATH
Rotor 0	DRV-SPINDLE



WIDELY USED IN



PHARMACEUTICALS



PAINTS & COATINGS



GREASE & LUBRICANTS



ADHESIVES



DETERGENTS



LABORATORIES &
RESEARCH

ARCHIT ASSOCIATES
A-242, Jahangir Puri,
Near Kali Mandir,
Delhi-110 033 (INDIA)



91-8130973844
9871406036
9810052732
9810052758



architasst@gmail.com

ALL PRODUCTS COMPLY WITH RELEVANT INDUSTRY STANDARDS



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844

E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC

SAFETY TODAY
SECURE TOMORROW



ALL SAFETY EQUIPMENT

COMPLETE PROTECTION. TOTAL SAFETY.

Your safety is our priority. We provide a wide range of high quality safety equipment to ensure protection in every work environment.



PREMIUM
QUALITY



RELIABLE
PROTECTION



COMFORTABLE
TO USE



COMPLIANT WITH
SAFETY STANDARDS

SAFETY HELMETS



- Strong & durable
- Adjustable & comfortable
- Available in multiple colors
- IS Certified

SAFETY HARNESS



- Full body protection
- Adjustable straps
- Shock absorbing lanyard
- IS 3521 Certified

SAFETY GOGGLES



- Impact resistant
- Anti-fog & anti-scratch
- Wide vision
- Comfortable fit

SAFETY JACKETS



- High visibility
- Reflective strips
- Lightweight & breathable
- Available in multiple sizes

WHY CHOOSE US?



Certified & Compliant

All products comply with national & international safety standards.



High Quality Products

Made from premium materials for maximum durability.



Comfort & Fit

Ergonomically designed for comfort during long hours of work.



Wide Range

One stop solution for all your safety needs.



Custom Branding

Logo printing & customized solutions available on request.

SAFETY SHOES



- Steel toe protection
- Anti-slip sole
- Oil & chemical resistant
- IS Certified

SAFETY GLOVES



- Cut, heat & impact resistant
- Excellent grip
- Comfortable & durable
- Suitable for all works

EAR PROTECTION



- Noise reduction
- Comfortable fit
- Reusable ear plugs
- High performance

RESPIRATORY PROTECTION



- Dust & particle protection
- Reusable & disposable options
- Comfortable & safe
- IS Certified

SAFETY BELTS



- Work positioning belt
- Adjustable & strong
- With hooks & lanyard
- IS 3521 Certified

FACE SHIELDS



- Full face protection
- Clear visibility
- Adjustable headgear
- Impact resistant

RAIN SUITS



- Waterproof material
- Durable & lightweight
- Comfortable fit
- Ideal for all weather

FALL ARRESTER



- Strong & durable
- Easy to install
- Shock absorbing
- IS 3521 Certified

INDUSTRIES WE SERVE

- Construction
- Manufacturing
- Oil & Gas
- Mining
- Power Plants
- Warehousing
- Laboratories
- And Many More...



WE ENSURE YOUR SAFETY, SO YOU CAN WORK WITH CONFIDENCE



REDUCES
WORKPLACE RISKS



ENHANCES
WORKER SAFETY



IMPROVES
PRODUCTIVITY



ENSURES LEGAL
COMPLIANCE



BUILT FOR EVERY
WORK ENVIRONMENT



SCAN TO VISIT
OUR WEBSITE

ARCHIT ASSOCIATES

A-242, Jahangir Puri, Delhi - 110033

Phone: 8130973844

Email: architasst@gmail.com

QUALITY SAFETY EQUIPMENT

FOR A SAFER WORKPLACE

Trusted by Professionals. Chosen for Protection.



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)
+91-8130973844, 9871406036 | architass@architass.com
www.architassociates.com
GSTIN NO. 07A0UPG0300F2ZC



SHOTCRETE PENETROMETER

PORTABLE DIGITAL PENETRATION RESISTANCE TESTER

FOR EARLY-AGE STRENGTH ASSESSMENT
OF SHOTCRETE

Shotcrete Penetrometer is a portable digital instrument used to determine the penetration resistance of freshly applied shotcrete. It provides a quick and reliable indication of early-age strength development and helps in deciding the safe time for removal of formwork, excavation or application of additional support.

FEATURES



Portable & Lightweight
Easy to carry and ideal for field testing



High Accuracy
Accurate measurement of penetration force



Digital Display
Large LCD with backlight for easy reading



Peak Hold Function
Captures maximum force during test



Data Storage & Transfer
Stores readings and transfers via USB



Interchangeable Needles
Quick change penetration needles



Rugged Carrying Case
Supplied in durable aluminium case



User Friendly Operation
Simple, reliable and maintenance free

TECHNICAL SPECIFICATIONS

Measuring Range	~ 1000 N
Resolution	0.1 N
Accuracy	± 0.5% Full Scale
Display	LCD Digital
Peak Hold	Yes
Data Storage	Yes
Data Output	USB
Power Supply	Rechargeable Battery
Operating Temperature	0 - 40°C
Carry Case	Aluminium Heavy Duty



COMPLETE WITH

- Digital Force Gauge
- Penetration Needle Holder
- Standard Penetration Needles (Set)
- Extension Rods - 2 Nos.
- USB Charging Cable
- Battery Charger
- Aluminium Carrying Case
- Instruction Manual



IDEAL FOR



TUNNEL
CONSTRUCTION



UNDERGROUND
MINING



METRO & RAIL
PROJECTS



ROCK
SUPPORT



SLOPE
STABILIZATION



HYDRO
PROJECTS

APPLICABLE STANDARDS

- EN 14488-2 (Method A)
- BS EN 14488-2
- ISO 14488-2

Test method for determining the penetration resistance of sprayed concrete (shotcrete) for estimation of early-age strength.

OPTIONAL ACCESSORIES

- Additional Needle Set
- Calibration Certificate
- Data Logging Software
- USB Data Cable
- Spare Battery



ACCURATE
RESULTS



RUGGED &
DURABLE



EASY TO
OPERATE



FIELD
PORTABLE



FAST & RELIABLE
TESTING



CALIBRATION
AVAILABLE



AFTER SALES
SUPPORT

ARCHIT ASSOCIATES - YOUR PARTNER IN TESTING EXCELLENCE

28
YEARS OF
EXPERIENCE



BUILT FOR LABS.
ENGINEERED FOR EXCELLENCE.
PRECISION. PERFORMANCE. TRUST.



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844 9810052732, 9810052758, 9871404034

E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



DIGITAL & ANALOG SIEVE SHAKER

ACCURATE • RELIABLE • EFFICIENT

Designed for precise particle size analysis in laboratories, quality control departments, R&D centers and educational institutions.

- Robust & Durable Construction
- Uniform Vibration for Accurate Results
- Low Noise & Maintenance Free
- Easy Operation & Safe Design
- Suitable for Dry Sieving
- Conforms to IS, ASTM Standards



DIGITAL MODELS



ANALOG MODELS

AVAILABLE SIZES

DIGITAL SIEVE SHAKER

20 CM DIA



Sieve Dia	200 mm
No. of Sieves	8 (2" Height)
Timer	0 - 30 Min
Power	230V, 50Hz

30 CM DIA



Sieve Dia	300 mm
No. of Sieves	8 (2" Height)
Timer	0 - 30 Min
Power	230V, 50Hz

45 CM DIA



Sieve Dia	450 mm
No. of Sieves	8 (2" Height)
Timer	0 - 30 Min
Power	230V, 50Hz

ANALOG SIEVE SHAKER

20 CM DIA



Sieve Dia	200 mm
No. of Sieves	8 (2" Height)
Timer	0 - 30 Min
Power	230V, 50Hz

30 CM DIA



Sieve Dia	300 mm
No. of Sieves	8 (2" Height)
Timer	0 - 30 Min
Power	230V, 50Hz

45 CM DIA



Sieve Dia	450 mm
No. of Sieves	8 (2" Height)
Timer	0 - 30 Min
Power	230V, 50Hz

KEY FEATURES



Strong & Sturdy Construction



Uniform Vibration for Consistent Sieving



Low Noise Operation



Digital Timer (0 - 30 Min)



Accommodates Standard Sieves



Safe & Easy Clamping System



Ideal for Laboratory & Industrial Use

IDEAL FOR

- Soil Testing Laboratories
- Cement & Concrete Testing
- Aggregate & Minerals Testing
- Quality Control Departments
- Educational & Research Institutes

ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844 9810052732, 9810052758, 9871404034

E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir,
Delhi-110 033 (INDIA)

☎ 91-8130973844 ✉ architassst@gmail.com

🌐 www.architassociates.com

GSTIN NO. 07A0UPG0300F22C

TRUSTED SUPPLIER OF
**MATERIAL TESTING &
LABORATORY INSTRUMENTS**



QUALITY
ASSURED
PRODUCTS



ACCURATE &
RELIABLE



LABORATORY &
SITE SOLUTIONS



CUSTOMIZED
TESTING
EQUIPMENT



TEST SIEVES

CONFORMING TO IS STANDARDS

FEATURES

- 🎯 High accuracy and precision performance
- 🛡️ Strong and durable construction
- 👍 Smooth finish and corrosion resistant
- ✅ Uniform aperture size as per IS Standards
- ✨ Easy to clean and maintain



BRASS SIEVE
20 cm DIA



GI SIEVE (POWDER COATED)
30 cm DIA



GI SIEVE (POWDER COATED)
45 cm DIA

TECHNICAL SPECIFICATIONS

Parameter	Brass Sieve 20 cm DIA	GI Sieve (Powder Coated) 30 cm DIA	GI Sieve (Powder Coated) 45 cm DIA
Type	Brass	GI (Powder Coated)	GI (Powder Coated)
Diameter	200 mm	300 mm	450 mm
Construction	Brass Frame with SS Wire Mesh	GI Frame with SS Wire Mesh	GI Frame with SS Wire Mesh
Finish	Polished Brass	Powder Coated	Powder Coated
Mesh Aperture	As per IS Standards	As per IS Standards	As per IS Standards
Conforming To	IS: 460 (Part 1) / IS: 460 (Part 2)		
Application	Particle Size Analysis in Laboratory, Quality Control & Material Testing		

WIDELY USED IN

- 🧪 CEMENT LABORATORIES
- 🏗️ CIVIL ENGINEERING LABORATORIES
- 🛡️ QUALITY CONTROL DEPARTMENTS
- 🔬 RESEARCH & INSTITUTES
- 🏗️ CONSTRUCTION INDUSTRY

MODELS AVAILABLE

BRASS SIEVE
20 cm DIA

GI SIEVE (POWDER COATED)
30 cm DIA

GI SIEVE (POWDER COATED)
45 cm DIA



QUALITY ASSURED
PRODUCTS



ACCURATE &
RELIABLE



CUSTOMIZED
TESTING EQUIPMENT



LABORATORY &
SITE SOLUTIONS



ARCHIT ASSOCIATES
A-242, Jahangir Puri,
Near Kali Mandir,
Delhi-110 033 (INDIA)



91-8130973844



architassst@gmail.com

ALL PRODUCTS COMPLY WITH RELEVANT INDUSTRY STANDARDS



ARCHIT ASSOCIATES

YOUR TRUSTED PARTNER IN LABORATORY & MATERIAL TESTING

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

+91-8130973844, 9871406036 | architassst@gmail.com

9810052732, 9810052758

www.architassociates.com



ACCURATE
RESULTS



RELIABLE
PERFORMANCE



QUALITY
ASSURED

TESTING EXCELLENCE IN EVERY RESULT

SILT CONTENT TEST APPARATUS (WITH SHAKER & WITHOUT SHAKER)

FOR DETERMINATION OF SILT CONTENT IN SOILS

Silt Content Test Apparatus is used to determine the percentage of silt content in soil by sedimentation method as per IS Standards.

The apparatus is available in two variants:

- With Shaker
- Without Shaker

PRINCIPLE

Soil sample is dispersed in water and transferred to a measuring cylinder. After thorough shaking (where applicable), the suspension is allowed to settle. The volume of silt in suspension at a specified time is measured to determine the silt content.

APPLICATION

- Determination of silt percentage in soil samples
- Useful in soil classification and engineering investigations
- Applicable for laboratory and field soil testing



WITHOUT SHAKER



WITH SHAKER

KEY FEATURES

- Available with shaker or without shaker
- Transparent glass bottle with rubber stopper and delivery tube
- Three measuring cylinders (1000 ml capacity) with scale
- Glass rod, thermometer & other accessories
- Sturdy wooden carrying case for safe storage & portability
- Easy to assemble, operate & clean
- Precision glassware for accurate results

TECHNICAL SPECIFICATIONS

Item	Specification
Bottle Capacity	4 Litres (Glass)
Measuring Cylinders	1000 ml (3 Nos.)
Temperature Range	0°C to 50°C
Shaker Travel	25 ± 1 mm (Only with Shaker)
Speed	120 ± 10 Strokes/min (Only with Shaker)
Dimension (L x W x H)	
Without Shaker	450 x 300 x 120 mm (Approx.)
With Shaker	450 x 320 x 280 mm (Approx.)
Net Weight	
Without Shaker	7 kg (Approx.)
With Shaker	15 kg (Approx.)

STANDARDS

IS: 2720 (Part - 4)
IS: 2720 (Part - 40)
BS 1377: Part 2
ASTM D 422

SCOPE OF SUPPLY (BOTH MODELS)

- Glass bottle 4 Litres with rubber stopper and delivery tubes - 1 No.
- Measuring cylinders 1000 ml - 3 Nos.
- Glass rod - 1 No.
- Rubber tubing - 1 No.
- Wooden carrying case - 1 No.
- Instruction manual - 1 No.

WITHOUT SHAKER MODEL

Designed for laboratories where manual shaking is preferred. Simple, compact and reliable for silt content determination.

Ideal For:

- Routine soil testing laboratories
- Field & site investigations
- Educational institutions

WITH SHAKER MODEL

Equipped with a constant speed shaker for uniform dispersion of soil sample ensuring repeatable and accurate results.

Ideal For:

- Advanced laboratories
- Research & quality control
- High volume testing environments



PRECISE RESULTS
Consistent & accurate
testing performance



SUPERIOR QUALITY
High grade components
for reliability



EASY TO OPERATE
User friendly design
& simple controls



TEST CERTIFICATES
Supplied with
calibration certificates



AFTER SALES SUPPORT
Prompt service &
technical support



SCAN TO VISIT
OUR WEBSITE

ARCHIT ASSOCIATES

A-242, Jahangir Puri, Delhi - 110033 (INDIA)

+91-8130973844, 9871406036

architassst@gmail.com

28
YEARS OF
EXPERIENCE



BUILT FOR LAST.
ENGINEERED FOR EXCELLENCE.
Precision. Performance. Perfection.



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844

E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



SKID RESISTANCE TESTER

For Measurement of Skid Resistance
of Road Surfaces



ACCURATE
RESULTS



ROBUST &
DURABLE



EASY TO
OPERATE



IDEAL FOR
LABORATORIES
& FIELD USE



CONFORMS TO
INTERNATIONAL
STANDARDS



DESCRIPTION

- The apparatus is mounted on a sturdy base having three adjustable feet so that the swing of the pendulum during the test moves through a truly circular path while operated directly on a road surface.
- The apparatus consists of a spring loaded rubber slider mounted on a pendulum arm.
- The tester measures the energy loss when a rubber slider edge is propelled over the surface under test.
- The slider lifting device is incorporated in the pendulum base assuring accurate adjustment operations. The height adjusting system is simple and reliable. The pointer, made from light alloy, has extremely low frictions granting high precision results.
- It operates a pendulum rotating about a spindle which is attached to a vertical pillar. At the end of the tubular arm a head of known mass is fitted with a rubber slider.
- The pendulum is released from a horizontal position so that it strikes the sample surface at a constant speed. The distance travelled by the head after hitting the sample is determined by the friction of the sample surface.
- The release mechanism of the pendulum arm has an original solution reducing the friction to minimum for better accuracy.
- This provides road engineers with a method of checking the resistance of wet and dry surfaces to slipping and skidding, both in the lab and on site.
- Lab testing kit for Stone Polishing Machine samples Optional extra.

TECHNICAL SPECIFICATIONS

Type	Pendulum Type
Test Method	British Pendulum Method
Scale Range	0 - 150 BPN (British Pendulum Number)
Least Count	0.5 BPN
Slider	Standard Rubber Slider
Pendulum Length	410 mm
Weight	Approx. 22 - 25 kg
Dimensions (L x W x H)	Approx. 700 x 300 x 600 mm
Accessories	Rubber Slider, Wrench, Cleaning Brush, User Manual

FEATURES

- Reliable and repeatable results
- Easy to use & maintain
- Robust construction for long life
- Smooth movement of pendulum
- Suitable for field and laboratory use

APPLICATION

- Highway & Road Construction
- Quality Control Laboratories
- Research & Development
- Airport Runways
- Parking Areas & Industrial Floors
- Road Maintenance Authorities

TEST PROCEDURE



STANDARDS



Conforms to Relevant
ASTM / BS / IS Standards

IDEAL FOR



Field Testing



Quality Control
Labs



Research &
Development



Construction
Projects

BENEFITS

- Improves Road Safety
- Ensures Quality Surfaces
- Simple, Fast & Accurate Testing

ARCHIT ASSOCIATES

A-242, Jahangir Puri,
Near Kali Mandir, Delhi-110 033 (INDIA)

91-8130973844

architasst@gmail.com

GSTIN NO. 07A0UPG0300F2ZC

**ARCHIT ASSOCIATES**A-242, Jahangir Puri, Near Kali Mandir,
Delhi-1100 033 (INDIA)

+91-8130973844

architasst@gmail.com Web www.architassociates.com

**TESTING EXCELLENCE
IN EVERY RESULT***Your Trusted Partner in Laboratory
& Material Testing*

With 28 years of rich experience, ARCHIT ASSOCIATES has been at the forefront of supplying high quality Laboratory and Material Testing Equipment to a wide range of industries, research institutes, quality control labs and educational institutions.

28 Years of
Industry
ExperienceWide Range of
Quality
ProductsExpert Advice &
After Sales
SupportPrecision,
Reliability
& Long Life

SLAKE / LIME REACTIVITY TEST APPARATUS

PRECISION TESTING. RELIABLE RESULTS.

Engineered for accurate determination of lime reactivity in water under controlled conditions. Built for laboratories that demand consistency, reliability and performance.

**PRECISE TEMPERATURE CONTROL**Upto 100°C with accuracy $\pm 0.01^\circ\text{C}$ **UNIFORM STIRRING**

Stirrer with Geared Motor for consistent results

**PREMIUM QUALITY**

Stainless Steel (SS 304), Insulated & durable construction

**EASY OPERATION**

User friendly controls for smooth testing

**TRUSTED PERFORMANCE**

Designed for laboratories that value accuracy

**SALIENT FEATURES**

- Stirrer with Geared Motor for uniform mixing
- Stainless Steel Container duly insulated
- Teflon Sealing Cap to minimize heat loss from the vessel
- Digital Temperature Indicator with Stirring Control Knob
- Stainless Steel Temperature Sensor
- Mounted on heavy duty stand for easy operation
- Temperature range upto 100°C with accuracy $\pm 0.01^\circ\text{C}$

TECHNICAL SPECIFICATIONS

- Temperature Range : 1 No. Stirrer with Geared Motor
- Temperature Accuracy : 1 No. Stainless Steel Container (Insulated)
- Stirring Speed : 1 No. Teflon Sealing Cap
- Container Material : 1 No. Digital Temperature Indicator with Stirring Control
- Insulation : 1 No. Stainless Steel Temperature Sensor
- Power Supply : 1 No. Heavy Duty Stand

BRASS TEST SIEVES

Conforms to IS 460 (Part 1) / ISO 3310-1

SIEVE SIZE (mm)	SIEVE SIZE (µm)
4.75 mm	4750 µm
2.36 mm	2360 µm
1.18 mm	1180 µm
600 µm	600 µm
300 µm	300 µm
150 µm	150 µm
75 µm	75 µm

Also available: Full set of Brass Sieves with Lid & Pan.

ACCESSORIES SUPPLIED**OTHER LIME TESTING ACCESSORIES****28 YEARS OF EXPERIENCE. THOUSANDS OF SATISFIED LABS. COUNTLESS SUCCESS STORIES.**R&D LABS &
UNIVERSITIESCEMENT &
CONSTRUCTION
INDUSTRYQUALITY CONTROL
LABORATORIESGOVERNMENT & PRIVATE
TESTING LABSRESEARCH & DEVELOPMENT
INSTITUTIONSCALL US TODAY
FOR EXPERT
SOLUTIONS9810052732
9810052758
9871406036**ARCHIT ASSOCIATES**A-242, Jahangir Puri, Near Kali Mandir,
Delhi-1100 033 (INDIA)
+91-8130973844
architasst@gmail.com**BUILT FOR LABS.
ENGINEERED FOR EXCELLENCE.**
Precision. Performance. Perfection.



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844

E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC

TESTING
EXCELLENCE IN
EVERY RESULT

SLUMP TEST APPARATUS

ACCURATE • RELIABLE • EASY TO USE

Used to determine the workability of concrete by measuring its slump in accordance with IS / ASTM Standards.

FEATURES & BENEFITS



Robust Construction
Heavy duty steel construction with rust proof coating.



Dimensional Accuracy
Precision engineered for accurate and repeatable results.



Easy to Use
Simple design for quick assembly and smooth operation.



Standard Compliance
Conforms to IS: 1199, ASTM C143 / C143M and other international standards.



Long Lasting
Durable, corrosion resistant and built for site & lab use.



Customized Solutions
We offer customized Slump Test Apparatus as per client requirements.



CUSTOMIZED SLUMP TEST APPARATUS

- ✓ Available in non-standard sizes as per project requirement
- ✓ Special molds for precast elements
- ✓ Customized base plate, handles & height
- ✓ Logo branding available
- ✓ Bulk order for projects, institutes & labs

**YOUR REQUIREMENT,
OUR SOLUTION.**



AVAILABLE MODELS (STANDARD)

SLUMP CONE - 100 mm



Top Dia.: 100 mm
Bottom Dia.: 200 mm
Height: 300 mm

SLUMP CONE - 150 mm



Top Dia.: 150 mm
Bottom Dia.: 300 mm
Height: 450 mm

SLUMP CONE - 200 mm



Top Dia.: 200 mm
Bottom Dia.: 400 mm
Height: 450 mm

SLUMP CONE - 250 mm



Top Dia.: 250 mm
Bottom Dia.: 500 mm
Height: 750 mm

APPLICATION

- ✓ Ready Mix Concrete
- ✓ Construction Sites
- ✓ Laboratory Testing
- ✓ Quality Control
- ✓ R&D and Educational Institutes



CONFORMS TO RELEVANT IS / ASTM STANDARDS

- IS 1199
- ASTM C143 / C143M
- EN 12350-2



ACCURATE
TESTING



DURABLE
CONSTRUCTION



RELIABLE
RESULTS



EASY TO
OPERATE



SCAN TO VISIT

ARCHIT ASSOCIATES

A-242, Jahangir Puri, Delhi - 110033

Phone: 8130973844

Email: architasst@gmail.com

QUALITY INSTRUMENTS

ACCURATE TESTING | RELIABLE RESULTS

BUILT FOR LABS, ENGINEERED FOR EXCELLENCE

TRUST FOR LONG-TERM ENGINEERING SOLUTIONS



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844 E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F22C



SPECIFIC GRAVITY BALANCE (BUOYANCY BALANCE)

Designed for accurate determination of Specific Gravity of materials using buoyancy method. Ideal for testing of Cube Moulds, Marshall Moulds and other industrial components. Equipped with two digital weighing balances and a top pan balance for high precision and reliable results.

FEATURES



High Accuracy

NABL calibrated digital balances ensure precise and reliable measurements.



Complete Set

Includes two digital weighing balances & one top pan balance for versatile use.



Heavy Duty Construction

Robust stand with lifting arrangement for easy handling of water tank.



Wide Application

Suitable for Specific Gravity testing of Cube Moulds, Marshall Moulds and other components.



Durable & Reliable

Built with quality materials for long service life and consistent performance.



COMPLETE SET INCLUDES

- ✓ NABL Calibrated Digital Weighing Balance
Cap.: 6 Kg x 0.1 g (Qty. 2 Nos.)
- ✓ Digital Top Pan Balance
Cap.: 15 Kg x 0.5 g (Qty. 1 No.)
- ✓ Bottom Hook for hanging arrangement
- ✓ Heavy Duty Stand with lifting arrangement
for water tank
- ✓ Water Tank (Stainless Steel)
- ✓ Hanger for specimen



TECHNICAL SPECIFICATIONS

Use	Determination of Specific Gravity by Buoyancy Method
Suitable For	Cube Moulds, Marshall Moulds & other components
Digital Weighing Balance (Qty. 2 Nos.)	Capacity: 6 Kg x 0.1 g (NABL Calibrated)
Digital Top Pan Balance (Qty. 1 No.)	Capacity: 15 Kg x 0.5 g (NABL Calibrated)
Weighing Method	Buoyancy (Water Displacement)
Construction	Heavy duty steel frame with corrosion resistant coating
Tank Material	Stainless Steel
Tank Lifting	Manual lifting arrangement with chain & crank
Accessories	Bottom hook, Hanger, Water tank
Power Supply (Balances)	230V AC, 50Hz

IDEAL FOR

- Civil Engineering Laboratories
- Quality Control Laboratories
- Road & Highway Testing
- Concrete & Construction Industry
- Research & Educational Institutes



CUBE MOULD



MARSHALL MOULD



ACCURATE
RESULTS



NABL
CALIBRATED



STRONG & DURABLE
CONSTRUCTION



EASY TO
OPERATE



RELIABLE &
LONG LASTING

ARCHIT ASSOCIATES – YOUR PARTNER IN TESTING EXCELLENCE





ARCHIT ASSOCIATES

YOUR TRUSTED PARTNER IN LABORATORY & MATERIAL TESTING

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

91-8130973844

architassst@gmail.com

www.architassociates.com

GSTIN NO. 07A0UPG0300F22C



ACCURATE
RESULTS



RELIABLE
PERFORMANCE



QUALITY
ASSURED

TESTING EXCELLENCE IN EVERY RESULT

TENSILE TESTING MACHINES

PRECISION • STRENGTH • RELIABILITY

Our Tensile Testing Machines are designed to deliver accurate and reliable results for a wide range of materials. Engineered for performance, built for durability.

KEY FEATURES

- High Precision & Accuracy
- Digital Display / Touch Screen Options
- Robust Construction for Long Life
- User Friendly & Easy Operation
- Over Load & Over Travel Protection
- Compliance with International Standards

FLOOR / TABLE TOP MODELS

FLOOR MODEL SINGLE LOAD CELL

Capacity : 20 Kg to 250 Kg



For Cable & Conductor

TABLE TOP MODEL WITH TOUCH SCREEN

Capacity : 20 Kg to 250 Kg



For PP & HDPE Ropes

3 KG TO 10 KG TENSILE TESTER

Capacity : 20 Kg to 250 Kg



For FIBC Belts

FLOOR MODEL DOUBLE LOAD CELL

Capacity : 50 Kg & 5000 Kg



TABLE TOP MODEL DOUBLE LOAD CELL

Capacity : Upto 250 Kg



TABLE TOP MODEL DOUBLE LOAD CELL

Capacity : Upto 500 Kg



ADVANCED & SPECIALIZED MODELS

SINGLE SCREW DOUBLE LOAD CELL

Capacity : 50 Kg & 5000 Kg



For FIBC Belt & Tape

SINGLE SCREW DOUBLE LOAD CELL

Capacity : 500 Kg & 5000 Kg



For FIBC Belt & Jumbo Bag Fabric

UTM



SINGLE LOAD CELL DOUBLE SCREW

Capacity : 2000 Kg to 10,000 Kg



For Cable & Conductor

3 LOAD CELL TENSILE TESTING MACHINE

Capacity : 50 Kg / 500Kg / 5000 Kg



For Tape, Fabric & Belt

DOUBLE LOAD CELL DUAL COLUMN TENSILE TESTER

Capacity : 500 Kg & 5000Kg



For FIBC Belt & Jumbo Bag Fabric

SPECIFICATIONS

BASIC DIGITAL TENSILE TESTING MACHINE

- 4 Line LCD Display unit / 4" HMI Touch Screen which shows Load, Elongation, Maximum Elongation & Elongation Percentage (%)
- Imported Guide Rods with Linear Motion Bearings for Vibration less & Frictionless Travel
- Over Load safety
- Over Travel safety
- Auto Reverse
- Auto Stop after completion of test
- Over Load safety
- Over Travel safety
- Speed Variation through AC Drive & Motor or Stopper Motor & Drive
- Unit Selection KG, Newton, KN
- Last 5 samples readings will be saved in controller

COMPUTERISED TENSILE TESTING MACHINE

- 4 Line LCD Display unit / 4" HMI Touch Screen which shows Load, Elongation, Maximum Elongation & Elongation Percentage (%)
 - Controller can save Data of last 5 samples
 - Speed Variation through AC Drive & Motor or Stopper Motor & Drive
 - Accuracy within tolerance limit of $\pm 1\%$
 - Force Measurement in KG, Newton, KN & MPa
 - Elongation through Potentiometer
 - Auto Reverse Facility
 - Auto Stop Facility
 - Over Load & Over Travel safety
 - Computer connection through USB
 - Speed Control through Software Program
 - Software calculates the results as per below:
- Maximum Load
 - Elongation at Maximum Load
 - Elongation Percentage at Maximum Load
 - Tensile Strength at Maximum Load
 - Break Load
 - Elongation at Break Load
 - Elongation Percentage at Break Load
 - Tensile Strength at Break Load
 - Yield Load
 - Elongation at Yield Load
 - Elongation Percentage at Yield Load
 - Tensile Strength at Yield Load
 - Modulus of Elasticity
 - Stiffness of Maximum Load (GRC)



PRECISE RESULTS
Consistent & accurate
testing performance



SUPERIOR QUALITY
High grade components
for long life



EASY OPERATION
Intuitive controls with
advanced technology



TEST CERTIFICATES
Supplied with calibration
certificate



AFTER SALES SUPPORT
Prompt service &
technical support



SCAN TO VISIT
OUR WEBSITE

ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near - 110033 (INDIA) 91-8130973844

architassst@gmail.com

28
YEARS OF
EXPERIENCE



BUILT FOR LABS.
ENGINEERED FOR EXCELLENCE.
Precision. Performance. Perfection.



ARCHIT ASSOCIATES

YOUR TRUSTED PARTNER IN LABORATORY & MATERIAL TESTING

📍 A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

☎ +91-8130973844 9871406036

✉ architasst@gmail.com

🌐 www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



SCAN TO VIEW
PRODUCT RANGE

TILE ABRASION TESTING MACHINE

REF. STD. IS 1237, 1706

The tile sample is pressed under a specified load against a grinding path, strewn evenly with an abrasive powder, revolving at the rate of 30 ± 1 rpm. At the end of 100 revolutions of the disc, the second parallel side of the tile is subjected to wear an equal number of revolutions of the disc. The wear of the tile is measured using thickness measuring device (which is supplied at extra cost).

The Abrasion Testing Machine consists of a replaceable grinding path, fitted on to a disc rotating about its vertical axis. The rotating disc and the grinding path are enclosed by a circular tray. A bracket holds the specimen. The loading is by a counterbalanced lever. A funnel is provided to facilitate charging the grinding path with abrasive powder.

Supplied complete with automatic pre-set revolution counter and one 7.5 kg weight, but without abrasive powder and thickness measuring device.



SALIENT FEATURES

- 🔍 Conforms to IS: 1237, 1706
- 🔧 Robust & Precision Engineering
- ⚖️ Counterbalanced Loading System
- 🔄 Rotating Disc with Replaceable Grinding Path
- 📊 Automatic Pre-set Revolution Counter
- 🔑 Funnel for Easy Charging of Abrasive Powder
- 🏗️ Sturdy Construction for Long Life
- 🛠️ Easy Operation & Maintenance

STANDARD ACCESSORIES



Replaceable
Grinding Path
(Disc)



Funnel



7.5 kg
Weight



Pre-set
Revolution
Counter



Tools &
Accessories



Instruction
Manual

TECHNICAL SPECIFICATIONS

Standards	IS: 1237, 1706
Grinding Disc Diameter	200 mm (Approx.)
Grinding Path Width	12.5 mm (Replaceable)
Grinding Path Material	Cast Iron
Speed of Disc	30 ± 1 rpm
Total Revolutions	100 Rev. on each of the two parallel sides
Loading	Counterbalanced Lever System
Standard Weight Supplied	7.5 kg
Motor	1/4 HP, Single Phase
Power Supply	230 Volts, 50 Hz, Single Phase, AC
Counter	Pre-set, 5 Digit, Auto Stop
Overall Dimensions (L x W x H)	800 x 550 x 750 mm (Approx.)
Net Weight (Approx.)	70 - 80 Kgs

OPERATING PRINCIPLE

- 1 The tile specimen is pressed under a specified load against a grinding path evenly strewn with abrasive powder, revolving at 30 ± 1 rpm.
- 2 After 100 revolutions, the second parallel side of the tile is subjected to equal number of revolutions.
- 3 The abrasive action wears the surface of the tile.
- 4 The wear of the tile is measured using a thickness measuring device (Optional).

APPLICATIONS

- 🏠 Ceramic & Vitrified Tile Industry
- 🔍 Quality Control Laboratories
- 🏢 Research & Development Centres
- 🏛️ Government Testing Laboratories
- 🎓 Educational Institutions
- 🏗️ Building Material Testing

OPTIONAL ACCESSORIES

- Thickness Measuring Device (Dial Gauge Type)
- Abrasive Powder (Carborundum / Silicon Carbide)
- Additional Weights
- Calibration Certificate



ACCURATE
RESULTS



RUGGED & DURABLE
CONSTRUCTION



PRECISE
PERFORMANCE



AUTO STOP
COUNTER



EASY TO OPERATE
& MAINTAIN



RELIABLE &
LONG LASTING

SCAN TO VISIT
OUR WEBSITE



ARCHIT ASSOCIATES

📍 A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

☎ +91-8130973844 | ✉ architasst@gmail.com | www.architassociates.com

28+
YEARS OF
EXPERIENCE



BUILT FOR LAB.
ENGINEERED FOR EXCELLENCE
Precision. Performance. Trust.

ARCHIT ASSOCIATES – YOUR PARTNER IN TESTING EXCELLENCE



ARCHIT ASSOCIATES

YOUR TRUSTED PARTNER IN LABORATORY & MATERIAL TESTING

📍 A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

☎ +91-8130973844, 9871406036 ✉ architasst@gmail.com

🌐 www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



ACCURATE
RESULTS



RELIABLE
PERFORMANCE



QUALITY
ASSURED

TESTING EXCELLENCE IN EVERY RESULT

TILE AUTOCLAVE

REF. STD.: IS 15477:2019 & ISO 10545-3

Tile Autoclave is used for determination of the effect of autoclave treatment on ceramic tiles. The unit subjects tiles to steam pressure and high temperature in saturated condition for a specified period. It is widely used in tile manufacturers, ceramic laboratories and quality control departments.

FEATURES

- ⚙️ **Complies with IS 15477:2019 & ISO 10545-3**
For accurate and reliable testing
- ⚙️ **Robust Construction**
Heavy duty chamber with durable finish
- 📊 **Digital Temperature Controller**
Accurate temperature control & indication
- 💧 **High Pressure & High Temperature**
Designed for saturated steam condition
- 💧 **Safety Pressure Relief Valve**
Ensures safe operation
- 👤 **Easy to Operate**
User friendly control panel
- 💧 **Drain Valve**
For easy water discharge
- 👤 **Three Phase Supply with MC**
Reliable electrical control with protection

TECHNICAL SPECIFICATIONS

Standard	IS 15477:2019, ISO 10545-3
Chamber Material	Stainless Steel
Chamber Construction	Heavy Duty
Chamber Capacity	Approx. 25 Litres
Working Pressure	2.0 ± 0.1 kg/cm ² (Saturated Steam)
Working Temperature	121 ± 2 °C
Temperature Controller	Digital PID Controller
Temperature Display	Digital
Pressure Gauge	Bourdon Type, 0 - 4 kg/cm ²
Heating Element	U Type Immersion Heater
Power Rating	3 MW (Approx.)
Electrical Supply	415V, 50 Hz, Three Phase, AC
Control	Start / Stop Push Button
Safety Device	Pressure Relief Valve
Drain	Drain Valve with Handle
Insulation	High Grade Mineral Wool
Overall Dimensions (L x W x H)	550 x 550 x 1150 mm (Approx.)
Net Weight (Approx.)	70 - 80 kg (Approx.)



STANDARD ACCESSORIES

- Pressure Gauge - 1 No.
- Safety Valve - 1 No.
- Specimen Holder - 1 No.
- Water Level Indicator - 1 No.
- Instruction Manual - 1 No.



APPLICATIONS



TILE
MANUFACTURERS



CERAMIC
LABORATORIES



QUALITY CONTROL
DEPARTMENTS



RESEARCH &
DEVELOPMENT



CONSTRUCTION
MATERIAL TESTING

OPTIONAL ACCESSORIES

- Additional Sample Rack
- Digital Timer
- Water Softener
- Automatic Water Level Controller
- Temperature Data Logger
- Calibration Certificate
- Spare Heating Element
- Insulated Outer Jacket



ACCURATE RESULTS
Consistent & reliable
test performance



RUGGED & DURABLE
Strong construction
for long life



EASY TO OPERATE
User friendly design
& controls



TEST CERTIFICATES
Supplied with
calibration certificates



AFTER SALES SUPPORT
Prompt service &
technical support



SCAN TO VISIT
OUR WEBSITE

ARCHIT ASSOCIATES

📍 A-242, Jahangir Puri, Delhi - 110033 (INDIA)

☎ +91-8130973844, +91 9871406036

✉ architasst@gmail.com

28
YEARS OF
EXPERIENCE



**BUILT FOR LABS.
ENGINEERED FOR EXCELLENCE.**
Precision. Performance. Trust.



ARCHIT ASSOCIATES

YOUR TRUSTED PARTNER IN LABORATORY & MATERIAL TESTING

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)
+91-8130973844 9871406036 architassst@gmail.com
www.architassociates.com
GSTIN NO. 07A0UPG0300F2ZC



SCAN TO VIEW
PRODUCT RANGE

LABORATORY TILE / CRAZING AUTOCLAVE

REF. STD. IS 13630 (PART 9)

The Laboratory Tile / Craze Autoclave is used for determining the resistance of ceramic and vitrified tiles against crazing under high-pressure steam conditions as specified in IS 13630 (Part 9).

The equipment is manufactured with a heavy-duty quadruple wall insulated pressure chamber to ensure excellent thermal efficiency, operator safety and long service life. The digitally controlled system provides precise temperature regulation while multiple safety mechanisms ensure reliable operation under high-pressure conditions.

SALIENT FEATURES

- Triple safety against accidental over pressure
Safety Valve / Pressure Control Switch /
Digital Temp. Indicator cum controller
- Simple in design and operation
- Easy heater replacement
- Calibrated Pressure Gauge
- Teflon coated boiling chamber
resistive against corrosion
- Hard chrome plated autoclave lid
resistive against corrosion
- Safety against Electrical Short circuit
- Digital Auto Tuning PID Temp. Controller
with Independent Soak Timer
- Cooling fan for faster cooling
- Sturdy construction with powder
coated finish



STANDARD ACCESSORIES



TECHNICAL SPECIFICATIONS

PARAMETER	SPECIFICATION
Standard	IS 13630 (Part 9)
Chamber Size	330 x 450 mm
Maximum Tile Size	12" x 16"
Chamber Volume	43 Liters
Maximum Pressure	7.5 kg/cm ² (735 kPa)
Chamber Construction	Heavy Duty Mild Steel with Teflon Coating
Wall Construction	Quadruple Wall
Insulation	Air Gap + Glass Wool
Autoclave Lid	Hard Chrome Plated
Pressure Gauge	Calibrated
Temperature Controller	Digital Auto Tuning PID
Timer	Independent Soak Timer
Safety Valve	Stainless Steel
Cooling	Forced Air Cooling Fan
Overall Dimensions (L x W x H)	800 x 1000 x 1070 mm
Electrical Supply	440V, 50Hz, Three Phase AC

PRINCIPLE OF OPERATIONS

- The tile specimens are placed inside the Teflon-coated pressure chamber and the lid is securely locked.
- Water is heated to generate saturated steam under controlled pressure.
- The Digital PID Controller accurately maintains the required temperature while the pressure regulator controls steam pressure automatically.
- Upon completion of the preset soak cycle, the heater switches OFF automatically.
- Pressure is safely released through the drain valve before opening the chamber. Multiple independent safety devices protect the operator against excessive pressure and temperature.

APPLICATIONS

- Ceramic Tile Manufacturers
- Quality Control Laboratories
- Research & Development Centres
- Government Testing Laboratories
- Building Material Testing Laboratories
- Educational Institutes

AVAILABLE MODELS

- 43 LITERS (STANDARD)
- CUSTOM CHAMBER SIZE
- OPTIONAL ACCESSORIES AVAILABLE
- BUILT TO MEET YOUR REQUIREMENT

OPTIONAL ACCESSORIES



ACCURATE
TEMPERATURE &
PRESSURE CONTROL

QUADRUPLE WALL
HEAVY DUTY
CONSTRUCTION

TRIPLE SAFETY
PROTECTION

RELIABLE
PERFORMANCE

LONG SERVICE
LIFE

EASY
MAINTENANCE

INDUSTRIAL GRADE
COMPONENTS

CALIBRATION
SUPPORT

PROMPT AFTER
SALES SERVICE



SCAN TO VISIT
OUR WEB SITE

ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)
+91-8130973844 | +91 9871406036 | architassst@gmail.com
www.architassociates.com

28+
YEARS OF
EXPERIENCE



BUILT FOR LABS.
ENGINEERED FOR EXCELLENCE.
Precision, Performance, Trust.



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844

E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



CERAMIC TILE FLATNESS & STRAIGHTNESS TESTING INSTRUMENT

Ref Std: IS 13630 Part 1 & ISO 10545-2

This equipment is designed to determine the flatness and straightness of ceramic wall tiles, in accordance with IS 13630 Part 1 and ISO 10545-2 standards. Testing is carried out by using a standard reference board made of high-grade aluminum alloy, customized to match the dimensions of the tile being tested. The apparatus provides a quick and reliable way to assess dimensional conformity of tiles during production and quality inspection.

FEATURES

- High-accuracy measurement with an error margin of ± 0.1 mm.
- Customizable aluminum alloy standard boards to match various tile sizes.
- Wide testing range for tiles from 60×60 mm up to 1600×1600 mm.
- Lightweight and durable construction for ease of handling and long operational life.
- Simple and efficient operation, suitable for both laboratory and production environments.

TECHNICAL SPECIFICATION

Parameter	Specification
Accuracy	± 0.1 mm
Suitable Tile Sizes	60×60 mm to 600×600 mm 60×60 mm to 800×800 mm 60×60 mm to 1000×1000 mm 60×60 mm to 1600×1600 mm
Standard Board Material	Aluminium Alloy (customized to sample size)
Measurement Capability	Flatness and Straightness of ceramic Tiles

APPLICABLE FOR

- Ceramic Wall & Floor Tiles
- Tile Manufacturing Units
- Quality Control Laboratories
- Inspection & Testing Centers
- Research & Development

STANDARD ACCESSORIES

- Aluminum Alloy Reference Board
- Dial Indicators (0-10 mm)
- Adjustable Measuring Carriage
- Adjustable Stops
- Carrying Case
- Instruction Manual

KEY BENEFITS

- ✓ Ensures compliance with IS 13630 Part 1 & ISO 10545-2 standards
- ✓ Quick and reliable dimensional conformity assessment
- ✓ Custom-built reference boards for different tile sizes
- ✓ Durable, precise and easy to use
- ✓ Essential instrument for maintaining product quality

AVAILABLE STANDARD SIZES (REFERENCE BOARD)

Reference Board Size (Approx.)

600×600 mm

800×800 mm

1000×1000 mm

1600×1600 mm



NOTE:

Other customized reference board sizes are available on request.



HIGH ACCURACY
MEASUREMENT



DURABLE &
RELIABLE



EASY TO USE &
MAINTAIN



IDEAL FOR LAB &
PRODUCTION



QUALITY CONTROL
ESSENTIAL

ARCHIT ASSOCIATES - YOUR PARTNER IN TESTING EXCELLENCE



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844

E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F22C

TESTING TODAY
BUILDING TOMORROW



ACCURATE



RELIABLE



DURABLE



TRUSTED

DIGITAL MOR TESTING MACHINE (MODULUS OF RUPTURE)

ACCURATE • PRECISE • RELIABLE

Used for determination of Modulus of Rupture (Flexural Strength) of cement concrete beams and other building materials. Microprocessor based digital unit with high precision load measurement and variable speed control.

KEY SPECIFICATIONS

	Suitable Sample Size	50 mm x 100 mm (Other sizes also available)
	Capacity (Load Range)	10 kg/cm ² to 750 kg/cm ²
	Load Accuracy	± 1% of Full Scale
	Speed Range	5 mm to 50 mm per minute (Variable Speed)
	Load Display	Digital (kgf / kN / kg/cm ²)
	Control System	Microprocessor based with Digital Display
	Power Supply	230V AC, 50Hz, Single Phase
	Standards	IS 515, IS 8781, IS 1786, ASTM C78 EN 12390-5 & Other International Standards



SCAN TO VISIT
OUR WEBSITE



VARIABLE SPEED
5 to 50 mm/min



HIGH ACCURACY
Precise & Reliable
Results



RUGGED & DURABLE
Construction



**MICROPROCESSOR
BASED CONTROL**



DIGITAL DISPLAY
Easy to Read

AVAILABLE SAMPLE SIZES (STANDARD MODELS)

50 mm x 100 mm
(STANDARD)



Capacity Range
10 to 750 kg/cm²

75 mm x 150 mm



Capacity Range
30 to 750 kg/cm²

100 mm x 100 mm



Capacity Range
10 to 750 kg/cm²

100 mm x 150 mm



Capacity Range
10 to 750 kg/cm²

150 mm x 150 mm



Capacity Range
10 to 750 kg/cm²

CUSTOM SIZES
(As per Requirement)



We also provide customized
fixtures & machines as per
client requirement.

* Other capacities and sample sizes available on request. * Special fixtures available for Block, Paver, Kerb Stone, RCC & Other Products.

SALIENT FEATURES

- ✓ Digital load indicator with high resolution
- ✓ Variable speed drive for uniform loading
- ✓ Automatic peak load hold & exact MOR calculation
- ✓ Strong, rigid & vibration free frame structure
- ✓ Overload protection & safety limit switch
- ✓ Compact, user friendly & low maintenance
- ✓ Suitable for laboratory & site quality control



PRECISE LOADING
SYSTEM



MICROPROCESSOR
CONTROL UNIT



ACCURATE MOR
RESULTS

IDEAL FOR TESTING OF

- Concrete Beams
- Precast Elements
- Kerb Stones
- Paving Blocks
- Fiber Reinforced Concrete
- Other Building Materials

CUSTOMIZATION AVAILABLE



Any Sample Size
(As per Requirement)



Customized Capacity
(As per Requirement)



Special Fixtures &
Accessories



Logo Branding
Available



Turnkey Lab Setup
Available



SCAN TO VISIT
OUR WEBSITE

ARCHIT ASSOCIATES

A-242, Jahangir Puri, Delhi - 110033

Phone: 8130973844

Email: architasst@gmail.com

QUALITY EQUIPMENT.
TRUSTED RESULTS.
STRONGER INFRASTRUCTURE.

For Better Testing Today,
For Better Tomorrow.



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844 | E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



ULTRASONIC THICKNESS GAGE (BASIC TYPE) CODE ISU-100D

ATTENTION: NOT SUITABLE FOR
CASTINGS AND ALUMINIUM WORKPIECES



- Measure the thickness from one side of objects, suitable for pipes, tanks, etc.
- Applicable material: steel, stainless steel, glass, copper, brass, polystyrene, titanium, etc.

SPECIFICATION (ON STEEL)

Measuring range		0.8 - 300 mm
Resolution		0.1 / 0.01 mm
Accuracy		±0.05mm (H<10mm) ±(0.05+H/1000)mm (10≤H<100mm) ±H/333mm (H≥100mm) H is the thickness to be measured in mm
Transducer	Frequency	5MHz
	Diameter (Ød)	10.8 mm
Minimum size of pipes		20 × 1.2 mm (diameter × wall thickness)
Velocity		1000-9999 m/s
Operation temperature		Main unit: -20-50°C Transducer: -20-60°C
Measuring unit		mm/inch
Power supply		2 × AAA batteries
Dimension		114 × 64 × 28 mm
Weight		200 g



LCD with backlight
for easy reading



High measurement accuracy
and stability



4.00mm block
for calibration



Easy to operate
and portable



Rugged enclosure for
protection & durability

APPLICATION EXAMPLE



STANDARD DELIVERY

Main unit	1 pc
Transducer	1 pc
Couplant	1 bottle
Battery (AAA)	2 pcs

PACKING

- Supplied in rugged protective enclosure with foam insert.
- Dimension (Enclosure): 192 x 120 x 60 mm (Approx.)
- Gross Weight: 0.65 kg (Approx.)



ACCURATE
MEASUREMENTS



RELIABLE
PERFORMANCE



QUALITY
ASSURED



PORTABLE &
DURABLE DESIGN



IDEAL FOR FIELD
INSPECTIONS



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844

E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



VERNIER CALIPER DIGITAL & ANALOG

PRECISE MEASUREMENT. RELIABLE PERFORMANCE.

Vernier Calipers are precision measuring instruments used for measuring internal, external, depth and step dimensions with high accuracy. Available in Digital and Analog models in various sizes to suit all industrial, laboratory and quality control applications.

DIGITAL VERNIER CALIPER



Make: INSIZE

ANALOG VERNIER CALIPER



Make: Mitutoyo

FEATURES

- High Accuracy**
Ensures precise and reliable measurement results.
- Hardened Stainless Steel**
Made from high quality stainless steel for long life and corrosion resistance.
- Four Way Measurement**
Measures external, internal, depth and step dimensions.
- Easy to Read**
Digital display for easy reading. Analog scale is clear and crisp.
- Smooth Operation**
Smooth sliding mechanism for comfortable and accurate operation.
- Data Output (Digital Models)**
With SPC data output port (for compatible models).

AVAILABLE SIZES & SPECIFICATIONS

Size (mm)	DIGITAL VERNIER CALIPER			ANALOG VERNIER CALIPER		
	Range (mm)	Accuracy (mm)	Resolution (mm)	Range (mm)	Accuracy (mm)	Resolution (mm)
150	0 - 150	± 0.02	0.01	0 - 150	± 0.02	0.02
300	0 - 300	± 0.03	0.01	0 - 300	± 0.03	0.02
600	0 - 600	± 0.05	0.01	0 - 600	± 0.05	0.02
1000	0 - 1000	± 0.07	0.01	0 - 1000	± 0.07	0.02
1600	0 - 1600	± 0.10	0.01	0 - 1600	± 0.10	0.02
2000	0 - 2000	± 0.12	0.01	0 - 2000	± 0.12	0.02

* Accuracy as per factory standard. (Excluding quantizing error)

DIGITAL VERNIER CALIPER - INSIZE

- Large LCD display
- Zero setting at any position
- Metric / Inch conversion
- Absolute / Incremental measurement
- Data output (for models with SPC port)
- Made of stainless steel
- Supplied in protective case



ANALOG VERNIER CALIPER - MITUTOYO

- Satin chrome finish on scale for glare-free reading
- High accuracy and repeatability
- Smooth sliding and fine adjustment
- Inside, outside, depth and step measurement
- Made of hardened stainless steel
- Supplied in protective case



APPLICATIONS

- Machining & Manufacturing
- Quality Control Laboratories
- Automotive & Aerospace
- Engineering Workshops
- Education & Training Institutes

STANDARD ACCESSORIES

- Protective Case
- Battery (Digital Models)
- Calibration Certificate
- Instruction Manual



OPTIONAL ACCESSORIES (DIGITAL MODELS)

- SPC Data Cable
- Wireless Data Transmitter
- USB Data Interface
- Bluetooth Data Output



PRECISE
MEASUREMENT



DURABLE &
RELIABLE



CORROSION
RESISTANT



EASY TO USE



WIDE RANGE
OF SIZES



GLOBAL
QUALITY
BRANDS

ARCHIT ASSOCIATES

A-242, Jahangir Puri,
Near Kali Mandir, Delhi-110 033 (INDIA)

91-8130973844

architasst@gmail.com

GSTIN NO. 07A0UPG0300F2ZC

ARCHIT ASSOCIATES - YOUR PARTNER IN TESTING EXCELLENCE



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844

E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07AOUPG0300F2ZC



DIGITAL & ANALOG VERNIER CALIPER

PRECISE MEASUREMENT. RELIABLE RESULTS.

High accuracy vernier calipers for inside, outside, depth and step measurements. Available in Digital and Analog models from trusted brands INSIZE and MITUTOYO. Ideal for laboratories, industries, workshops and quality control applications.

Available Sizes (All Models)

150 mm | 300 mm | 600 mm | 1000 mm |
1600 mm | 2000 mm

DIGITAL VERNIER CALIPER (INSIZE)



- Resolution: 0.01 mm / 0.0005"
- Accuracy: ± 0.02 mm
- Large LCD Display
- Zero Setting at Any Position
- Metric / Inch Conversion
- Hardened Stainless Steel
- Smooth Movement & Long Life



DIGITAL VERNIER CALIPER (MITUTOYO)



- Resolution: 0.01 mm / 0.0005"
- Accuracy: ± 0.02 mm
- ABSOLUTE Encoder Technology
- No Origin Reset Needed
- Large LCD Display
- IP67 Coolant & Dust Proof
- High Quality Stainless Steel



ANALOG VERNIER CALIPER (INSIZE / MITUTOYO)

- Reading: 0.02 mm / 0.001"
- Accurate Vernier Scale
- Hardened Stainless Steel
- Fine Adjustment
- Satin Chrome Finish
- Available in INSIZE & MITUTOYO



STANDARD SPECIFICATIONS

Model Type	150 mm / 6"	300 mm / 12"	600 mm / 24"	1000 mm / 40"	1600 mm / 60"	2000 mm / 80"
Measuring Range	0 - 150 mm	0 - 300 mm	0 - 600 mm	0 - 1000 mm	0 - 1600 mm	0 - 2000 mm
Resolution (Digital)	0.01 mm / 0.0005"	0.01 mm / 0.0005"	0.01 mm / 0.0005"	0.01 mm / 0.0005"	0.01 mm / 0.0005"	0.01 mm / 0.0005"
Resolution (Analog)	0.02 mm / 0.001"	0.02 mm / 0.001"	0.02 mm / 0.001"	0.02 mm / 0.001"	0.02 mm / 0.001"	0.02 mm / 0.001"
Accuracy (Digital)	± 0.02 mm	± 0.02 mm	± 0.02 mm	± 0.02 mm	± 0.02 mm	± 0.02 mm
Accuracy (Analog)	± 0.03 mm	± 0.03 mm	± 0.04 mm	± 0.04 mm	± 0.04 mm	± 0.06 mm
Depth Measurement	✓	✓	✓	✓	✓	✓
Inside / Outside / Step	✓	✓	✓	✓	✓	✓
Material	Hardened Stainless Steel					

KEY FEATURES

- ✓ High Accuracy & Repeatability
- ✓ Smooth Sliding Mechanism
- ✓ Inside, Outside, Depth & Step Measurement
- ✓ Metric / Inch Conversion (Digital)
- ✓ Strong, Durable & Long Life

APPLICATIONS

- ✓ Engineering & Manufacturing
- ✓ Quality Control & Inspection
- ✓ Research & Development
- ✓ Education & Training Institutes
- ✓ Laboratories & Workshops

BRANDS AVAILABLE



Mitutoyo
WORLD CLASS PRECISION
YOU CAN TRUST

ARCHIT ASSOCIATES

A-242, Jahangir Puri,
Near Kali Mandir, Delhi-110 033 (INDIA)

91-8130973844

architasst@gmail.com

GSTIN NO. 07AOUPG0300F2ZC

ARCHIT ASSOCIATES - YOUR PARTNER IN TESTING EXCELLENCE



ARCHIT ASSOCIATES

Trusted Supplier of
MATERIAL TESTING & LABORATORY INSTRUMENTS

A-242, Jahangir Puri, Near Kali Mandir,
Delhi-110 033 (INDIA)

architasst@gmail.com

91-8130973844, 9871406036,
9810052732, 9810052758

www.architassociates.com

GSTIN NO. 07A0UPG0300F22C

TRUSTED SUPPLIER OF
MATERIAL TESTING &
LABORATORY INSTRUMENTS



QUALITY
ASSURED
PRODUCTS



ACCURATE &
RELIABLE



LABORATORY &
SITE SOLUTIONS



CUSTOMIZED
TESTING
EQUIPMENT

VIBRATING MACHINE FOR CEMENT MORTAR

IS : 10080

The machine is capable for compaction of cement mortar cubes as per IS: 4031 and designed as per Standard IS : 10080.

The complete frame assembly consists of device to hold a 7.06 cm cube mould to feed the sample in the mould.

The frequency of vibration is 12000 ± 400 R.P.M. Vibration per minute.

Supplied complete with one 7.06 cm cube mould and a DIGITAL Timer.

Complete with One No. poking rod to meet IS : 10080 requirements.




TECHNICAL SPECIFICATIONS

Standard	IS : 10080
Mould Size	7.06 cm Cube
Frequency of Vibration	12000 ± 400 R.P.M. (Vibration per minute)
Drive	Electric Motor
Timer	Digital Timer
Supplied With	+ 1 No. 7.06 cm Cube Mould + 1 No. Poking Rod
Power Supply	230V AC, 50Hz

COMPLETE WITH



7.06 cm Cube Mould

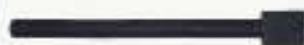

Poking Rod

APPLICATION



Used for compaction of cement mortar cubes as per IS: 4031 and IS: 10080.

FEATURES

- ✓ Sturdy and rigid construction
- ✓ Exact frequency as per standard
- ✓ Easy to operate and maintain
- ✓ Low noise and vibration
- ✓ Long service life

WIDELY USED IN



CEMENT
INDUSTRY



CONSTRUCTION
QUALITY CONTROL



RESEARCH
LABORATORIES



ENGINEERING
COLLEGES



TESTING
LABORATORIES

A-242, Jahangir Puri,
Near Kali Mandir,
Delhi-110 033 (INDIA)



91-8130973844
9871406036
9810052732
9810052758



architasst@gmail.com

ALL PRODUCTS COMPLY WITH RELEVANT INDUSTRY STANDARDS



ARCHIT ASSOCIATES

YOUR TRUSTED PARTNER IN LABORATORY & MATERIAL TESTING

📍 A-242, Jahangir Puri, Near Kali Mandir,
Delhi-110 033 (INDIA)

☎ 91-8130973844 ✉ architassst@gmail.com

🌐 www.architassociates.com

GSTIN NO. 07AOUPG0300F2ZC



ACCURATE
RESULTS



RELIABLE
PERFORMANCE



QUALITY
ASSURED

CONCRETE VIBRATING TABLE

FOR COMPACTING CONCRETE SPECIMENS

Concrete Vibrating Table is used for compacting concrete in moulds, cubes, cylinders, beams and other precast elements. It removes entrapped air and provides a uniform, dense and high quality concrete specimen.

KEY FEATURES

- ✓ Robust and rigid construction
- ✓ High compaction efficiency
- ✓ Uniform vibration for consistent results
- ✓ Heavy-duty vibration motor
- ✓ Easy to operate and maintain
- ✓ Suitable for laboratory and field use



TECHNICAL SPECIFICATIONS

PARAMETER	SPECIFICATION
Table Size	500 x 500 mm / 750 x 750 mm / 1000 x 1000 mm (custom options available)
Top Plate	Mild Pteel, machined flat, with anti-corrosive powder coating
Vibration Motor	Heavy-duty motor, 1 HP / 1.5 HP, 2800 RPM
Vibration Frequency	3600 vibrations per minute (adjustable)
Control Panel	On/Off switch with vibration timer (optional)
Power Supply	230V, 50 Hz single phase / 415V, 3 phase
Overall Construction	Heavy-duty welded steel frame with spring mounting
Finish	Anti-corrosive powder coated
Applications	Compacting of concrete cubes, cylinders, beams, pavers, tiles and other precast elements

APPLICATIONS

- Concrete testing laboratories
- Precast concrete factories
- Research & development laboratories
- Construction quality control
- Educational institutions



STRONG
CONSTRUCTION



HIGH
EFFICIENCY



CONSISTENT
RESULTS



LOW
MAINTENANCE

NOTE Specifications are subject to change without prior notice.



91-8130973844



architassst@gmail.com



A-242, Jahangir Puri,
Near Kali Mandir,
Delhi-110 033 (INDIA)

ARCHIT ASSOCIATES – YOUR PARTNER IN TESTING EXCELLENCE



ARCHIT ASSOCIATES

YOUR TRUSTED PARTNER IN LABORATORY & MATERIAL TESTING

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

91-8130973844, 9871406036 | architassst@gmail.com

www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



ACCURATE
RESULTS



RELIABLE
PERFORMANCE



QUALITY
ASSURED

TESTING EXCELLENCE IN EVERY RESULT

STAINLESS STEEL WATER DISTILLATION UNIT

4 LITRE / HOUR CAPACITY

High quality Stainless Steel Water Distillation Unit designed to produce pure, pyrogen free distilled water, ideal for laboratories, hospitals, clinics and research applications.

KEY FEATURES

- Made of high grade stainless steel (SS 304) for long life & corrosion resistance
- Produces high purity distilled water, pyrogen free
- Efficient heating system for faster distillation & low power consumption
- Sturdy construction with elegant finish
- Easy to install, operate & maintain
- Suitable for continuous operation



CAPACITY
4
LITRE/HOUR

APPLICATIONS

- Laboratories & Research Institutes
- Hospitals & Clinics
- Pharmaceutical Industries
- Educational Institutions
- Quality Control Labs



PURE & SAFE
DISTILLED WATER



RUST PROOF
& DURABLE



ENERGY EFFICIENT
PERFORMANCE



BUILT FOR
LONG LIFE

MANESTY
TYPE STILL
LABORATORY GRADE

TECHNICAL SPECIFICATIONS

Capacity	4 Litre / Hour (Approx.)
Construction Material	Stainless Steel
Type	Manesty Type
Water Consumption	Approx. 80 to 100 Lt./Hr.
Power Supply	230V AC, 50 Hz
Power Rating	3.0 kW (Approx.)
Distilled Water Outlet	1 No.
Cooling Water Inlet	1 No.
Drain Outlet	1 No.
Overall Dimensions (Approx.)	450 mm (L) x 250 mm (W) x 900 mm (H)

STANDARD ACCESSORIES INCLUDED

- SS Distillation Unit (Manesty Type)
- Stainless Steel Condenser
- Level Control Arrangement
- Inlet & Outlet Connections
- Instruction Manual

OPTIONAL ACCESSORIES

- Wall Mounting Kit
- Automatic Water Level Controller
- Tubing & Fittings
- Extra Condenser



PRECISE RESULTS
Consistent & accurate
performance



SUPERIOR QUALITY
High grade components
for reliability



EASY TO OPERATE
User friendly design
& simple controls



TEST CERTIFICATES
Supplied with
calibration certificates



AFTER SALES SUPPORT
Prompt service &
technical support



SCAN TO VISIT
OUR WEBSITE

ARCHIT ASSOCIATES

A-242, Jahangir Puri, Delhi - 110033 (INDIA)

91-8130973844, 9871406036

architassst@gmail.com

28
YEARS OF
EXPERIENCE



BUILT FOR LABS.
ENGINEERED FOR EXCELLENCE.
Precision. Performance. Perfection.



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir, Delhi-110 033 (INDIA)

Phone: 91-8130973844, E-mail: architasst@gmail.com

Web: www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



SCAN TO VIEW
PRODUCT RANGE

COMPLETE RANGE OF WATER TESTING EQUIPMENT & ACCESSORIES

ACCURATE • RELIABLE • EASY TO USE

High quality instruments and solutions for accurate water quality analysis in laboratories, industries, wastewater treatment plants and field applications.



HIGH ACCURACY
& RELIABILITY



WIDE RANGE OF
APPLICATIONS



EASY TO USE



PRECISE &
FAST RESULTS



RUGGED &
DURABLE DESIGN



COMPLETE RANGE
OF ACCESSORIES



LAB pH METER

LAB TDS METER

LAB ORP METER

POCKET TYPE METERS



POCKET
pH METER



POCKET
TDS METER



POCKET
ORP METER



POCKET
EC METER



POCKET
SALINITY METER



COMPLETE
TESTING KIT
Portable • Comprehensive • Reliable

DISPOSABLE TESTERS



pH
TEST KIT
pH 4.0 - 9.0
100 TESTS



TDS
TEST KIT
0 - 2000 µg/L
100 TESTS



CHLORINE
TEST KIT
0 - 5 ppm
100 TESTS



HARDNESS
TEST KIT
0 - 500 ppm
100 TESTS



IRON
TEST KIT
0 - 10 ppm
100 TESTS



NITRATE
TEST KIT
0 - 50 ppm
100 TESTS



REAGENTS & SOLUTIONS

ACCESSORIES



pH ELECTRODE



TDS PROBE



ORP PROBE



CALIBRATION
SOLUTIONS
(pH 4.01 / 7.00 / 9.18)



BUFFER SOLUTIONS
(pH / ORP / TDS)



STANDARD
SOLUTIONS



KCL SOLUTION &
STORAGE BOTTLE

APPLICATIONS



DRINKING WATER
TESTING



WASTEWATER
TREATMENT



INDUSTRIAL WATER
TESTING



BOILER & COOLING
WATER



POOL & SPA
WATER



RO / DI WATER
TESTING



ENVIRONMENTAL
MONITORING



ACCURATE
MEASUREMENTS
Reliable results you
can trust



PREMIUM QUALITY
CONSTRUCTION
Built for laboratory &
field use



USER FRIENDLY
OPERATION
Intuitive design for
every user



AFTER SALES SUPPORT
& CALIBRATION
Expert support for long
life performance

ARCHIT ASSOCIATES

A-242, Jahangir Puri,
Near Kali Mandir, Delhi-110 033

91-8130973844

architasst@gmail.com



ARCHIT ASSOCIATES

A-242, Jahangir Puri, Near Kali Mandir,
Delhi-110 033 (INDIA)

91-8130973844

architasst@gmail.com

www.architassociates.com

GSTIN NO. 07A0UPG0300F2ZC



TRUSTED SUPPLIER OF
MATERIAL TESTING &
LABORATORY INSTRUMENTS



QUALITY
ASSURED
PRODUCTS



ACCURATE &
RELIABLE



LABORATORY &
SITE SOLUTIONS



CUSTOMIZED
TESTING
EQUIPMENT

PREMIUM WEIGHING BALANCES

High Accuracy • Reliable • Advanced Technology

SARTORIUS WEIGHING BALANCE



- High precision electronic balance
- Auto calibration
- RS232 / USB interface
- Ideal for laboratory & research use

MODEL	CAPACITY	READABILITY	ACCURACY
CP Series	160 g	0.1 mg	±0.1 mg
CP Series	320 g	0.1 mg	±0.1 mg
CP Series	620 g	0.1 mg	±0.1 mg
CP Series	2200 g	1 mg	±1 mg
CP Series	6200 g	1 mg	±1 mg

sartorius

ACZET WEIGHING SCALE



- High accuracy weighing scale
- Stable & durable performance
- RS232 interface
- Ideal for lab & field use

MODEL	CAPACITY	READABILITY	ACCURACY
ACN 220G	220 g	0.1 mg	±0.1 mg
ACN 320G	320 g	0.1 mg	±0.1 mg
ACN 520G	520 g	0.1 mg	±0.1 mg
ACN 2200G	2200 g	1 mg	±1 mg
ACN 3200G	3200 g	1 mg	±1 mg
ACN 6200G	6200 g	1 mg	±1 mg

aczet

SACALATECH WEIGHING SCALE



- Heavy duty industrial scale
- Stainless steel platform
- RS232 interface
- Easy operation & long life

MODEL	CAPACITY	READABILITY	ACCURACY
SAB 302	300 g	0.01 g	±0.02 g
SAB 602	600 g	0.01 g	±0.02 g
SAB 1002	1100 g	0.01 g	±0.02 g
SAB 1502	1500 g	0.01 g	±0.02 g
SAB 3KE	3100 g	0.1 g	±0.1 g
SAB 6KE	6100 g	0.1 g	±0.2 g

SACALATECH
TECHNOLOGY MEETS PRECISION

COMPREHENSIVE WEIGHING SOLUTIONS

ANALYTICAL BALANCE (SCALETEC - AXIS RANGE)



MODEL	CAPACITY	READABILITY	ACCURACY
APN 802200	80220 g	0.0001 mg	±0.0001 mg
APN 2200	220 g	0.1 mg	±0.1 mg
APN 2205	220 g	0.1 mg	±0.1 mg



Auto Internal
Calibration



High Repetitive
Compliance

PRECISION BALANCE (SCALETEC - AXIS RANGE)



MODEL	CAPACITY	READABILITY	ACCURACY
ACA 320G	320 g	0.001 g	±0.001 g
ACA 520G	520 g	0.001 g	±0.001 g
ACA 820G	820 g	0.001 g	±0.002 g
ACA 3200G	3200 g	0.01 g	±0.01 g
ACA 6200G	6200 g	0.01 g	±0.02 g



Auto Internal
Calibration



Large Graphic
Display



Database &
Statistics

COMPACT BALANCE (SCALETEC - AXIS RANGE)



MODEL	CAPACITY	READABILITY	ACCURACY
ATN 120G	120 g	0.1 mg	±0.1 mg
ATN 100G	100 g	0.1 mg	±0.1 mg
ATN 2200	220 g	0.1 mg	±0.1 mg



Compact
Design



Fast
Stabilization



User
Friendly

HIGH CAPACITY BALANCE (SCALETEC - AXIS RANGE)



MODEL	CAPACITY	READABILITY	ACCURACY
SAB 2200CL	2200 g	0.01 g	±0.02 g
SAB 3000CL	3000 g	0.01 g	±0.03 g
SAB 5200CL	5200 g	0.01 g	±0.03 g
SAB 6200CL	6200 g	0.01 g	±0.03 g



High Capacity



Strong Build



Reliable Accuracy

KEY FEATURES

- ✓ High Accuracy & Repeatability
- ✓ Multiple Weighing Units
- ✓ RS232 / USB / LAN Connectivity
- ✓ GLP / GMP / ISO Compliant
- ✓ Counting, Percentage, Statistics
- ✓ Internal / External Calibration

WIDELY USED IN



Laboratories



Pharmaceutical
Industry



Research &
Development



Quality Control
Labs



Education &
Institutes

ARCHIT ASSOCIATES
A-242, Jahangir Puri,
Near Kali Mandir,
Delhi-110 033 (INDIA)

91-8130973844

architasst@gmail.com

ALL PRODUCTS COMPLY WITH RELEVANT INDUSTRY STANDARDS